


Quality Products for Dentistry Worldwide
Made in USA



Equipment · **Components** · **Tubing**

Dear Colleagues,

As in previous years, we continue to grow and improve our product line to meet the needs of the dental community worldwide. The 2017-2018 updated catalog features several new products, identified by  symbol throughout the catalog, as well as a new line of Orthodontic Equipment.

Since 1996, our goal and mission have remained the same – to establish a long-term relationship and to provide you with quality, American made products . As always, we strive to provide you with competitive prices accompanied with professional and reliable service.

We are honored for your continued patronage and look forward to working with you for many years to come.

GENERAL ORDERING INFORMATION

For pricing quotations or orders, please include the following information in your correspondence:

- Catalog number
- Description
- Quantity
- Color, if applicable
- Method of shipping preferred

TERMS

- For Domestic clients – Advance payment or net 30 days with approved credit.
- For International clients – Advance payment by wire transfer or confirmed, irrevocable letter of credit, at sight.
- Credit cards are acceptable, 4% transaction fee applies.
- All prices are subject to change without notice.

RETURN/EXCHANGE MERCHANDISE

- All returns must have prior company approval.
- Indicate catalog number, invoice number and reason for return.
- Unless defective, a re-stocking fee of 15% may be charged.
- Carefully package all returning items to prevent shipping damage.
- Return shipping must be freight prepaid.
- Custom orders are final and may not be returned for credit.

WARRANTY

DentalPartsUSA (a division of U.S. MEDLINK) warrants its products to be free from defect in material and workmanship only. No other warranties are expressed or implied. Warranty periods are indicated in each section of the catalog. Written notice of breach must be given to DentalPartsUSA within this period. Warranty is limited to repair parts and replacement of items by DentalPartsUSA. This warranty does not cover damage resulting from improper installation or maintenance, accident or misuse. This warranty does not cover damage that results from using cleaning, disinfecting or sterilization chemicals and processes. No claim for labor or consequential damages will be allowed.

Contact Information

DentalPartsUSA
418 North Glendale Avenue, Unit J
Glendale, California 91206 U.S.A.

Telephone: 1-818-502-9954
Facsimile: 1-818-502-0501
Website: www.dentalpartsusa.com

DENTAL EQUIPMENT	7 - 37
DENTAL COMPONENTS	41 - 125
DENTAL TUBING	129 - 139
NUMERICAL INDEX PAGES	141



Dental Chairs / Stools / Color Chart	7
Orthodontic Equipment	11
Over The Patient Delivery Units	14
Beside Patient Delivery Units Beneath Shelf Mount Delivery Units	16
Cabinet Mount Delivery Units	17
Doctor's Mobile Delivery Units	21
Handpiece Control Heads Panel Mount Handpiece Controls	24
Wall Mount Handpiece Control Heads	28
Vacuum Arm Systems	30
Cuspidors	33
Operatory Lights	34
Operatory Products / Accessories	35

DentalPartsUSA strives to offer dental products which are made in USA. The dental equipment featured in this section are made in American facilities that have been producing quality products for domestic and international dental practices for decades. We offer affordable dental units that are simple, yet durable and attractive. Compared to expensive, electronically driven dental units, our units are known for the practical, pneumatic features requiring minimal maintenance for years of continuous operation. The patient chairs range from deluxe to standard to economy models, which can be customized with numerous optional features. We provide a warranty in material and workmanship of up to 5 years on dental chairs and units, and 3 years on stools. The California Series dental patient chairs have a 10-year warranty on motors. Ask us about customizing the units and chairs to fit your needs.

Standard Upholstery Color Chart

Boltalex® Upholstery 100% Vinyl Face, Poly/Cotton Backing with PreFixx(R) protective coating for increased stain resistance.

Also available - ULTRA LEATHER Upholstery (additional charge applies)



California Series Dental Patient Chairs



Model USM-CA 550

Retail Price \$4,200

Model USM-CA 450

Retail Price \$3,900

Model USM-CA 350

Retail Price \$3,500

Standard Features

	USM-550	USM-450	USM-350
All steel construction covered with powder coated plastic eliminating chipping	✓	✓	✓
Chair frame made of steel	✓	✓	
Chair frame covered with durable plastic		✓	
Bottom heavy - no bolt down required	✓	✓	✓
Electro-mechanical with 800 lbs of lift capacity (primary lift motor made in Germany, secondary back-rest motor made in USA)	✓	✓	✓
Standard articulating headrest included	✓	✓	✓
Ultra thin seat back	✓	✓	✓
Memory foam cushions, ideal for longer procedures	✓		
Asepsis seamless cushions		✓	
Removable/reversible cushions			✓
Hand controls operation on both sides of chair	✓	✓	✓
Auto return and adjustable preset (Auto return only on USM- CA 350)	✓	✓	✓
Wide back-rest	✓	✓	✓
Winged arm covers for elbow support			✓
Drop down arms for easy access	✓	✓	
Available in 110V and 220V	✓	✓	✓
Over 150 color choices	✓	✓	✓
10-YEAR WARRANTY	✓	✓	✓

Optional Features for all model chairs

Chair foot control	add \$460
Push button articulating headrest	add \$440
Chair swivel	add \$425

Doctor's & Assistant's Stools



Model DR-96B Basic Doctor's Stool Retail Price \$342	Model DR-96T Traditional Doctor's Stool Retail Price \$739	Model AT-96B Basic Assistant's Stool Retail Price \$466	Model N-965 Basic Operator's Stool Retail Price \$337
---	---	--	--

Standard Features	DR-96B	DR-96T	AT-96B	N-965
Stool: Triple Lever Adjustment		✓		
Stool: Single Lever Adjustment	✓		✓	✓
Back: Oval Back	✓			
Back: D-Shaped		✓		
Seat: Contoured	✓			
Seat: Round		✓	✓	
Seat: Round Vacuum-formed seat				✓
Arm: T-Stem			✓	
Cylinder: Standard - Metal Base	✓			
Cylinder: Standard - Plastic Base		✓	✓	
Cylinder: Chrome pneumatic				✓
Foot Ring: Fixed			✓	
Base: Metal - Fully Powder Coated	✓			
Base: Plastic Base		✓	✓	
Base: Aluminum - Taupe				✓
Casters: Carpet - Beige	✓	✓	✓	
Casters: Hard floor - Taupe				✓
Cushioning: Standard Foam	✓	✓	✓	✓
Upholstery: Standard Vinyl	✓	✓	✓	✓

Dental Equipment Package

- Asepsis Patient Chair
- 2 HP Automatic Air Driven Delivery Unit with syringe
- Economy Doctor Stool
- Halogen Operatory Light (or optional - LED Operatory Light)
- Assistant Delivery (or with Cuspidor Delivery System)
- All products are Made-in-USA (Halogen Light only – assembled in USA)
- Available in White or Surf powder coat
- Warranty Period: 5 years - (Except: Halogen Light – 1 year, LED Light – 2 years)
- **Factory ISO 9000, ETL, UL and CSA Certified



Model INT-2001

Retail Price \$7,520

Standard Dental Chair, Asepsis Upholstery

- Electro-mechanical Cantilever system
- Metal construction with HiGlos powder-coat inside and out
- Asepsis quick-change upholstery - choice of 10 standard colors
- 7-way foot control with two programmable preset positions and automatic exit/return
- Wheelchair accessible double articulating headrest
- Synchronized seat and back tilt
- Left/right access with swing up arms
- Locking swivel, 60 degree
- Low voltage
- Available in 110V and 220V
- Made in USA
- Five-year warranty



Model 2000

Retail Price \$5,240

Economy Dental Chair, Asepsis Upholstery

- Electro-mechanical Cantilever system
- Metal construction with HiGlos powder-coat inside and out
- Asepsis upholstery - choice of 10 standard colors
- 5-way foot control with automatic exit/return
- Wheelchair accessible double articulating headrest
- Synchronized seat and back tilt
- Fixed arms
- Thermal overload protection
- Low voltage
- Available in 110V and 220V
- Made in USA
- Five-year warranty



Model 1000

Retail Price \$4,827

Orthodontic Chairs



shown with double articulating headrest, plush upholstery & trim

2000-060 Wide Back Asepsis Upholstery Retail Price \$3,378	2000-060SA Sling Back Asepsis Upholstery Retail Price \$3,492	2000-060T Plush Ultraleather Upholstery & Trim Retail Price \$3,986	2000-060S Sling Back Plush Ultraleather Upholstery & Trim Retail Price \$3,986	2000-060E ESense Plush Narrow Back Ultraleather Uphols. & Trim Retail Price \$3,986
---	--	--	---	--

Standard Features

- Electromechanical backrest
- Pedestal base with 5" manual height adjustment
- Foot operation switch with automatic exit (return), mounted on pedestal
- Non-articulating headrest
- All metal construction with HiGlos powdercoat inside and out
- Synchronized seat and back tilt
- No armrests (optional add on)
- Warranty: 5 years
- Quick change upholstery

Optional Features for Chairs above:

Catalog No.	Description	Upgrade Charge
7000-021	Double Articulating Headrest in lieu of standard	\$114
2000-243	Stationary Arms	\$228
2000-244	Moveable Arms (arms swing up)	\$404
7000-96-1.1	Asepsis Upholstered Arm Covers, per set	\$108
2000-226	Chair Mount Adapter	\$286
1014	Slings for ESense Upholstery (For chair model 2000-060E only)	\$116

Standard Ortho Cart



shown with
2HP unit &
VAC package

CPU Ortho Cart



shown with
3HP unit &
VAC package



Model USM-SP200

Standard Ortho Cart, without delivery unit

Retail Price: \$3,120

Model USM-SP400

CPU Ortho Cart, without delivery unit

Retail Price: \$3,660

Standard Features

	USM SP200	USM SP400
Choice of Casters/Mobile (CA) or Fixed (FB) Base	✓	✓
Five drawer module, with four shallow drawers and one deep drawer	✓	✓
Rear access door	✓	
Access door for CPU		✓
Keyboard pullout		✓
Choose from 6 standard laminate tops and 10 standard cabinet colors	✓	✓
Warranty: 1 year	✓	✓

Optional Features

FPM103G	Flat panel monitor mount (monitor not included)	\$ 408
1007	Install of Panel Mount in Cart (each install)	\$ 66

Standard Work Surface Laminates



Standard Cabinetry Colors



Ortho Panel Mount HP Units



Model USM-4062 2HP Manual Panel Mount Delivery Unit Retail Price: \$972	Model USM-4063 3HP Automatic Panel Mount Delivery Unit Retail Price: \$1,678
--	---

Standard Features	USM-4062	USM-4063
Panel mount 2 handpiece manual control unit with:	✓	
Panel mount 3HP automatic control unit		✓
4" non-swivel holder bar	✓	
8" non-swivel holder bar		✓
3-way syringe with autoclavable tip	✓	✓
Standard disc-type, chrome plated foot control	✓	✓
Junction box accessories (no box), includes air/water shut off valves, regulators & filters	✓	✓
Warranty: 5 years	✓	✓

Optional Features for Panel Mount Delivery Units:

Catalog No.	Description	Upgrade Charge
4091-7	Vacuum package with HVE & SE valves & tubing	\$212
4091	Vacuum package (HVE, SE, SYR)	\$344
4003-193-1	Panel mount clean water bottle (bottle water only)	\$186
4095-4	Panel mount clean water system (city/bottle switch)	\$370
4064	Swivel holder bar	\$152
4003-276	Additional 4" holder bar	\$ 36
4003-275	Additional 8" holder bar	\$ 54
5025-1	Junction box with lid (no accessories)	\$144
4003-240	1-1/4" Spiral flex vacuum tubing, per foot	\$ 8
4003-123	2" Spiral flex vacuum tubing, per foot	\$ 10
4000-055	Additional umbilical per foot	\$ 32

Over The Patient Delivery Units

Doctor's Delivery Units with Pneumatic Brake Flex Arm



Post NOT included



Top Post Mount Units include 28\"/>

Model USM-6000

3HP Automatic Euro 2" dia. Post Mount

Retail Price \$3,800

Model USM-6010

3HP Automatic Euro Top Post Mount

Retail Price \$3,800

Standard Features

	USM-6000	USM-6010
Asepsis designed automatic selection for 3 handpieces	✓	✓
Side-mounted recessed control adjustments	✓	✓
3-way syringe with quick change autoclavable tip	✓	✓
Master On/Off toggle with indicator	✓	✓
Pressure adjustment for each handpiece	✓	✓
Easy to read pressure gauge	✓	✓
Air coolant flow adjustment	✓	✓
Non-retracting water coolant adjustment & purge	✓	✓
Asepsis handpiece tubing, sterling	✓	✓
Stainless steel instrument tray, 9¾" x 13½"	✓	✓
Arm-mounted tray holder, 30" lateral movement, 360° rotation	✓	✓
Asepsis flex arm system with pneumatic brake, 52" reach, 27" vertical travel	✓	✓
Unit mount post, 28"		✓
Standard foot control	✓	✓
Utility center assembly, 14" x 9" x 5⅜" with master valves and regulators	✓	✓
7' umbilical tubing	✓	✓

Optional Features

	Retail Price	USM-6000	USM-6010
Delete utility center	deduct \$400	✓	✓
Additional umbilical tubing, per foot	add \$11	✓	✓
Wet/dry foot control in lieu of standard	add \$100	✓	✓
Chip blower foot control in lieu of standard	add \$100	✓	✓
Wet/dry chip blower foot control in lieu of standard	add \$150	✓	✓
Unit mount post, 40" with bushing	add \$165	✓	
Light post extension with butterfly clamp, post cap, light cord and sleeve	add \$250		✓
Chair bracket (specify model of chair)	add \$290	✓	✓

Over The Patient Delivery Units

Doctor's Delivery Units with Pneumatic Brake Flex Arm



Model USM-3300
3HP Automatic 2" dia. Post Mount

Model USM-3310
3HP Automatic Top Post Mount

Retail Price \$3,295

Retail Price \$3,295

Standard Features

	USM-3300	USM-3310
Asepsis designed automatic selection for 3 handpieces	✓	✓
Side-mounted recessed control adjustments	✓	✓
3-way syringe with quick change autoclavable tip	✓	✓
Master On/Off toggle with indicator	✓	✓
Pressure adjustment for each handpiece	✓	✓
Easy to read pressure gauge (protected under cover)	✓	✓
Air coolant flow adjustment	✓	✓
Non-retracting water coolant adjustment & purge	✓	✓
Asepsis handpiece tubing, sterling (specify straight or coiled)	✓	✓
Stainless steel instrument tray, 9¾" x 13½"	✓	✓
Asepsis flex arm system with pneumatic brake, 52" reach, 27" vertical travel	✓	✓
Unit mount post, 28"		✓
Standard foot control	✓	✓
Utility center assembly, 14" x 9" x 5⅜" with master valves and regulators	✓	✓
7' umbilical tubing	✓	✓

Optional Features

	Retail Price		
Delete utility center	deduct \$400	✓	✓
Delete handpiece station with tubing	deduct \$50	✓	✓
Additional handpiece station with tubing (dry only)	add \$150	✓	✓
Silcryn™ Tubing	add \$20	✓	✓
Additional umbilical tubing, per foot	add \$11	✓	✓
Wet/dry foot control in lieu of standard	add \$100	✓	✓
Chip blower foot control in lieu of standard	add \$100	✓	✓
Wet/dry chip blower foot control in lieu of standard	add \$150	✓	✓
Unit mount post, 40" with bushing	add \$165	✓	
Light post extension with butterfly clamp, post cap, light cord and sleeve	add \$250		✓
Chair bracket (specify model of chair)	add \$290	✓	✓

Beside The Patient Delivery Unit

Hygiene System with Telescoping Arm and Pivoting Tray Mount Assembly



Model USM-5300
2HP Automatic 2" dia. Post Mount

Retail Price \$1,910

Beneath Shelf Mount Delivery Unit

Doctor's Delivery Unit



Filter/Regulators must be ordered separately

Model USM-250
2HP Automatic Beneath Shelf Mount

Retail Price \$815

Model USM-200 (not shown)
2HP Manual Beneath Shelf Mount

Retail Price \$735

Standard Features

	USM-5300	USM-250	USM-200
Automatic selection for 2 handpieces (1 wet, 1 dry)	✓	✓	
Manual selection for 2 handpieces (1 wet, 1 dry)			✓
3-way syringe with quick change autoclavable tip	✓	✓	✓
1 autoclavable HVE and saliva ejector valve with tubing	✓		
Vacuum canister with removable screen	✓		
Master On/Off toggle		✓	
Pressure adjustment for each handpiece	✓	✓	✓
Easy to read pressure gauge	✓	✓	✓
Air coolant flow adjustment	✓		
Non-retracting water coolant adjustment	✓	✓	✓
Asepsis handpiece tubing, sterling (specify straight or coiled)	✓	✓	✓
Telescoping arm with breakaway feature, 29" to 40" horizontal reach	✓		
Standard foot control	✓	✓	✓
7' umbilical tubing	✓	✓	✓

Optional Features

	Retail Price	USM-5300	USM-250	USM-200
Utility center assembly, 14" x 9" x 5 ³ / ₈ "	add \$400	✓	✓	✓
Water coolant to both handpieces	add \$40	✓	✓	✓
Silcryn™ Tubing	add \$15	✓	✓	✓
Additional umbilical tubing, per foot	add \$11	✓	✓	✓
Unit mount post, 28" with bushing	add \$125	✓		
Unit/light post, 40" with bushing	add \$165	✓		
Stainless steel instrument tray, 9 ³ / ₄ " x 13 ¹ / ₂ " 3	add \$60	✓		

Cabinet Mount Delivery Units

Doctor's Rear Delivery Units

Cabinet shown for display only, not included in price.



Model USM-3305
3HP Automatic Rear Delivery
Slide-Out Mount Control Head

Retail Price: \$2,090

Model USM-3350
3HP Automatic Rear Delivery
Slide-Out Mount Control Head with Vacuum

Retail Price: \$2,250

Standard Features	USM 3305	USM 3350
Asepsis designed automatic selection for 3 handpieces	✓	✓
Side-mounted recessed control adjustments	✓	✓
3-way syringe with quick change autoclavable tip	✓	✓
Master On/Off toggle with indicator	✓	✓
Pressure adjustment for each handpiece	✓	✓
Easy to read pressure gauge (under cover)	✓	✓
Air coolant flow adjustment	✓	✓
Non-retracting water coolant adjustment & purge	✓	✓
Asepsis handpiece tubing, sterling	✓	✓
Vacuum canister with removable 1 ⁷ / ₈ " screen		✓
1 autoclavable HVE and Saliva Ejector valve with tubing		✓
Side-out mount with 15° pivot left and right	✓	✓
Standard foot control	✓	✓
Master valve assembly with filter and regulator	✓	✓
7' umbilical tubing	✓	✓

Optional Features	Retail Price	USM 3305	USM 3350
Delete master valve assembly with filter and regulator	deduct \$280	✓	✓
Delete handpiece station with tubing	deduct \$50	✓	✓
Additional handpiece station with tubing	add \$150	✓	✓
Additional syringe on doctor or assistant's side	add \$120	✓	✓
Additional umbilical tubing, per foot	add \$11	✓	✓
Wet/dry foot control in lieu of standard	add \$100	✓	✓
Chip blower foot control in lieu of standard	add \$100	✓	✓
Wet/dry chip blower foot control in lieu of standard	add \$150	✓	✓
Silcryn™ Tubing	add \$30	✓	✓

Cabinet Mount Delivery Units

Doctor's Rear Delivery Units

Cabinet shown for display only, not included in price.



Model USM-4350
3HP Automatic Rear Delivery
Duo Swing System with Vacuum

Retail Price: \$2,600

Model USM-5100
Rear Delivery Assistant's Swivel
Vacuum System

Retail Price: \$790

Standard Features

	USM 4350	USM 5100
Pivoting arm mounted asepsis designed automatic selection for 3 handpieces	✓	
Side-mounted recessed control adjustments	✓	
3-way syringe with quick change autoclavable tip	✓	✓
Master On/Off toggle with indicator	✓	
Pressure adjustment for each handpiece	✓	
Easy to read pressure gauge (under cover)	✓	
Air coolant flow adjustment	✓	
Non-retracting water coolant adjustment & purge	✓	
Asepsis handpiece tubing, sterling	✓	
Pivoting arm mounted vacuum canister with removable 1 ⁷ / ₈ " screen	✓	✓
1 autoclavable HVE and Saliva Ejector valve with tubing	✓	✓
Standard foot control	✓	
Master valve assembly with filter and regulator	✓	
7' umbilical tubing	✓	✓

Optional Features

	Retail Price	USM 4350	USM 5100
Delete master valve assembly with filter and regulator	deduct \$280	✓	
Delete handpiece station with tubing	deduct \$50	✓	
Additional handpiece station with tubing	add \$150	✓	
Additional syringe on doctor or assistant's side	add \$120	✓	
Additional umbilical tubing	add \$11	✓	✓
Wet/dry foot control in lieu of standard	add \$100	✓	
Chip blower foot control in lieu of standard	add \$100	✓	
Wet/dry chip blower foot control in lieu of standard	add \$150	✓	
Silcryn™ Tubing	add \$30	✓	
Silcryn™ Tubing	add \$15		✓

Cabinet Mount Delivery Units

Side Delivery Unit



Cabinet shown for display only,
not included in price.

Model USM-3360

Side Delivery with Pneumatic Brake Flex Arm

Retail Price: \$3,400

Standard Features

	USM 3310
Asepsis designed automatic selection for 3 handpieces	✓
Side-mounted recessed control adjustments	✓
3-way syringe with quick change autoclavable tip	✓
Master On/Off toggle with indicator	✓
Pressure adjustment for each handpiece	✓
Easy to read pressure gauge (under cover)	✓
Air coolant flow adjustment	✓
Non-retracting water coolant adjustment & purge	✓
Asepsis handpiece tubing, sterling (specify straight or coiled)	✓
Asepsis flex arm system with pneumatic brake, 9" rigid arm, 39" reach, 19" vertical travel	✓
2" wall mount bracket	✓
Side-out mount with 15° pivot left and right	✓
Standard foot control	✓
Master valve assembly with filter and regulator	✓
7' umbilical tubing	✓

Optional Features

	Retail Price	
Delete master valve assembly with filter and regulator	deduct \$280	✓
Delete handpiece station with tubing	deduct \$50	✓
Additional handpiece station with tubing (dry only)	add \$150	✓
2" dia. post mount rigid arm, 24" (in lieu of 9" arm)	N/C	✓
Additional umbilical tubing, per foot	add \$11	✓
Wet/dry foot control in lieu of standard	add \$100	✓
Chip blower foot control in lieu of standard	add \$100	✓
Wet/dry chip blower foot control in lieu of standard	add \$150	✓
Silcryn™ Tubing	add \$20	✓

**See additional selection of Arm Systems and Brackets in the "Arms & Trays" section.

Cabinet Mount Rear Delivery Units

Filter/Regulators must be ordered separately



Model USM-4700
Rear Delivery, 3HP Automatic
Under-Counter Mount

Retail Price: \$1,435

Model USM-4710
Rear Delivery, 3HP Automatic
Under-Counter Mount, With Vacuum

Retail Price: \$1,630

Standard Features

	USM 4700	USM 4701
Automatic selection for 3 handpieces	✓	✓
Slide out access for service	✓	✓
3-way syringe with quick change autoclavable tip	✓	✓
Master On/Off toggle with indicator	✓	✓
Pressure adjustment for each handpiece	✓	✓
Easy to read pressure gauge	✓	✓
Non-retracting water coolant adjustment and purge	✓	✓
Standard foot control	✓	✓
Asepsis handpiece tubing, sterling (specify straight or coiled)	✓	✓
1 each HVE & SE with tubing		✓
Vacuum canister with removable screen		✓
7' umbilical tubing	✓	✓

Optional Features

	Retail Price	USM 4700	USM 4701
Utility center assembly (less housing & cover)	add \$280	✓	✓
Delete handpiece station with tubing	deduct \$50	✓	✓
Additional umbilical tubing	add \$11	✓	✓
Wet/dry foot control in lieu of standard	add \$100	✓	✓
Chip blower foot control in lieu of standard	add \$100	✓	✓
Wet/dry chip blower foot control in lieu of standard	add \$150	✓	✓
Silcryn™ Tubing	add \$30	✓	✓

Doctor's Mobile Delivery Units

Introducing
Mobile Carts from
Beaverstate Dental
Equipment!



Model A-3100 3HP Automatic Swing Mobile Unit	Model A-3150 3HP Automatic Duo Swing Mobile Unit w/Vacuum	Model A-3101 3HP Automatic Cart	Model A-3151 3HP Automatic Cart with Vacuum
Retail Price \$2,820	Retail Price \$3,400	Retail Price \$2,670	Retail Price \$2,790

Standard Features	A-3100	A-3150	A-3101	A-3151
Asepsis designed automatic selection for 3 handpieces	✓	✓	✓	✓
3-way syringe with quick change autoclavable tip	✓	✓	✓	✓
Dual 3-way syringe with quick change autoclavable tip		✓		✓
Master On/Off toggle with indicator	✓	✓	✓	✓
Pressure adjustment for each handpiece	✓	✓	✓	✓
Easy to read pressure gauge (protected under cover)	✓	✓	✓	✓
Non-retracting water coolant adjustment, air coolant adjustment and purge	✓	✓	✓	✓
Asepsis handpiece tubing, sterling (specify straight or coiled)	✓	✓	✓	✓
16" x 22" Starboard™ Work Surface, Dolphin Gray			✓	✓
Painted H-frame with height adjustment, 27" to 35"	✓	✓	✓	✓
Vacuum canister with removable screen		✓	✓	✓
1 autoclavable HVE and saliva ejector valve, tubing and holder		✓	✓	✓
Standard foot control	✓	✓	✓	✓
Utility center assembly, 14" x 9" x 5 ³ / ₈ ", with valves and regulators	✓	✓	✓	✓
7' umbilical tubing	✓	✓	✓	✓

Optional Features	Retail Price	A-3100	A-3150	A-3101	A-3151
Delete utility center	deduct \$400	✓	✓	✓	✓
Delete handpiece station with tubing	deduct \$50	✓	✓	✓	✓
Additional handpiece station with tubing (dry only)	add \$150	✓	✓	✓	✓
Additional syringe for doctor or assistant's side	add \$120	✓	✓	✓	✓
Add QD Kit (1/4" and 3/8" QD's with bracket & hardware)	add \$110	✓	✓	✓	✓
Quick disconnect kit (1/4" with bracket and hardware)	add \$38	✓	✓	✓	✓
Delete storage organizer assembly	deduct \$300			✓	✓
Additional umbilical tubing, per foot	add \$11	✓	✓	✓	✓
Wet/dry foot control in lieu of standard	add \$100	✓	✓	✓	✓
Chip blower foot control in lieu of standard	add \$100	✓	✓	✓	✓
Wet/dry chip blower foot control in lieu of standard	add \$150	✓	✓	✓	✓
Silcryn™ Tubing	add	\$20	\$30	\$20	\$30

Doctor's Mobile Delivery Units



Model USM-4100

3HP Automatic Doctor's Mobile Unit

Retail Price \$2,365

Model USM-4150

3HP Automatic Duo Mobile Unit with Vacuum

Retail Price \$2,630

Model USM-4250

3HP Automatic Duo Swing Mobile Unit with Vacuum

Retail Price \$4,000

Standard Features

	USM-4100	USM-4150	USM-4250
Asepsis designed automatic selection for 3 handpieces	✓	✓	✓
Side mounted control adjustments	✓	✓	✓
Pivoting arms, 360° swing for easy left/right conversion			✓
3-way syringe with quick change autoclavable tip	✓	✓	✓
Master On/Off toggle with indicator	✓	✓	✓
Pressure adjustment for each handpiece	✓	✓	✓
Easy to read pressure gauge (protected under cover)	✓	✓	✓
Non-retracting water coolant adjustment, air coolant adjustment and purge	✓	✓	✓
Asepsis handpiece tubing, sterling (specify straight or coiled)	✓	✓	✓
16" x 22" Starboard™ Work Surface, Dolphin Gray			✓
Painted H-frame with height adjustment, 27" to 35"	✓	✓	✓
Vacuum canister with removable screen		✓	✓
1 autoclavable HVE and saliva ejector valve, tubing and holder		✓	✓
Standard foot control	✓	✓	✓
Utility center assembly, 14" x 9" x 5 ³ / ₈ ", with valves and regulators	✓	✓	✓
7' umbilical tubing	✓	✓	✓

Optional Features

	Retail Price	USM-4100	USM-4150	USM-4250
Delete utility center	deduct \$400	✓	✓	✓
Delete handpiece station with tubing	deduct \$50	✓	✓	✓
Additional handpiece station with tubing (dry only)	add \$150	✓	✓	✓
Additional syringe for doctor or assistant's side	add \$120	✓	✓	✓
Add QD Kit (1/4" and 3/8" QD's with bracket & hardware)	add \$110	✓	✓	✓
Quick disconnect kit (1/4" with bracket and hardware)	add \$38	✓	✓	✓
Delete storage organizer assembly	deduct \$300			✓
Additional umbilical tubing, per foot	add \$11	✓	✓	✓
Wet/dry foot control in lieu of standard	add \$100	✓	✓	✓
Chip blower foot control in lieu of standard	add \$100	✓	✓	✓
Wet/dry chip blower foot control in lieu of standard	add \$150	✓	✓	✓
Silcryn™ Tubing	add	\$20	\$30	\$30

Doctor's Mobile Delivery Units

Assistant's Mobile Delivery Units



USM-4100A 3HP Automatic Doctor's Mobile Unit	USM-4150A 3HP Automatic Doctor's Mobile Unit with Vacuum	USM-4500A Assistant's Mobile Unit	USM-4550A Assistant's Mobile Unit with Vacuum
Retail Price \$2,730	Retail Price \$3,265	Retail Price \$1,260	Retail Price \$1,610

		USM-4100A	USM-4150A	USM-4500A	USM-4550A
Standard Features					
Automatic selection for 3 handpieces		✓	✓		
3-way syringe with quick change autoclavable tip		✓	✓		✓
Master On/Off toggle with indicator		✓	✓		
Pressure adjustment for each handpiece		✓	✓		
Easy to read pressure gauge		✓	✓		
Non-retracting water coolant adjustment, air coolant adjustment & purge		✓	✓		
Asepsis handpiece tubing, sterling (specify straight or coiled)		✓	✓		
16" x 22" Starboard™ Work Surface, Dolphin Gray		✓	✓	✓	✓
Storage organizer assembly		✓	✓	✓	✓
Painted H-frame with height adjustment, 27" to 35"		✓	✓	✓	✓
Vacuum canister with removable screen			✓		✓
1 autoclavable HVE and saliva ejector valve, tubing and holder			✓		✓
Standard foot control		✓	✓		
Utility center assembly, 14" x 9" x 5 ³ / ₈ " with valves and regulators			✓	✓	
7' umbilical tubing		✓	✓		✓
Optional Features					
Delete utility center	deduct \$400	✓	✓		
Delete handpiece station with tubing	deduct \$50	✓	✓		
Additional handpiece station with tubing (dry only)	add \$150	✓	✓		
Additional syringe for assistant's side	add \$120	✓	✓		
QD Kit (1/4" & 3/8" QD's w/bracket & hardware)	add \$110	✓	✓		✓
Additional QD Kit (1/4" with mounting bracket)	add \$38	✓	✓		
Delete storage organizer	deduct \$300		✓	✓	✓
Additional umbilical tubing, per foot	add \$11	✓	✓		✓
Wet/dry foot control in lieu of standard	add \$100	✓	✓		
Chip blower foot control in lieu of standard	add \$100	✓	✓		
Wet/dry chip blower foot control in lieu of standard	add \$150	✓	✓		
Silcryn™ Tubing	add	\$20	\$30		\$15

Handpiece Control Heads

Filter/Regulators must be ordered separately



Model USM-330

3HP Automatic Control Head Only

Retail Price \$1,530

Model USM-300

3HP Automatic Control Head Only

Retail Price \$1,465

Standard Features

	USM-330	USM-300
Automatic selection for 3 handpieces	✓	✓
Side-mounted recessed control adjustments	✓	✓
3-way syringe with quick change autoclavable tip	✓	✓
Master On/Off toggle with indicator	✓	✓
Pressure adjustment for each handpiece	✓	✓
Easy to read pressure gauge	✓	✓
Non-retracting water coolant adjustment & purge	✓	✓
Air coolant flow adjustment	✓	✓
Water coolant On/Off and flow control	✓	✓
Asepsis handpiece tubing, sterling (specify straight or coiled)	✓	✓
Standard foot control	✓	✓
7' umbilical tubing	✓	✓

Optional Features

	Retail Price	USM-330	USM-300
Utility center assembly, 14" x 9" x 5 ³ / ₈ ", with valves and regulators	add \$400	✓	✓
Delete handpiece station with tubing	deduct \$50	✓	✓
Additional handpiece station with tubing (dry only)	add \$150	✓	✓
Add stainless steel tray, 10" x 6 ¹ / ₂ " and pad	add \$50	✓	✓
Add stainless steel tray, 9 ³ / ₄ " x 13 ¹ / ₂ " and pad	add \$60	✓	✓
Add plastic tray, 9 ³ / ₄ " x 13 ¹ / ₂ "	add \$30	✓	✓
Additional umbilical tubing, per foot	add \$11	✓	✓
Wet/dry foot control in lieu of standard	add \$100	✓	✓
Chip blower foot control in lieu of standard	add \$100	✓	✓
Wet/dry chip blower foot control in lieu of standard	add \$150	✓	✓
Silcryn™ Tubing	add \$20	✓	✓

Panel Mount Handpiece Controls

Filter/Regulators must be ordered separately



Model USM-200P
2HP Manual Control

Retail Price \$750

Model USM-220P
2HP Automatic Control

Retail Price \$860

Model USM-400P
3HP Automatic Control

Retail Price \$1,155

Standard Features

	USM-200P	USM-220P	USM-400P
Automatic selection for 3 handpieces			✓
Automatic selection for 2 handpieces		✓	
Manual selection for 2 handpieces	✓		
3-way syringe with quick change autoclavable tip	✓	✓	✓
Instrument holder bar, 6"	✓	✓	
Instrument holder bar, 8"			✓
Master On/Off toggle with indicator			✓
Pressure adjustment for each handpiece	✓	✓	✓
Easy to read pressure gauge	✓	✓	✓
Non-retracting water coolant, adjustment and purge			✓
Water coolant adjustment, for 3 handpieces			✓
Water coolant adjustment, 1 handpiece	✓	✓	
Asepsis handpiece tubing, sterling (specify straight or coiled)	✓	✓	✓
Standard foot control	✓	✓	✓
7' umbilical tubing		✓	✓
7' of 1/4" supply line	✓		

Optional Features

	Retail Price	USM-200P	USM-220P	USM-400P
Master valve assembly with filters and regulators (less housing & cover)	add \$280	✓	✓	✓
Water coolant to both handpieces	add \$50	✓	✓	
Additional supply line, per foot	add \$2	✓		
Additional umbilical tubing, per foot	add \$11		✓	✓
Wet/dry foot control in lieu of standard	add \$100			✓
Chip blower foot control in lieu of standard	add \$100			✓
Wet/dry chip blower foot control in lieu of standard	add \$150			✓
Silcryn™ Tubing	add	\$15	\$15	\$20

Panel Mount Handpiece Controls

Filter/Regulators must be ordered separately



Model-4500P

3HP Automatic Control

Retail Price \$1,550

Model USM-4550P

3HP Automatic Control with Vacuum

Retail Price \$1,730

Standard Features

	USM-4500P	USM-4550P
Automatic selection for 3 handpieces	✓	✓
Slide out service access	✓	✓
3-way syringe with quick change autoclavable tip	✓	✓
Instrument holder bar, 12"	✓	
Instrument holder bar, 16"		✓
Master On/Off toggle with indicator	✓	✓
Pressure adjustment for each handpiece	✓	✓
Easy to read pressure gauge	✓	✓
Non-retracting water coolant, adjustment and purge	✓	✓
HVE and saliva ejector valves, tubing and holder		✓
Vacuum canister with removable screen		✓
Asepsis handpiece tubing, sterling (specify straight or coiled)	✓	✓
Standard foot control	✓	✓
7' umbilical tubing	✓	✓

Optional Features

	Retail Price	USM-4500P	USM-4550P
Master valve assembly with filters and regulators (less housing & cover)	add \$280	✓	✓
Additional umbilical tubing, per foot	add \$11	✓	✓
Wet/dry foot control in lieu of standard	add \$100	✓	✓
Chip blower foot control in lieu of standard	add \$100	✓	✓
Wet/dry chip blower foot control in lieu of standard	add \$150	✓	✓
Silcryn™ Tubing	add	\$20	\$30

Vacuum Canister Horizontal Surface Mounting Bracket

VCF-339-00 - \$25



Panel Mount Handpiece Control Heads

Filter/Regulators must be ordered separately



Model-4501P
3HP Automatic Control

Retail Price \$1,520

Model USM-4551P
3HP Automatic Control with Vacuum

Retail Price \$1,720

Standard Features

	USM-4501P	USM-4551P
Automatic selection for 3 handpieces	✓	✓
Slide out service access	✓	✓
3-way syringe with quick change autoclavable tip	✓	✓
Dual 3-way syringe with quick change autoclavable tip		✓
Dr.'s side has 10-1/2" cabinet mounted pivoting bar	✓	✓
Assistant's side has 8-1/2" cabinet mounted pivoting bar		✓
Master On/Off toggle with indicator	✓	✓
Pressure adjustment for each handpiece	✓	✓
Easy to read pressure gauge	✓	✓
HVE and saliva ejector valves, tubing and holder		✓
Vacuum canister with removable screen		✓
Non-retracting water coolant, adjustment and purge	✓	✓
Asepsis handpiece tubing, sterling (specify straight or coiled)	✓	✓
Standard foot control	✓	✓
7' umbilical tubing	✓	✓

Optional Features

	Retail Price	USM-4501P	USM-4551P
Master valve assembly with filters and regulators (less housing & cover)	add \$280	✓	✓
Delete handpiece station with tubing	deduct \$50	✓	✓
Additional umbilical tubing, per foot	add \$11	✓	✓
Wet/dry foot control in lieu of standard	add \$100	✓	✓
Chip blower foot control in lieu of standard	add \$100	✓	✓
Wet/dry chip blower foot control in lieu of standard	add \$150	✓	✓
Silcryn™ Tubing	add	\$20	\$35



**Vacuum Canister Horizontal
Surface Mounting Bracket**
VCF-339-00 - \$25

Wall Mount Handpiece Control Heads

Filter/Regulators must be ordered separately



Model USM-220A
2HP Automatic Control Head

Retail Price \$690



Model USM-222A (not shown)
2HP Automatic Control Head, Both HP's Wet

Retail Price \$740

Model USM-400A
3HP Automatic Control Head

Retail Price \$1,020



Model USM-200A
2HP Manual Control Head

Retail Price \$610



Model USM-202A (not shown)
2HP Manual Control Head, Both HP's Wet

Retail Price \$650

Model USM-100A
1HP Control

Retail Price \$510



Model USM-101A
Lab HP Control

Retail Price \$320



Model USM-050A
Lab HP Control Kit

Retail Price \$340

	USM-220A	USM-400A	USM-200A	USM-100A	USM-101A	USM-050A
Standard Features						
Automatic selection for 2 handpieces	✓					
Automatic selection for 3 handpieces		✓				
Manual selection for 2 handpieces			✓			
Single handpiece control with coolant				✓		
Single handpiece control without coolant					✓	✓
3-way syringe with quick change autoclavable tip	✓	✓	✓	✓		
Master On/Off toggle		✓				
Pressure adjustment for each handpiece	✓	✓	✓	✓	✓	✓
Easy to read pressure gauge	✓	✓	✓	✓	✓	✓
Non-retracting water coolant adjustment for 1 HP (1 wet/1 dry)	✓		✓			
Water coolant On/Off and flow control				✓		
Asepsis handpiece tubing, sterling (specify straight or coiled)	✓	✓	✓	✓	✓	✓
Standard foot control	✓	✓	✓	✓	✓	✓
7' umbilical tubing		✓		✓		
7' of 1/4" supply line	✓		✓		✓	✓
Optional Features						
Utility center assembly, 14" x 9" x 5 3/8" with filters and regulators	add \$400	✓	✓	✓	✓	✓
Additional air supply tubing, per foot	add \$2	✓		✓	✓	✓
Additional umbilical tubing, per foot	add \$11		✓	✓		
Silcryn™ Tubing	add	\$15	\$20	\$15	\$10	\$5

Pivot Mount (1/2" dia.)
ASF-104-02 - \$35



Horizontal Surface Mount Bracket
ASF-106-01 - \$20



Utility Center and Master Valves
VAF-144-017 - \$400



Vacuum Arm Systems



Model USM-5200A

Telescoping Arm System with Vacuum Package, 1 HVE, 1 SE, 1 Syringe

Retail Price \$950

Model USM-5220A

Telescoping Arm System with Vacuum Package, 1 HVE, 1 SE, No Syringe

Retail Price \$800

Model USM-5222A

Telescoping Arm System with Vacuum Package, 2 HVEs, 1 SE, No Syringe

Retail Price \$870

		USM-5200A	USM-5220A	USM-5222A
Standard Features				
Telescoping arm with breakaway feature, 31" to 42" horizontal reach		✓	✓	✓
Holder bar with 3 holders		✓		✓
3-way syringe with quick change autoclavable tip		✓		
Flow adjustment for 3-way syringe		✓		
Vacuum canister with removable screen		✓	✓	✓
1 autoclavable HVE and saliva ejector valve with tubing		✓	✓	
2 autoclavable HVEs and saliva ejector valve with tubing				✓
Standard foot control		✓	✓	✓
2" dia. post mount (post not included, if needed specify length)		✓	✓	✓
10' umbilical tubing		✓	✓	✓
Optional Features				
Additional umbilical, per foot	Retail Price			
	add \$10	✓	✓	✓
Silcryn tubing	add \$15	✓	✓	✓
Tray mount assembly with SS tray (see image below)	add \$175	✓	✓	✓
Wall mount bracket in lieu of post mount	N/C	✓	✓	✓



Tray mount assembly with SS tray

Vacuum Arm Systems



Model USM-5035A

Horizontal Folding Arm, 32",
2" dia. post mount with Vacuum Package

Retail Price \$765

Model USM-5125A

Horizontal Folding Arm, 24",
2" dia. wall mount with Vacuum Package

Retail Price \$720

Model USM-5000A [shown]

Horizontal Folding Arm, 24",
2" dia. post mount with Vacuum Package

Retail Price \$740

Model USM-5130A

Horizontal Folding Arm, 32",
2" dia. wall mount with Vacuum Package

Retail Price \$750

		USM-5000A	USM-5035A	USM-5125A	USM-5130A
Standard Features					
Folding arm, 24"		✓		✓	
Folding arm, 32"					✓
Holder bar with 3 holders		✓	✓	✓	✓
3-way syringe with quick change autoclavable tip		✓	✓	✓	✓
Flow adjustment for 3-way syringe		✓	✓	✓	✓
Vacuum canister with removable screen		✓	✓	✓	✓
1 autoclavable HVE and saliva ejector valve with tubing		✓	✓	✓	✓
Standard foot control 2" dia. post mount (post not included, if needed specify length)		✓	✓	✓	✓
Wall mount bracket				✓	✓
5' umbilical tubing		✓	✓	✓	✓
Optional Features	Retail Price				
1/4" air or water quick disconnect	add \$38	✓	✓	✓	✓
With flow control	add \$62				
Additional HVE valve with tubing holder & 4-position bar holder	add \$110	✓	✓	✓	✓
Additional umbilical, per foot	add \$8	✓	✓	✓	✓
Silcryn tubing	add \$15	✓	✓	✓	✓

Vacuum Arm Systems



Model USM-5230A

Hygiene System, 2" dia. post mount with Vacuum Telescoping Arm

Retail Price \$1,015

Model USM-5525A

Hygiene System, 1HP, 2" dia. post mount with Folding Arm and Vacuum Package

Retail Price \$1,000

Standard Features	USM-5230A	USM-5525A
1 (dry only) handpiece		✓
Folding arm, 24"		✓
Tray mount with holder bar and 3 holders		✓
3-way syringe with quick change autoclavable tip	✓	✓
Flow adjustment for 3-way syringe	✓	✓
1/4" Female water quick disconnect	✓	
Asepsis handpiece tubing, sterling		✓
Plastic instrument tray, 9-3/4" x 13-1/2"		✓
Vacuum canister with removable screen	✓	✓
1 saliva ejector valve with tubing	✓	✓
Standard foot control		✓
2" dia. post mount (post not included, if needed specify length)	✓	✓
5' umbilical tubing	✓	
7' umbilical tubing		✓
Optional Features	Retail Price	
Water quick disconnect	add \$62	✓
With flow control	add \$125	
Unit mount post, 28" with bushing	add \$165	✓
Unit/light post, 40" with bushing	add \$9	✓
Jacketed tubing, per foot	add \$15	✓
Silcryn™ Tubing	add \$15	✓

Cuspidors



Model USM-3200

Asepsis Design Cuspidor

Retail Price \$1,015

Model USM-3250

Asepsis Design Cuspidor with Vacuum Accessories and Water Quick Disconnect

Retail Price \$1,530

Standard Features

	USM-3200	USM-3250
Easy clean asepsis design	✓	✓
Ceramic bowl	✓	✓
Quick release autoclavable cup filler and bowl rinse spout	✓	✓
Adjustable timed bowl rinse	✓	✓
Gravity drain with dual venting	✓	✓
Vacuum canister with removable screen		✓
1 autoclavable HVE and saliva ejector valve with tubing		✓
1/4" Female water quick disconnect		✓
Low voltage electric vacuum switch		✓
2" dia. Post mount	✓	✓
6' umbilical tubing	✓	✓

Optional Features

	Retail Price		
Timed cup filler	add \$190	✓	✓
Additional umbilical tubing	add \$11	✓	✓
3-Way Syringe with tubing and holder	add \$130	✓	✓
Silcryn™ Tubing	add \$10	✓	✓

Celux Operatory Lights

- Exclusive ball and socket pivot allows unrestricted multi axis positioning
- Lamp head may be adjusted vertically, horizontally or diagonally
- Sharply defined rectangular light pattern, 3" x 7" at 28", 3½" x 4½" at 29"
- 3800° to 4500° Kelvin temperature
- 4500° to 5100° Kelvin temperature
- Easy reach ON/OFF/intensity switch
- 2-way switch for low (12,000 LUX) and high (24,000 LUX) intensity
- 2-way switch for low (18,000 LUX) and high (23,000 LUX) intensity
- Glass reflector with dichroic coating filters out infrared rays, providing cool light
- 12V—55 Watt tungsten halogen H-3 bulb with twin filaments
- 24 V, 100 Watt halogen bulb with single filament (designed for use in Aerolight with 1000 hour rating)
- Easy bulb replacement
- Easy open safety shield allows for cleaning
- Slim profile minimizes heat
- Available in 110 Volt or 220 Volt
- Available in 110 Volt
- Made in Taiwan



Catalog No.	Description	Retail Price
CLX-100	Celux Dental Operatory Light	\$1,250

LED Operatory Light

- Available in 4 configurations: Post Mount, Ceiling Mount, Wall Mount, Cabinet Mount
- Sharply defined rectangular light pattern, 4" x 6" at 29"
- Easy reach feather touch on/off switch
- 3 axis head rotation
- "T" handles for easy positioning
- 5000 low setting to 5500 high setting, Kelvin temperature
- Low voltage 12 VDC power supply
- Low energy consumption 4-12 watts
- Dimming feature from low to high/high to low
- Accommodates ceiling heights from 8ft to 10ft

Catalog No.	Description	Retail Price
LED-105	Post Mount	\$2,200
LED-105CL	Ceiling Mount	\$2,790
LED-105W	Wall Mount	\$3,230
LED-105CA	Cabinet Mount	\$3,230



Notes:

- When ordering LED light in Ceiling Mount configuration, provide ceiling height measured from floor to hard mounting surface.
- For total mounting bracket concealment in drop-ceiling configuration, allow 7-1/4" inch depth between hard surface mount to drop ceiling.
- Ceiling mount light is shipped with domed cover to completely conceal mounting bracket.



Autoscaler® Ultrasonic Dental Scaler



- Autoscaler® Ultrasonic Dental Scaler with Automatic Tuning
- Made in USA.
- Special Periodontal and Routine Scaling
- Precision Water Control
- Warranty: 6 years on Circuit Boards, 2 years on all other components.

Catalog No.	Description	Retail Price
US-SEI-01	Table Top Model, AW 25C	\$665.00
US-SEI-02	Table Top Model, AW 30C	\$665.00
US-INS-01	Universal Curette, each	\$96.00
US-INS-02	Slimline Style Perio, each	\$96.00

Specifications for table top models AW25C & AW30C

Height: 2.2" Weight: 8.6" Depth: 7.9"
 Weight: 7 pounds
 Insert Stroke: 0.001/0.003 elliptical motion

Output Power

Operator adjustable from front panel up to 70 watts. Watts measured at AC primary. Linear Power Control.

Water Input Pressure Needed

Minimum 20 PSIG to MAX 80 PSIG pressure range input city or bottle water supply. Water flow requirement from insert is 18cc per minute MINIMUM FLOW to 40cc per minute.

AC Current Draw and Input Voltage

0.8 amps @ 110/125 VAC 50/60Hz fuse with AGC 2.0 amp fast blow style

Special Order Model

0.4 amps @ 220/230 VAC 50/60 Hz. Replace fuse with AGC 2.0 amp fast blow style

*** Built-in model Autoscaler also available

Autoscaler Is Supplied With (approx. measurement)

- 1 - 8.0' Detachable handpiece assembly
- 1 - 8.0' Foot Switch
 (not supplied when ordered with Air Micro Switch)
- 1 - 8.0' Water Tubing with 1/4 Male NON Shut Off quick release attached
- 1 - Complete instructions and internet web support

Frequency Generated

2 models offered: 25,000Hz CPS and 30,000Hz CPS (cycles per second)

Accepts DENTSPLY® CAVITRON® Brand (TFI, P, SLI and FSI inserts)

Optional Feature

Air Micro Switch – allows the Autoscaler® to be interfaced into the operatory air system, bypassing the need for a manual footswitch

Required air pressure range to function properly is 20 PSIG to 60 PSIG. Manual foot switch not included with this option.

Solmetex® Hg5 Series Amalgam Separators



Provide outstanding performance capturing amalgam and sediment from dental waste streams.

- | | |
|---|--|
| <p>Includes:</p> <ul style="list-style-type: none"> ISO 11143 certified at greater than 99% Easy to operate & virtually maintenance-free Affordable & cost-efficient Flexible design for hard-to-fit applications Functional for wet or dry vacuum systems Hg5 Series is clear Supported by responsive experts Environmentally friendly | <ul style="list-style-type: none"> NO Contracts with hidden fees NO Electricity NO Daily maintenance NO Decanting NO Pumps NO Timers |
|---|--|

Catalog No.	Description	Retail Price
Hg5-001	<p>Hg5 amalgam separator</p> <ul style="list-style-type: none"> - Accommodates up to 10 operatories - Dimensions 29½”h x 10”w x 8”d - 1½” inlet and outlet piping - Fernco couplers for easy install - Floor/wall mount brackets provided 	\$860.00
Hg5-Mini	<p>Mini amalgam separator</p> <ul style="list-style-type: none"> - Accommodates up to 4 operatories - Dimensions 21½”h x 11”w x 8”d - ¾” inlet and outlet piping - Fernco couplers for easy install - Floor/wall mount brackets provided - Designed specifically to fit under a counter 	\$860.00
Hg5HV-001	<p>High volume amalgam separator</p> <ul style="list-style-type: none"> - Accommodates up to 20 operatories - Dimensions 18” w x 28”h x 22½”d - 3” inlet and outlet piping - Wall mount brackets provided 	\$2,860.00
Hg5-002CR	Collection container with recycle kit	\$338.00
Hg5-002C	Collection container only	\$197.00

Vacuum Pumps And Air Compressors

We offer numerous standard and customized, quality models of vacuum pumps and air compressors to meet your needs. Please contact our office with your requirements for a competitive quotation.

EcoVac Dry
Vacuum System



Liquid Ring
Vacuum System



Dental Air Dryers



Tech-West Ultra-Clean Oil-less & Oil-free Air Compressors



DENTAL EQUIPMENT	7 - 37
DENTAL COMPONENTS	41 - 125
DENTAL TUBING	129 - 139
NUMERICAL INDEX PAGES	141



Arms / Trays	41
Syringes	51
Foot Controls	55
Vacuum Components	59
Holders / Holder Bars	69
Valves / Filters / Regulators	77
Unit Accessories	91
Lighting / Reflectors / Bulbs	93
Fittings	99
Service Parts	115
Handpieces / HP Accessories / HP Connectors	119

We offer a wide selection of American made, quality dental components that are currently installed in numerous well-known domestic and international brands. Compared to other suppliers, you will find our parts excel in workmanship and compete in price. Ask us about customizing parts for your application.

2" Dia. Ring Mount Flexible Arm System

- Asepsis 2" dia. Post Mount Flex Arm System with Pneumatic Brake
- Control Mount Bracket
- 24" Rigid Arm and 24" Flex Arm
- Adjustable for weights from 3 to 15 lbs.

Catalog No.	Description	Retail Price
ASF-100-01	24" Rigid w/24" Flexible Arm	\$940.00



2" Dia. Ring Mount Flexible Arm System

- Asepsis 2" dia. Ring Mount Flex Arm System with Pneumatic Brake
- Control Mount Bracket
- 9" Rigid Arm and 24" Flex Arm
- Requires 2" Wall Mount (item ASF-103-02)
- Short rigid arm ideal for side delivery. Adjustable for weights from 3 to 15 lbs.

Catalog No.	Description	Retail Price
ASF-100-02	9" Rigid w/24" Flexible Arm	\$940.00
ASF-103-02	Wall Mount for 2" dia Post Mounted Arms	\$135.00



Top Post Mount Flexible Arm System

- Asepsis Top Post Mount Flex Arm System with Pneumatic Brake
- Control Mount Bracket
- 24" Rigid Arm and 24" Flex Arm
- Adjustable for weights from 3 to 15 lbs.

Catalog No.	Description	Retail Price
ASF-100-03	24" Rigid w/24" Flexible Arm with Arm bushing (ACF-594-01)	\$940.00



Control Mount Bracket

Catalog No.	Description	Retail Price
ASF-106-350	9" Bracket only	\$160.00



Manual 2" Dia. Post Mount Flexible Arm System

- Manual 2" Post Mount Flex Arm System with Friction Brake and Tray.
- Measures 15" deep x 13½" wide
- Holds 9¾" x 13½" tray
- 63" reach, 27" vertical travel with light spring

Catalog No.	Description	Retail Price
ASF-100-04	28" Rigid w/24" Flexible Arm	\$950.00
ASF-100-04TP	28" Rigid w/24" Flexible Arm, top post mount	\$950.00



Horizontal Folding Arm Systems

- Horizontal Folding Arms
- 3 position holder bar with 3 handpiece holders
- 2" dia. post mount or wall mount
- 24" or 32" horizontal reach

Catalog No.	Description	Retail Price
ASF-100-05	24" Wall Mount	\$240.00
ASF-100-06	24" Post Mount	\$240.00
ASF-100-07	32" Wall Mount	\$250.00
ASF-100-08	32" Post Mount	\$250.00



Assistant's Telescoping Arm

- Lightweight aluminum extrusion
- Breakaway, includes 2" post mount clamp
- Accepts hanger bars/trays with 1/2" pin
- Compatible with wall and surface mounts
- 24" - 40" Extension Arm, 8" Rigid Arm
- Leveling adjustment
- Powder-coat white finish
- Use VCF-339-02 to mount C.V. canister on end of arm

Catalog No.	Description	Retail Price
ASF-101-01	Assistant's Telescoping Arm	\$310.00



Post Mount Folding Arm / Holder Assembly

- Includes 12" folding arm
- Post mount bracket with swivel for 2" post
- Holder bar, 3 universal cavities
- Accepts additional folding arm sections
- Holder bar rotates 360°
- Powder-coat white finish

Catalog No.	Description	Retail Price
ASF-101-02	Post mount folding arm/holder assembly	\$230.00



Post Mount Telescoping Arm / Holder Assembly

- Includes ASF-101-01 Telescopic Arm
- Post Mount Clamp
- CV Mounting Block
- Holder Bar, 3 Universal Cavities
- Holder bar rotates 360°
- Powder-coat white finish

Catalog No.	Description	Retail Price
ASF-101-03	Post mount telescoping arm/holder assembly	\$459.00



Folding Arm Sections

- Folding Arm Sections only.
- 1/2" dia. pin mounts to a 2" dia. post mount clamp or wall mount
- 1/2" ID female end with bronze bushing
- Available in white (W), gray (G), or sterling (S) powdercoat finish
- Sections 12" and 16" also available in anodized (A)

Catalog No.	Description	Retail Price
ASF-102-01	6" Section	\$62.00
ASF-102-02	8" Section, Dual Pivot Arm	\$64.00
ASF-102-03	12" Section	\$70.00
ASF-102-04	16" Section	\$76.00



NEW

Wall Mounts For 1/2" Dia. Pin



103-01



103-11



103-12

- Designed for mounting Folding and Telescoping Arm Systems to 2" dia. post
- Accepts 1/2" O.D. pin.

Catalog No.	Description	Retail Price
ASF-103-01	Vertical Mount for 1/2" dia. pin, sterling, white or gray	\$58.80
ASF-103-11	Horizontal Mount for 1/2" dia. pin	\$58.80
ASF-103-12	Under Counter Mount for 1/2" dia. pin	\$58.80

Wall Mount For 2" Dia. Post Mount Arms

- Used to mount Post Mount Arm Systems
- Includes mounting screws

Catalog No.	Description	Retail Price
ASF-103-02	Wall Mount for 2" dia. post mount arms	\$135.00



Counter Top Mount

- Designed for mounting Folding and Telescoping Arm Systems to any flat surface
 - Can be drilled in the countertop from top or bottom (allow a 1-inch rise)
- Includes mounting screws
- Accepts 1/2" O.D. pin

Catalog No.	Description	Retail Price
ASF-103-03	Counter Top Mount	\$50.00



Cabinet Mount Bracket (X-Ray Column)

- Designed for use with x-ray column
- Accepts 2" diameter post

Catalog No.	Description	Retail Price
ASF-103-04	Cabinet Mount Bracket	\$185.00



Bracket For 1/2" Dia. Pin Mount

- Used with horizontal folding arms
- 1/2" pin protrudes from bottom of mount
- Includes mounting screws

Catalog No.	Description	Retail Price
ASF-104-01	Bracket for 1/2" dia. pin mount	\$40.00



1/2" Dia. Pivot Mount

- Designed for use with USM-200 and USM-400 series Delivery Units to mount onto horizontal folding arms
- 2" aluminium mount with 1/2" dia. pin
- includes mounting screws

Catalog No.	Description	Retail Price
ASF-104-02	1/2" dia. Pivot Mount	\$35.00



Rear Delivery Mounts

- Designed for use with rear delivery unit
- Mounts to existing cabinets
- Includes mounting screws

Catalog No.	Description	Retail Price
ASF-105-01	Single Mount	\$150.00
ASF-105-02	Double Mount	\$190.00



Horizontal Surface Mount Bracket

- Designed for mounting USM-200 and USM-400 series handpiece controls and other horizontal surfaces
- Includes mounting screws

Catalog No.	Description	Retail Price
ASF-106-01	Horizontal Surface Mount Bracket	\$20.00



Post Mount Clamps

- Mounts horizontal folding arms with 1/2" dia. pin to 1 3/4" or 2" dia. post

Catalog No.	Description	Retail Price
ASF-107-01	2" dia. Single Clamp	\$60.00
ASF-107-02	2" dia. Dual Clamp	\$80.00
ASF-107-11	1 3/4" dia. Single Clamp	\$70.00
ASF-107-12	1 3/4" dia. Dual Clamp	\$80.00



Clamp-On Water System & Accessory Adapters

- Adapts water systems or accessories to posts

Catalog No.	Description	Retail Price
ASF-107-03	2" Post Adapter, available in White powder coat (W) or Anodized aluminum (A)	\$64.00
ASF-107-04	2- 3/8" Post Adapter, White powder coat	\$70.00



Vacuum Canister Post Brackets

- For mounting vacuum canisters on a 2" O.D. post

Catalog No.	Description	Retail Price
VCF-342-01	Post Mount Adapter, anodized with gray bracket	\$64.00
VCF-342-02	Post Mount Adapter, white powder-coated, with gray bracket	\$64.00



Pivot Riser Mounting Kit

- Mounts above or below cabinet or cart surface and receives 1/2" dia. pin arm for special configuration of arms and trays
- Includes four 1/4" x 1" mounting screws

Catalog No.	Description	Retail Price
ASF-108-01	Pivot Riser Mounting Kit	\$90 .00



Top Post Monitor Mount

- 75 mm and 100 mm mounting patterns
- Fits into 2" dia. post
- Slotted to accommodate monitor cable inside post

Catalog No.	Description	Retail Price
ASF-109-01	Top Post Monitor Mount	\$360.00



Top Post Mount Brackets, 2"

- For mounting folding and telescoping arm assemblies
- Accepts 1/2" O.D. pin
- Select from White or Gray scratch resistant powder coat finish
- Fits radius style and most vertical posts
- Fits 2" diameter post

Catalog No.	Description	Retail Price
ASF-108-02	Top Post Mount	\$136.50



Flat Panel Monitor Mount, Vertical Post Mount

- Utilizes 1/2" pin for use wth folding rigid arms & mounts
- Fits either 75mm or 100mm VESA patterns
- Includes 2" vertical post mount
- Tilts up & down 30°
- Keyboard tray includes 8" dual pivot arm
- Monitor & keyboard have quick adjustment knob
- Swivels left to right 55°

Catalog No.	Description	Retail Price
ASF-109-02	Flat Panel Monitor Mount, Vertical Post Mount, white (W) or gray (G)	\$320.25
ASF-109-02K	Flat Panel Monitor Mount, Vertical Post Mount with full size Keyboard mount, white (W) or gray (G)	\$595.00



Flat Panel Monitor Mount, Top Post

Catalog No.	Description	Retail Price
ASF-109-04	Flat Panel Monitor Mount, Top Post, white (W) or gray (G)	\$340.00
ASF-109-04K	Flat Panel Monitor Mount, Top Post Mount with full size Keyboard mount, white (W) or gray (G)	\$880.00



Wall Mount Board

- Designed for extra security when mounting Wall Mounts
- Measures 12" x 12"

Catalog No.	Description	Retail Price
ASF-110-01	Wall Mount Board	\$154.00



Ring Mount Flexible Arm Systems

- Flexible Counterbalanced Arm System with pneumatic brake
- Provides 27" vertical travel using the 24" flexible arm
- Accommodates arm mounted handpiece control systems with tubing bundles concealed and protected within the arm
- Adjustable for weights from 3 to 13½ lbs
- Designed for use with 2" dia. post

Catalog No.	Description	Retail Price
ASF-111-01	26" Rigid w/24" Flexible Arm, Horizontal Reach 66"	\$800.00
ASF-111-02	26" Rigid w/18" Flexible Arm, Horizontal Reach 40"	\$800.00
ASF-111-03	22" Rigid w/24" Flexible Arm, Horizontal Reach 44"	\$800.00
ASF-111-04	22" Rigid w/18" Flexible Arm, Horizontal Reach 38"	\$800.00



Top Post Mount Flexible Arm Systems

- Flexible Counterbalanced Arm System with pneumatic brake
- Provides 27" vertical travel using the 24" flexible arm
- Adjustable for weights from 3 to 13½ lbs.
- Designed for use with 2" dia. post
- Accommodates arm mounted handpiece control systems with tubing bundles concealed and protected within the arm

Catalog No.	Description	Retail Price
ASF-112-01	26" Rigid w/24" Flexible Arm, Horizontal Reach 66"	\$800.00
ASF-112-02	26" Rigid w/18" Flexible Arm, Horizontal Reach 40"	\$800.00
ASF-112-03	22" Rigid w/24" Flexible Arm, Horizontal Reach 44"	\$800.00
ASF-112-04	22" Rigid w/18" Flexible Arm, Horizontal Reach 38"	\$800.00
ASF-112-05	12" Rigid w/24" Flexible Arm, Horizontal Reach 32"	\$800.00
ASF-112-06	12" Rigid w/18" Flexible Arm, Horizontal Reach 26"	\$800.00



Flexible Arm

- Flexible Counterbalanced Arm System with pneumatic brake
- Provides 27" vertical travel using the 24" flexible arm
- Accommodates arm mounted handpiece control systems with tubing bundles concealed and protected within the arm
- Adjustable for weights from 3 to 13½ lbs.
- Designed for use with either top post or ring mounted rigid arm

Catalog No.	Description	Retail Price
ASF-121-02	24" Flexible Arm	\$520.00
ASF-121-03	18" Flexible Arm	\$520.00



2" Dia. Ring Mount Rigid Arm

- Accommodates all lengths of Flexible Arms with tubing bundles concealed within the arm
- Designed for use with 2" dia. post

Catalog No.	Description	Retail Price
ASF-122-01	26" Ring Mount Rigid Arm	\$280.00
ASF-122-02	22" Ring Mount Rigid Arm	\$280.00



Top Post Mount Rigid Arm

- Accommodates all lengths of Flexible Arms with tubing bundles concealed within the arm
- Designed for use with 2" dia. control post with 1½" I.D. bearing or Wall Mount Cat No. ASF-133-01

Catalog No.	Description	Retail Price
ASF-123-01	26" Top Post Mount Rigid Arm	\$280.00
ASF-123-02	22" Top Post Mount Rigid Arm	\$280.00
ASF-123-03	12" Top Post Mount Rigid Arm	\$280.00



Wall Mount For Top Mount Arms

- Designed for mounting all Top Mount Arms to any flat vertical surface and provides leveling adjustment
- Mounting plate 6" x 6"

Catalog No.	Description	Retail Price
ASF-133-01	Wall Mount III for Top Mount Arms	\$175.00



Wall Mount For Post Mount Arms

- Designed for mounting a 2" dia. post to any flat vertical surface and provides leveling adjustment
- Mounting plate 6" x 6"

Catalog No.	Description	Retail Price
ASF-134-01	Wall Mount IV for Post Mount Arms	\$185.00



Utility Shelf Assembly

- 12" x 12" Utility Shelf pivots on mounting bracket to provide easy access
- 12" x 12" storage surface
- Powder coated finish
- Mounting bracket

Catalog No.	Description	Retail Price
ASF-155-01	2" dia. Post Mount or H-frame mount	\$170.00
ASF-155-03	Wall Mount	\$170.00



Cabinet/Cart Tray Holder Assembly

- Mounts above or below cabinet or cart surface
- Stainless Steel Tray, 9-3/4" x 13-1/2"
- 12" pivoting arm
- Mounting arm with screws

Catalog No.	Description	Retail Price
ASF-156-01	Cabinet/Cart tray holder assembly	\$330.00



Swing Tray

- For use with Flexible or Telescoping Arm Systems
- Swing Tray has a 13" pivot arm with 1/2" pin
- 12" x 12" work surface

Catalog No.	Description	Retail Price
ASF-160-01	Swing Tray	\$102.00



Tub & Tray Holder

- For use with Flexible or Telescoping Arm Systems
- Pivots on a rigid arm to provide easy access
- Powder coated finish
- Utility shelf, 9 1/2" x 12 1/2"
- 12" rigid arm

Catalog No.	Description	Retail Price
ASF-161-01	Utility Shelf	\$150.00



Tray Mount with Holder Bar

- For use with horizontal arm with 1/2" dia. pin mount
- Holds 9-3/4" x 13-1/2" tray

Catalog No.	Description	Retail Price
ASF-161-02	Tray mount with holder bar	\$165.00



SS Tray Holder with Holder Bar

- Designed for use with Folding or Telescoping Arm Systems
- May be used as a holder for plastic or stainless steel trays
- Stainless steel tray holder, 9 3/4" x 13 1/2"
- Wrap around holder bar, 1/2" dia.

Catalog No.	Description	Retail Price
ASF-162-01	Tray Holder with Hanger Bar	\$198.00



Tray Holder

- Designed for use with Folding or Telescoping Arm Systems
- May be used as a holder for plastic or stainless steel trays
- Stainless steel tray holder, 9 3/4" x 13 1/2"
- Mounted bracket with 1/2" O.D. pin

Catalog No.	Description	Retail Price
ASF-163-01	Tray Holder	\$164.00



Stainless Steel Tray

Catalog No.	Description	Retail Price
ASF-164-00	SS Tray, 10 ⁵ / ₈ " x 15 ¹ / ₈ "	\$63.30
ASF-164-01	SS Tray, 9 ³ / ₄ " x 13 ¹ / ₂ "	\$39.50
ASF-164-02	Non-slip pad 11 ³ / ₄ " x 6 ¹ / ₄ "	\$5.00



Plastic Molded Tray

- Molded plastic with non-slip bottom
- Available colors: black, white, and gray

Catalog No.	Description	Retail Price
ASF-165-01	Tray, 9 ⁵ / ₈ " x 13 ³ / ₈ "	\$18.00
ASF-165-02	Tray, 6 ³ / ₈ " x 9 ¹ / ₄ "	\$9.00



Plastic Utility Tray

- Molded plastic with non-slip bottom
- Available in white, beige and gray colors

Catalog No.	Description	Retail Price
ASF-171-01	2" Post Mount with post adapter	\$120.00
ASF-171-02	Utility Tray post adapter for smaller posts (1 ³ / ₈ " dia. to 1 ⁷ / ₈ " dia.)	\$4.00
ASF-171-03	8" Extension arm	\$33.00
ASF-171-04	Wall/Cabinet mount	\$33.00



Post Mount Tray Table

- Rugged molded plastic
- Swivels for optimum placement
- Post mount, adaptor selection included
- Approximate size 14" x 10" oblong
- Surf finish

Catalog No.	Description	Retail Price
ASF-172-01	Post Mount Tray Table, 14" x 10" oblong, surf finish	\$102.00



Aluminum Utility Tray, Post Mount

- Clamps onto 2" diameter post
- 10" x 10" surface area
- Sturdy, heavy gauge aluminum
- Tough powder-coated finish

Catalog No.	Description	Retail Price
ASF-172-02	Post Mount Utility Tray, aluminum, powder coated	\$127.00



U-Type Mobile Stand

- Satin Chrome finish
- Adjustable height, 23" to 32"
- 2" diameter with 3" x 3" mounting plate
- 15° angled upright
- Twin wheel casters
- Stand Size 18"W x 18"D

Catalog No.	Description	Retail Price
ACF-581-01	U-Type Mobile Stand	\$370.00



H-Type Mobile Stand

- Satin chrome finish
- Stable "H" frame base
- Adjustable vertical height from 25" to 41"
- 2" diameter post with 3" x 3" mounting plate
- Twin wheel casters for carpet or hard floor
- Stand size - 20"W x 18"D

Catalog No.	Description	Retail Price
ACF-582-01	H-Type Mobile Stand	\$370.00



Wheel Casters

- Smooth rolling, dual wheel
- Caster 2" diameter, pin 7/16" diameter x 1" long
- Non-marking

Catalog No.	Description	Retail Price
ACF-582-03	Soft wheel for hard floor caster, gray (G) and surf (S)	\$8.00/ea
ACF-582-04	Soft wheel for hard floor caster, gray (G) and surf (S)	\$40.00/pkg of 5
ACF-582-05	Hard wheel for carpet, gray (G)	\$12.00/ea
ACF-582-06	Cart Caster Insert	\$8.00/ea



Post / Cart Frame End Cap

- 2" Round End Cap. Fits inside most 2" O.D. posts, designed for all carts

Catalog No.	Description	Retail Price
ACF-584-01	Post or Cart Frame End Cap, gray	\$3.70



Reducing Sleeve

- Delrin spacer for enlarging diameter of post to accept 2" diameter clamp

Catalog No.	Description	Retail Price
ACF-584-03	Reducing Sleeve, 2" to 1.5", pkg of 2	\$19.50
ACF-584-04	Reducing Sleeve, 2" to 1.7", pkg of 2	\$19.50
ACF-584-05	1/16" Nyliner to expand to a 2" post	\$19.50



Post Clamp

- Designed to adapt light post to control post

Catalog No.	Description	Retail Price
ACF-592-01	Post Clamp	\$94.00



Light Post Extension Kits

- Scratch resistant powder coat finish
- Available in White (W), Gray (G) or Metal (M)
- Specify bushing - P&C or Ritter

Catalog No.	Description	Retail Price
ACF-592-02	24" Post, includes bushing, clamp, 8' extension cord,	\$180.00
ACF-592-03	30" Post, includes bushing, clamp, 8' extension cord,	\$185.00



COMPONENTS

Arms / Trays

Light Posts (Stainless Steel)

- All posts are 2" dia. durable stainless steel
- Available in 5 lengths

- Designed to be chair adapter mounted
- Light bushings sold separately - see item ACF-594-01

Catalog No.	Description	Retail Price
ACF-593-01M	Light Post, 24", Stainless Steel	\$136.00
ACF-593-03M	Light Post, 30", Stainless Steel	\$145.00
ACF-593-04M	Light Post, 40", Stainless Steel	\$190.00
ACF-593-06M	Light Post, 48", Stainless Steel	\$210.00
ACF-593-07M	Light Post, 60", Stainless Steel	\$219.00



Post Cap (Stainless Steel)

- Post Cap, 2" with O-Ring, Stainless Steel
- Use with Light Posts
- 2" Diameter

Catalog No.	Description	Retail Price
ACF-593-01MC	Post Cap, SS	\$28.00



Unit/Light Post

- Powder coat finish
- 2" diameter

Catalog No.	Description	Retail Price
ACF-593-01	Light Post, 24"	\$140.00
ACF-593-02	Light Post, 28" with clamp	\$170.00
ACF-593-06	Light Post, 48"	\$150.00
ACF-593-07	Light Post, 60" with light bearing	\$171.00



Bushing

- Designed to adapt light post to control post

Catalog No.	Description	Retail Price
ACF-594-01	Bushing, 1.832 OD x 1.700 ID (Fits Celux, Pelton, Adec, Belmont)	\$25.00
ACF-594-02	Bushing, used inside ACF-594-01 above, 1.700 OD x 1.380 ID (Fits Celux, Pelton, Adec, Belmont)	\$27.50
ACF-594-03	Bushing, 1.825 OD x 1.380 ID (Fits DentalEz, Faro, Forest, Ritter)	\$30.00



Chair Adapters

- Durable polyester powder coat finish
- Available for a wide variety of chairs and equipment
- Specify model needed

Catalog No.	Description	Retail Price
ACF-600-01	Chair Adapters for various models (A-dec, Belmont, DentalEz, Pelton & Crane, Marus)	\$300.00



Floor Plates

- Powder coated metal

Catalog No.	Description	Retail Price
ACF-601-01	Floor Plate, 6"	\$230.00
ACF-601-02	Floor Plate, 8"	\$240.00



Autoclavable Valve Core Syringe

- Complete handle and head assembly are fully autoclavable
- Quick-Change tip disconnect
- Accepts standard autoclavable tips or disposable tips
- Patented quick-disconnect turns air and water off when syringe is removed
- Precise control of air, water and spray
- Syringe with tubing also available

Catalog No.	Description	Retail Price
SYR-STSD-AS	Autoclavable Syringe with tip (includes Quick Disconnect, QD holder)	\$153.00



Autoclavable Continental Syringe

- Complete handle and head assembly are fully autoclavable
- Quick-Change tip disconnect
- Accepts standard autoclavable tips or disposable tips
- Patented quick-disconnect turns air and water off when syringe is removed
- Precise control of air, water and spray
- Syringe with tubing also available

Catalog No.	Description	Retail Price
SYR-STSD-AC	Autoclavable Continental Syringe with tip (includes Quick Disconnect, QD holder)	\$153.00



Precision Comfort™ Syringe

- Patented ergonomic design & precision touch button technology
- More control & delivery of air and/or water
- Softer buttons
- Increased angle & button size
- Syringe with tubing also available
- Made in USA

Catalog No.	Description	Retail Price
SYR-STSD-N	Precision Comfort with indexing tip	\$71.00



Precision Comfort Syringe™
Trademark of DCI International

Quick-Clean® Syringes

- Available in Standard, Continental and Valve Core styles
- Smooth contours for quick, thorough cleaning and a sleek, contemporary look
- Use standard autoclavable tips
- Precise control of air, water and spray
- Syringe with tubing also available
- Made in USA

Catalog No.	Description	Retail Price
SYR-STSD	Quick-Clean® Standard Syringe with tip	\$71.00
SYR-VCSD	Quick-Clean® Valve-Core Style Syringe with tip	\$75.00
SYR-CSD	Quick-Clean® Continental Syringe with tip	\$71.00
SYR-STST	Standard Syringe with tip (made in Taiwan)	\$50.00



Quick-Clean® - Registered Trademark of DCI International

Quick-Change Autoclavable Syringe Tips, Universal

Catalog No.	Description	Retail Price
SYR-233-01	Syringe tips, made in Taiwan	\$4.00
SYR-234-01	Syringe Tips, pkg of 5, made in USA	\$23.50
SYR-234-02	Syringe Tips, pkg of 50, made in USA	\$200.00



Autoclavable Indexing Syringe Tips

Catalog No.	Description	Retail Price
SYR-234-03	Indexing tips, pkg of 5	\$27.50
SYR-234-05	Indexing tips, pkg of 50	\$255.00



A-dec® Style Syringe Tips, Autoclavable

Catalog No.	Description	Retail Price
SYR-235-01	Syringe tip, each	\$4.90
SYR-235-02	Syringe tips, pkg of 5	\$24.00
SYR-235-03	Syringe tips, pkg of 50	\$205.00

* Note: These parts are not made by or in any way associated with A-dec®.



Seal-Tight® Disposable Air/Water Syringe Tips

Catalog No.	Description	Retail Price
SYR-236-00	Seal-Tight® Disposable Air/Water Syringe Tips, pkg of 200	\$48.00
SYR-236-0A	A Style Adapter (A-dec)	\$14.40
SYR-236-0M	M Style Adapter (Engle-Marco)	\$14.40
SYR-236-0D	D Style Adapter (Dental-EZ)	\$14.40
SYR-236-0P	P Style Adapter (Press Ring)	\$14.40

NEW



Repair Kit For Autoclavable Standard Syringe

• Syringe Repair Kit includes o-rings, springs, tool

Catalog No.	Description	Retail Price
SYR-031-00	Syringe Repair Kit for Autoclavable Valve Core Syringes	\$8.80



Repair Kit For Autoclavable Continental Syringe

• Syringe Repair Kit includes o-rings, springs, tool

Catalog No.	Description	Retail Price
SYR-031-01	Repair Kit For Autoclavable Continental Syringe	\$9.80



Repair Kit For Standard & Continental Syringes

• Syringe Repair Kit includes o-rings, springs, tool

Catalog No.	Description	Retail Price
SYR-031-02	Repair Kit For Standard & Continental Syringes	\$8.00



Repair Kit For Valve Core Style Syringes


• Syringe Repair Kit includes o-rings, springs, tool

Catalog No.	Description	Retail Price
SYR-031-03	Repair Kit For Valve Core Style Syringes	\$8.20



Repair Kit For Precision Comfort Syringes

Catalog No.	Description	Retail Price
SYR-031-04	Repair Kit for DCI Precision Comfort Syringe	\$8.00

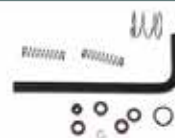

Syringe Replacement Cartridge

• Standard and Continental Style Button Cartridge with spring

Catalog No.	Description	Retail Price
SYR-241-01	Standard Syringe Replacement Cartridge, plated brass	\$9.00



Syringe Replacement Cartridge

• Valve Core Style Button Cartridge with spring

Catalog No.	Description	Retail Price
SYR-241-03	Valve Core Syringe Replacement Cartridge	\$9.20



Syringe Replacement Cartridge

- Precision Comfort Syringe Cartridge with spring

Catalog No.	Description	Retail Price
SYR-241-04	Precision Comfort Syringe Replacement Cartridge with spring, plated brass	\$9.20



Replacement Cartridge For Autoclavable Valve Core Syringe

Catalog No.	Description	Retail Price
SYR-241-05	Autoclavable Button Assembly, Syringe Button with Spring (U.S. Patent No. 5,199,871)	\$9.20



Replacement Cartridge For Autoclavable Continental Syringe

Catalog No.	Description	Retail Price
SYR-241-06	Replacement Syringe Button with spring for Continental Autoclavable Syringe	\$11.00



Syringe Service Tool

Catalog No.	Description	Retail Price
SYR-242-01	Syringe Service Tool, for removal of cartridge, tip adapter and barbs	\$2.50



3-In-1 Syringe Service Tool

Catalog No.	Description	Retail Price
SYR-242-02	3-in-1 Syringe Service Tool, includes Valve Core tool, Quick-Change adapter tool, Button Removal tool	\$8.00



3-In-1 Syringe Service Tool

Catalog No.	Description	Retail Price
SYR-242-03	Syringe Button Removal Tool	\$12.00



Flow Adjustment Block

Provides air/water metering to syringe. Includes 1/16" & 1/8" barbs, washers and sleeves. Available in plastic or plated brass

Catalog No.	Description	Retail Price
SYR-243-01	Flow Adjustment Block	\$18.90



Syringe O-Ring Kit, Universal

Includes all o-rings and service tool.

Catalog No.	Description	Retail Price
SYR-245-02	Autoclavable Syringe O-ring kit	\$15.50



Parts & Tools For A-Dec® Syringes



Catalog No.	Description	Retail Price
SYR-246-01	Tool to fit A-Dec® Syringes (ref: A-dec #98.0010.00)	\$16.00
SYR-246-02	Syringe Button Screw (ref: A-dec #23.1193.01), pkg of 5	\$26.80
SYR-246-03	Brass Swivel Sleeve (ref: A-dec #23.0886.00), pkg of 5	\$3.80
SYR-246-04	Valve Core (ref: A-dec #026.033.01), pkg of 10	\$6.00
SYR-246-05	Conical Button Spring (ref: A-dec #013.038.01), pkg of 10	\$6.90
SYR-246-06	Button Assy to fit A-Dec® Syringes (ref: A-dec 23.1232.00)	\$14.00
SYR-246-07	Repair kit for A-Dec® Soft Touch Syringe with 2 buttons	\$26.00
SYR-246-08	Service Kit to Fit A-Dec® Lever Syringes (ref: A-dec #23.1093.00)	\$22.00

*Not made by A-Dec® or under license or supervision of A-Dec®



Syringe Water Heaters

- Compact 4oz. reservoir type water heater with manual or air-actuated on/off control
- Temperature pre-set at 105° for comfort and quick recovery
- May be base or side mounted
- Not recommended with clean water systems
- Reservoir-type heater, 200watt, 1.67 amp
- 10-32 x 1/16" barbs with sleeves

Catalog No.	Description	Retail Price
SYR-007-06	110V Automatic with on/off auto switch	\$320.00
SYR-007-07	220V Automatic with on/off auto switch	\$320.00
SYR-007-08	110V Manual with on/off switch	\$320.00
SYR-007-09	220V Manual with on/off switch	\$320.00



Heated Syringe Upgrade

- Heat water to existing syringes
- Includes low voltage transformer, heated tubing(s) & air-electric switch
- Ultra-fast recovery
- Straight or coiled
- gray tubing

Catalog No.	Description	Retail Price
SYR-007-10	Heated Syringe Upgrade, complete kit	\$202.00
SYR-007-10P	2V Power Supply only, 110V/12Vac 40VA/3.3 Amp	\$63.40
SYR-006-xx	Replacement Heated Syringe Tubing available tubing colors: gray (14), sand (54), sterling (72)	\$69.00



2-Hole Sheathed Syringe Tubing

(See tubing section for all syringe tubings)

Catalog No.	Description	Retail Price
2S-005-xx	7' straight tubing with 5' sheath, each available in gray (14), sand (54), sterling (72)	\$13.42



3-Hole Sheathed Syringe Tubing

Catalog No.	Description	Retail Price
3S-005-xx	3-Hole Sheathed Syringe Tubing, price per 7 ft length available in gray (14), surf (92)	\$8.74



Standard Disc Type

- Die cast base, chrome plated cover
- 7' polyurethane 2-hole tubing (4-hole tubing with relay)
- Light touch, variable speed, high volume
- Available in two styles



Color Code (xx)

- Gray (14)
- Sterling (54)
- Black (01)
- Sand (72)
- Surf (92)



Catalog No.	Description		Retail Price
FCF-251-01-x-SI	Without Relay	Style I	\$75.00
FCF-251-01E-x-SII	Without Relay	Economy, Style II	\$68.00
FCF-251-01S-x-SII	Without Relay	Standard, Style II	\$77.00
FCF-251-02-x-SI	With Relay	Style I	\$110.00
FCF-251-02S-x-SII	With Relay	Standard, Style II	\$107.00



Chip Blower Type

Includes all the features of the standard type above plus toe operated momentary valve for chip blower, 7' polyurethane 3-hole tubing (5-hole tubing with relay)

Catalog No.	Description		Retail Price
FCF-252-01-x-SI	Without Relay	Style I	\$170.00
FCF-252-01-x-SII	Without Relay	Style II	\$130.00
FCF-252-02-x-SI	With Relay	Style I	\$187.00
FCF-252-02-x-SII	With Relay	Style II	\$160.00



Scaler Disc Type

Includes all the features of the standard type above plus toe operated momentary valve for on/off control of ultrasonic scaler (3-way, momentary), 7' polyurethane 3-hole tubing (5-hole tubing with relay)

Catalog No.	Description		Retail Price
FCF-252-03-x-SII	Without Relay		\$130.00
FCF-252-04-x-SII	With Relay		\$160.00



Wet/Dry Only Type

Includes all the features of the standard type above plus toggle valve for on/off control of water coolant, 7' polyurethane 3-hole tubing (4-hole tubing with relay).

Catalog No.	Description		Retail Price
FCF-253-01-x-SI	Without Relay	Style I	\$170.00
FCF-253-01-x-SII	Without Relay	Style II	\$130.00
FCF-253-02-x-SI	With Relay	Style I	\$187.00
FCF-253-02-x-SII	With Relay	Style II	\$160.00



Wet/Dry with Chip Blower Type

Includes all the features of the standard type above plus toggle valve for on/off control of water coolant and toe operated momentary valve for chip blower, 7' polyurethane 4-hole tubing (5-hole tubing with relay).

Catalog No.	Description		Retail Price
FCF-261-01-x-SI	Without Relay	Style I	\$180.00
FCF-261-01-x-SII	Without Relay	Style II	\$165.00
FCF-261-02-x-SI	With Relay	Style I	\$197.00
FCF-261-02-x-SII	With Relay	Style II	\$195.00



A-Dec® & Midmark Replacement Foot Control

- Foot controls - wet/dry & scaler button
- A-DEC® & Midmark replacement (A-DEC® #38.1805.00)

Catalog No.	Description	Retail Price
FCF-270-01	A-dec® & Midmark Replacement Foot Control, Wet/Dry & Scaler Button, dark surf color only.	\$195.00



Chip Blower Kit

Designed for adding chip blower feature to existing foot controls (not for use with standard disc type). Includes chip blower valve and mounting hardware.

Catalog No.	Description	Retail Price
FCF-263-01-SI	2 Way chip blower Style I	\$70.00
FCF-263-01-SII	2 Way Chip blower Style II	\$24.00
FCF-263-02-SI	3 Way chip blower Style I	\$50.00
FCF-263-02-SII	3 Way chip blower Style II	\$24.00

Note: Order foot control tubing separately, see TUBING.



Wet/Dry Kit

Designed for adding wet/dry feature to existing foot controls (not for use with standard disc type). Includes wet/dry toggle valve and mounting hardware.

Catalog No.	Description	Retail Price
FCF-264-01-SI	Toggle wet/dry kit Style I	\$35.00
FCF-264-01-SII	Toggle wet/dry kit Style II	\$24.00

Note: Order foot control tubing separately, see TUBING.



Foot Control Toggle Only

- Replacement Toggle

Catalog No.	Description	Retail Price
FCF-265-01-SI	FC toggles for style I Foot Controls, plastic	\$50.00/pkg of 10
FCF-265-02-SI	FC toggles for style I Foot Controls, metal	\$58.00/pkg of 10
FCF-265-03K-SII	FC Toggle Kit for style II Foot Controls, gray	\$4.50/ea



Foot Control Service Kit (Drive Air Valve Kit)

- Includes plunger, poppet, o-rings, and springs


Catalog No.	Description	Retail Price
FCF-852-01-SI	Complete kit for style I Foot Controls	\$15.00



Foot Control Components, Standard

- Replacement parts for standard Foot Controls

Catalog No.	Description	Retail Price
FCF-852-02	Foot Control Components	\$15.00

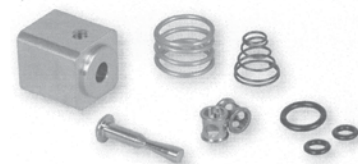

Foot Control Repair Kit

- For Style I wet/dry disc-type foot controls
- Includes plunger, poppet, o-rings, & springs

Catalog No.	Description	Retail Price
FCF-854-01-SI	Wet/dry valve kit, toggle	\$21.00



Foot Control Repair Kit

- For Style II standard disc-type foot controls

Catalog No.	Description	Retail Price
FCF-855-01-SII	Repair Kit for Style II, Disc Type	\$7.80



Foot Control Repair Kit

- For Style II chip blower disc-type foot controls

Catalog No.	Description	Retail Price
FCF-856-01-SII	Repair Kit for Style II, Chip Blower Disc-Type	\$8.90



Foot Control Repair Kit

- For Style II wet/dry disc-type foot controls

Catalog No.	Description	Retail Price
FCF-857-01-SII	Repair Kit for Style II, Wet/Dry Disc-Type	\$11.50



Foot Control Ring

- Fits all Style II Foot Controls

Catalog No.	Description	Retail Price
FCF-858-01-SII	Foot Control Ring, Gray	\$9.00



Foot Control Ring

- Fits all Style II Foot Controls

Catalog No.	Description	Retail Price
FCF-859-01-SII	Signal Relay Repair Kit	\$6.50



Marus Foot Control Repair Kit

- For Repairing Older Marus Style Foot Controls
- 55R527, 55R531, 55R535 Foot Controls
- Marus OEM # 85R186

Catalog No.	Description	Retail Price
FCF-859-02	Marus Foot Control Repair Kit	\$44.00



Replacement Parts For A-Dec® Foot Controls

Catalog No.	Description	Retail Price
FCF-860-01	Parts to fit A-dec® Foot Control III, Service Kit, (A-dec® ref #90.0593.00)	\$24.00
FCF-860-02	Parts to fit A-dec® Foot Control II, Complete Repair Kit, (A-dec® ref #90.0312.00)	\$16.50
FCF-860-03	Parts to fit A-dec® Foot Control I, Complete Rebuild Kit, (A-dec® ref #90.0010.00)	\$17.00
FCF-860-04	Parts to fit A-dec® Foot Control, Lever Foot Control Service Kit, (A-dec® ref #90.0313.00)	\$22.00



** Parts not made by A-Dec®

See Holder Section for various mounting options

Arm/Surface Mount Assistant's Package

- Deluxe HVE and SE with asepsis tubing
- 3-Way contour syringe & flow control block
- Extra port for additional HVE valve
- Integrated vacuum canister with locking lid
- 1/2" mounting pin allows a variety of mounting solutions
- 7' umbilical tubing
- 1/4" female quick disconnect with auto shut-off
- See Arms/Trays section for 1/2" pin mounting options
- Available in 4 tubing colors – gray, sand, sterling, surf

Catalog No.	Description	Retail Price
VCF-200-01P	Arm/Surface Mount Assistant's Package	\$595.00



Swing Mount Assistant's Package

- Under counter mounting
- Includes QD autoclavable SE & HVE with aseptic tubing
- Syringe with 6' aseptic tubing & flow control block
- Solids canister with screen
- 1 Piece molded swinging holder
- 10' Umbilical


Catalog No.	Description	Retail Price
VCF-200-01S	Swing Mount Assistant's Package	\$620.00



Wall & Cabinet Mount Assistant's Packages

- Vacuum canister with locking lid
- 8" Flush mount holder bar & holders
- Choice of Economy or Deluxe QD valves - 2 autoclavable HVEs and 1 autoclavable SE
- Smooth aseptic suction tubing, 7'

Catalog No.	Description	Retail Price
VCF-300-01E	Economy Package, 3 position, with 2 HVE's & SE, No Syr, with holder bar & holders	\$245.00
VCF-300-01D	Deluxe Package, 3 position, with 2 HVE's & SE, No Syr, with holder bar & holders	\$255.00

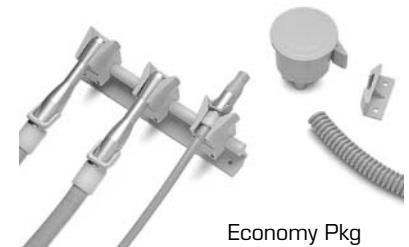

Economy Pkg

Precision Comfort Assistant's Packages

- Precision Comfort® QD HVE's & SE include 7' of tubing
- Syringe includes 7' of tubing & flow control block
- Vacuum canister, bracket mounted
- Smooth vacuum tubing, 7'

Catalog No.	Description	Retail Price
VCF-200-01	Precision Comfort Syringe, HVE & SE	\$290.00
VCF-200-02	Precision Comfort Syringe, 2 HVE's & SE	\$355.00
VCF-200-03	2 Precision Comfort HVE's & SE	\$190.00



Deluxe Assistant's Package

- Deluxe HVE's & SE with metal lever and 7' of tubing
- Syringe includes 7' of tubing & flow control block
- Vacuum canister, bracket mounted
- Smooth vacuum tubing 7'

Catalog No.	Description	Retail Price
VCF-201-01	Quick Clean Syringe, Deluxe HVE & SE	\$290.00
VCF-201-02	Quick Clean Syringe, 2 Deluxe HVE's & SE	\$355.00
VCF-201-03	2 Deluxe HVE's & SE	\$190.00



Economy Assistant's Package

- Economy HVE's & SE with plastic gray lever and 7' of tubing
- Syringe includes 7' of tubing & flow control block
- Vacuum canister, bracket mounted
- Smooth vacuum tubing 7'

Catalog No.	Description	Retail Price
VCF-202-01	Quick Clean Syringe, Economy HVE & SE	\$280.00
VCF-202-02	Quick Clean Syringe, 2 Economy HVE's & SE	\$345.00
VCF-202-03	2 Economy HVE's & SE	\$180.00



Deluxe Autoclavable QD HVE's

- All anodized aluminum finish lever
- Autoclavable valve features finger tip on/off and flow control with unique lever operated design
- Swivel quick disconnect allows 360° rotation of entire valve and adapts to any 1/2" I.D. central vacuum tubing
- After valve is installed in tubing, it can be snapped apart at the swivel to allow autoclaving of valve assembly (lever, body, and o-rings)
- Universal valve has one piece adapter added for use with universal type handpiece holders
- Accepts 7/16" O.D. tips
- Available with Standard Barb (S), Curved Barb (C), or Tapered Barb (T)

Catalog No.	Description	Retail Price
VCF-302-00	Contour Vacuum Valve, Slider HVE	\$68.00
VCF-302-01	Deluxe Autoclavable Vacuum Lever, Standard	\$45.00
VCF-302-02	Deluxe Autoclavable Vacuum Lever, Universal (specify delrin or aluminum barb)	\$47.00



Deluxe Autoclavable QD Saliva Ejector

- Anodized aluminum finish lever
- Features finger tip on/off and flow control with unique lever operated design
- Swivel quick disconnect works with universal holders allowing 360° rotation of entire valve and adapts to any 3/8" O.D. saliva ejector tubing
- After valve is installed in tubing, it can be snapped apart at the base to allow autoclaving of valve assembly.

Catalog No.	Description	Retail Price
VCF-303-00	Deluxe Autoclavable Saliva Ejector, No Tip Required	\$40.00
VCF-303-01	Deluxe Saliva Ejector Valve with lever, push on tip	\$42.00



Shown with tapered barb, gray tip.
Also available with straight barb and black tip.

Elite SS Saliva Ejector Valve, Lever Action



- Includes same features as the autoclavable saliva ejectors but in stainless steel material

Catalog No.	Description	Retail Price
VCF-305-00	Hard Tip SE Valve, SS, Lever Action with Viton® O-Rings, no SE-screen	\$22.30
VCF-305-01	Hard Tip SE Valve, SS, Lever Action with Viton® O-Rings, straight white delrin barb, no SE-screen	\$26.50
VCF-305-02	SE Valve, SS, Lever Action with Viton® O-Rings, straight black delrin barb	\$30.20
VCF-305-03	SE Valve, SS, Lever Action with Viton® O-Rings, straight white delrin barb	\$30.20

All valves are also available with 6' tubing

DCI Precision Comfort® Autoclavable Premium Valves

- Ergonomic single piece metal lever
- New O-ring placement fits Evacu-Tip & Tip-A-Dilly™
- No rubber tip required

Catalog No.	Description	Retail Price
Extended Vacuum Valves used with Standard Holder		
VCF-306-01	Single Valve w/ Quick-Disconnect	\$49.00
VCF-306-02	Replacement Quick-Disconnect Swivel for 1/2" Tubing	\$13.00
VCF-306-03	Replacement Lever & Spool Assembly	\$20.00

Short Vacuum Valves used with Vacuum Holder		
VCF-307-01	Single Valve w/ Quick-Disconnect	\$46.00
VCF-307-02	Replacement Quick-Disconnect Swivel for 1/2" Tubing	\$13.00
VCF-307-03	Replacement Lever & Spool Assembly	\$20.00


Saliva Ejectors used with Standard Holder		
VCF-308-01	Single Valve w/ Quick-Disconnect	\$43.00
VCF-308-02	Replacement Quick-Disconnect Swivel for 1/2" Tubing	\$11.00
VCF-308-03	Replacement Lever & Spool Assembly	\$20.00



Economy Autoclavable Vacuum Lever Valves

- Autoclavable valve features finger tip on/off and flow control with unique lever operated design
- Swivel quick disconnect allows 360° rotation of entire valve and adapts to any 1/2" I.D. central vacuum tubing
- After valve is installed in tubing, it can be snapped apart at the swivel to allow autoclaving of valve assembly (lever, body, and o-rings)
- Universal valve has one piece adapter added for use with universal type handpiece holders
- Accepts 1/16" O.D. tips
- Also available with 6' tubing
- Specify delrin or aluminum barb

Catalog No.	Description	Retail Price
VCF-301-01	Autoclavable Vacuum Lever, Standard	\$34.00
VCF-301-02	Autoclavable Vacuum Lever, Universal	\$34.00

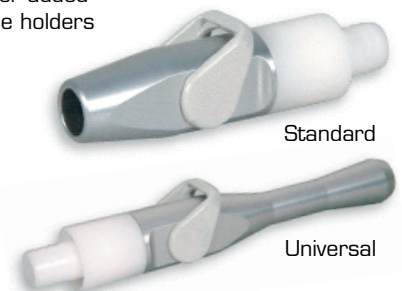

Economy Autoclavable Vacuum Lever Valves

- With Lever or Traditional Rotary On/Off
- 360° swivel quick-disconnect
- Lightweight comfortable body
- Hangs by tip in syringe holders
- Easy to disassemble for cleaning

Catalog No.	Description	Retail Price
VCF-309-00	HVE with lever	\$40.00
VCF-309-01	HVE with Rotary On/Off	\$36.00




NEW

Economy Autoclavable Saliva Ejectors

- Swivel quick disconnect works with universal holders allowing 360° rotation of entire valve and adapts to any 3/8" O.D. saliva ejector tubing
- After valve is installed in tubing, it can be snapped apart at the base to allow autoclaving of valve assembly. SE with lever features finger tip on/off and flow control with unique lever operated design
- Also available with 6' tubing
- Shown here with white delrin, tapered barb, gray tip
- Also available with straight barb, with aluminum barb, with black tip

Catalog No.	Description	Retail Price
VCF-311-01	Saliva Ejector Valve with lever	\$25.00
VCF-311-02	Saliva Ejector Valve without lever	\$15.00

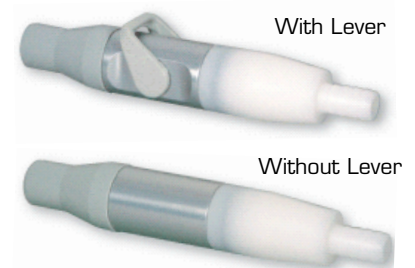

Saliva Ejector, Lever Action

- Saliva ejector lever action valve with push-on tip
- Available with straight or tapered barb

Catalog No.	Description	Retail Price
VCF-311-03	SE, Lever Action, straight barb, gray	\$24.00
VCF-311-04	SE, Lever Action, tapered barb, gray	\$24.00



Saliva Ejector, Rotary On/Off

- Saliva ejector rotary action valve with push-on tip
- Available with straight or tapered barb

Catalog No.	Description	Retail Price
VCF-312-01	SE, Rotary Action, straight barb, gray	\$18.00
VCF-312-02	SE, Rotary Action, tapered barb, gray	\$18.00



Saliva Ejector Tip Assembly

- Push-on tip without shut-off

Catalog No.	Description	Retail Price
VCF-313-01-G	SE with push-on tip, gray	\$19.00
VCF-313-01-B	SE with push-on tip, black	\$11.90



Saliva Ejector Adapter

- Adapts saliva ejector tip for use in central vacuum valve

Catalog No.	Description	Retail Price
VCF-314-01	Saliva Ejector Adapter, plastic	\$7.00
VCF-314-02	Saliva Ejector Adapter, aluminum	\$8.00
VCF-314-03	Adapter only, plastic	\$4.00
VCF-314-04	Adapter only, aluminum	\$5.00



Replacement A-dec® Style Vacuum Valves

Catalog No.	Description	Retail Price
VCF-315-01	Autoclavable Extended Vacuum Valve	\$70.00
VCF-315-02	Autoclavable Short Vacuum Valve	\$64.00
VCF-315-03	Autoclavable Saliva Ejector	\$64.00

** These parts are manufactured to fit A-dec® equipment, not under license or supervision of A-dec®.



Precision Comfort® Conversion QD Swivels

Catalog No.	Description	Retail Price
VCF-316-01	A-dec® HVE QD Swivel	\$20.60
VCF-316-02	Pelton & Crane HVE QD Swivel	\$20.60
VCF-316-03	Cattani HVE QD Swivel, 26mm Swivel, 18mm Barb	\$24.70
VCF-316-04	Cattani SE QD Swivel, 20mm Swivel, 11mm Barb	\$20.60
VCF-316-05	Durr HVE QD Swivel, 28mm Swivel, 22mm Barb	\$24.70
VCF-316-06	Siemens HVE QD Swivel, 19mm Swivel, 15mm Barb	\$24.70

** These parts are manufactured to fit A-dec® equipment, not under license or supervision of A-dec®.



Tubing Adapter – 1/2" ID x 3/8" OD Barb

- Connects 3/8" ID vacuum tubing to the spare HVE port on Vacuum canisters

Catalog No.	Description	Retail Price
VCF-320-01	Tubing Adapter, 1/2" ID tube x 3/8" OD barb, black	\$8.80



SE Tubing Adapter – 7/16" OD x 7/32" Barb

- Allows connection of SE and tubing
- Inserts into vacuum canister port or HVE valve

Catalog No.	Description	Retail Price
VCF-320-02	SE Tubing Adapter, 7/16" OD x 7/32" barb, white	\$5.20



Autoclavable Saliva Ejector Tip (325°F)

- For use with saliva ejector assemblies

Catalog No.	Description	Retail Price
VCF-321-01	Push-On tip, gray	\$2.90
VCF-321-02	Push-On tip, black	\$2.90



Saliva Ejector Check Valve

- Prevents back-flow
- Easy to install barbed ends
- Long lasting plastic
- Non-restricting

Catalog No.	Description	Retail Price
VCF-322-01	SE Check Valve	\$14.00



Saliva Ejector to HP Holder Adapter

- Slips over your old SE base to allow use in an auto holder

Catalog No.	Description	Retail Price
VCF-322-02	SE to HP Holder Adapter	\$4.50



Saliva Ejector Screen

- Fits push-on or screw-on style saliva ejector assemblies

Catalog No.	Description	Retail Price
VCF-323-01	Plastic Screen	\$0.60/ea
VCF-323-02	Metal Screen	\$9.00/ea



Lever

- Allows for on/off and flow control with use of finger tip
- Both levers are autoclavable

Catalog No.	Description	Retail Price
VCF-325-01	Saliva Ejector Lever	\$3.50
VCF-325-02	Central Vacuum Lever	\$4.00



COMPONENTS

Vacuum Components

Plug

- Stops suction and allows saliva ejector and central vacuum tubings to be placed in universal holder after valve has been removed from quick disconnect base

Catalog No.	Description	Retail Price
VCF-326-01	Plug	\$12.00



Central Vacuum Canister – Bottom Port

- Standard Vacuum Canister, 1/2" Drain with Locking Lid, large
- Removable solids collector
- Two high velocity ports for 1/2" I.D. tubing
- One saliva ejector port for 3/8" I.D. tubing
- Extra port plugs
- Assembly includes body, lid with o-ring and screen

Catalog No.	Description	Retail Price
VCF-332-01	Canister Assembly, bottom port, gray	\$24.00



Central Vacuum Canister – Side Port

- Standard Vacuum Canister, 1/2" Drain with Locking Lid, large
- Removable solids collector
- Two high velocity ports for 1/2" I.D. tubing
- One saliva ejector port for 3/8" I.D. tubing
- Assembly includes body, lid with o-ring and screen

Catalog No.	Description	Retail Price
VCF-333-01	Canister Assembly, side port, gray	\$24.00



Central Vacuum Canister

- Standard Vacuum Canister, 1/2" Drain with Push Lid, small
- Features removable solids collector
- Two high velocity ports for 1/2" I.D. tubing
- One saliva ejector port for 3/8" I.D. tubing
- Optional water drain port for 1/8" I.D. tubing
- Assembly includes body, lid with o-ring, screen, vertical & horizontal brackets

Catalog. No.	Description	Retail Price
VCF-335-01	Canister Assembly, push lid, small, gray	\$22.00



Central Vacuum Canister

- Standard Vacuum Canister, 1/2" Drain with Twist Lid, small
- Removable solids collector
- Two high velocity ports for 1/2" I.D. tubing
- One saliva ejector port for 3/8" I.D. tubing
- Optional water drain port for 1/8" I.D. tubing
- Assembly includes body, lid with o-ring, screen, vertical & horizontal brackets

Catalog. No.	Description	Retail Price
VCF-336-01	Canister Assembly, twist lid, small, gray	\$22.00



Vacuum Screens

- Replacement Screen for Central Vacuum Canister

Catalog No.	Description	Retail Price
VCF-337-01	Covered Screen, clear, small - 1 ⁷ / ₈ "	\$11.45/pkg of 25
VCF-332-02	Covered Screen, clear, large - 2 ¹ / ₈ "	\$11.45/pkg of 25



Solids Collector

- For use with central vacuum canisters
- Green fits #VCF-335-01 and #VCF-336-01 canister assemblies
- Blue fits #VCF-332-01 & #VCF-333-01 canister assemblies

Catalog No.	Description	Retail Price
VCF-331-04	Solids Collector 1 ⁷ / ₈ ", green	\$12.50/pkg of 25
VCF-331-04B	Solids Collector 1 ⁷ / ₈ ", green	\$34.00/pkg of 100
VCF-333-02	Solids Collector 2 ¹ / ₈ ", blue	\$12.50/pkg of 25
VCF-333-02B	Solids Collector, 2 ¹ / ₈ ", blue	\$34.00/pkg of 100



Solids Collector

Catalog No.	Description	Retail Price
VCF-333-05	Fits A-dec® 300/500, (A-dec® #11.1284.00)	\$49.45/pkg of 144
VCF-333-06	Fits A-dec® Cascade & Performer Screen, Disposable (A-dec® #11.1191.00)	\$50.00/pkg of 144



Solids Collector - Pelton & Crane®

Catalog. No.	Description	Retail Price
VCF-333-07	Solids Collector, Pelton & Crane® Spirit Disposable Screen (#9431388)	\$50.00



Central Vacuum Canister Plugs

Catalog. No.	Description	Retail Price
VCF-334-01	Saliva ejector valve port	\$0.26
VCF-334-02	Central vacuum valve port	\$0.34



Vacuum Canister Mounting Brackets

- Replacement Brackets for Central Vacuum Canisters

Catalog No.	Description	Retail Price
VCF-338-01	Bracket only, horizontal	\$3.00
VCF-338-02	Bracket only, vertical	\$3.00



Vacuum Canister Mounting Brackets

Catalog No.	Description	Retail Price
VCF-339-00	Horizontal surface mounting bracket, allows 1-1/14" umbilical tubing	\$25.00
VCF-339-01	Wall or Counter Mounting Bracket	\$16.50
VCF-339-02	Under Counter Mounting Bracket	\$14.00



Vacuum Canister Post Brackets

- Allows either central vacuum canister, #VCF-332-01 or #VCF-333-01, to be mounted to a 2" O.D. post.

Catalog No.	Description	Retail Price
VCF-342-01	Post Mount Adapter, anodized with gray bracket	\$64.00
VCF-342-02	Post Mount Adapter, white powder-coated, with gray bracket	\$64.00



In-line Manual On/Off Valve

- Manual operated vacuum on/off valve may be added in-line to central vacuum tubing
- Enables easy removal of vacuum canister lid without shutting vacuum down

Catalog No.	Description	Retail Price
VCF-345-01	1/2" x 1/2" Barbed fittings	\$19.00
VCF-345-02	5/8" x 5/8" Barbed fittings	\$24.00
VCF-345-03	1/2" x 5/8" Barbed fittings	\$24.00



Autoclavable Dry Oral Cup (Plastic)

- Fits all vacuum valves with 7/16" dia. tips

Catalog No.	Description	Retail Price
VCF-351-01	Dry Oral Cup	\$9.30
VCF-351-02	Liners for Dry Oral Cup, 1000/pkg	\$170.00



Combo-Tiptm Premium Disposable Evacuator Tips

- Vented and Non-Vented Ends on the same Evacuator tip
- Unique "S" Tip with no sharp edges for maximum comfort
- Latex Free

Catalog No.	Description	Retail Price
VCF-352-01	Evacuator Disposable Tips, white	\$14.00/ pkg of 100



Comfort Plus Saliva Ejectors

- Rounded edges, smooth surfaces and perfectly placed suction ports
- Enhance patient comfort while protecting mucosal tissue
- Latex Free

Catalog No.	Description	Retail Price
VCF-352-02	SE Disposable Tips, white	\$14.00/pkg of 100



Cuspidor Strainer

- Molded white plastic
- Available in two sizes

Catalog No.	Description	Retail Price
VCF-353-01	Cuspidor Screen with white lid, small, 1-5/8" Fits: A-dec Cascade and Performer cuspidor bowls	\$8.70/pkg of 25
VCF-353-02	Cuspidor Screen with gray lid, large, 2-1/2" Fits: A-dec, Belmont, Dabi Atlante, Dansereau, DCI, Den-Tal-Ez, Forest, Knight, Marus, MDT, Proma Ritter and Schein	\$8.70/pkg of 25



Cuspidor Bowl

- Quality porcelain.
- Bowl 8" OD x 6¾" ID
- Stem 2½" OD x 1½" ID

574-02



574-03



574-05

574-06



Catalog No.	Description	Retail Price
ACF-574-01	Cuspidor Bowl, with 2 o-rings and screen	\$96.50
ACF-574-02	Cuspidor Bowl only	\$91.50
ACF-574-03	Screen only, A-dec® style	\$5.00
ACF-574-05	Cup Filler Spout, stainless steel	\$24.00
ACF-574-06	Cuspidor Bowl Rinse Spout, stainless steel	\$28.00
ACF-575-01	Cuspidor Time Delay Kit	\$83.00

575-01



Cuspidor Vacuum Drain Kit

- Converts any cuspidor to vacuum drain
- Accepts 5/8" drain line, adaptable
- Fully automatic, silent running
- In-line filter included

Catalog No.	Description	Retail Price
VCF-354-01	Cuspidor Vacuum Drain Kit	\$190.00
VCF-354-02	Workstation Vacuum Drain Kit	\$189.50
VCF-354-03	Diaphragm only	\$9.30
VCF-354-04	Replacement Screen	\$6.70



Vacuum Drain Assembly Kit

- Select cuspidor drain kit or a sink drain
- Eliminates the need for conventional plumbing
- Allows installation where a gravity drain is not available or practical
- Includes flush toggle for remote operation

Catalog No.	Description	Retail Price
VCF-355-01	Cuspidor Drain Kit	\$250.00
VCF-355-02	Workstation Vacuum Drain Kit	\$235.00



NEW

Air/Water Saliva Ejector

- Air/water operated saliva ejector
- Vacuum generator
- Saliva ejector tip with tubing, 6'
- Toggle valve
- ½" I.D. drain tubing, 9'
- Installation hardware

Catalog No.	Description	Retail Price
VCF-361-01	Air/Water SE Assembly, gray	\$160.70



Air/Water Saliva Ejector Venturi

- Air/Water Saliva Ejector Venturi
- Plated brass construction

Catalog No.	Description	Retail Price
VCF-362-01	Air/Water Saliva Ejector Venturi	\$49.50



Air Venturi

- Venturi with insert for air only saliva ejector
- Simple to connect
- Plastic outer with brass ports

Catalog No.	Description	Retail Price
VCF-362-02	Air/ Venturi, plastic	\$32.50



Air Vacuum Venturi

- Air Vacuum Venturi For SE and HVE
- Venturi only
- Plated brass construction.

Catalog No.	Description	Retail Price
VCF-362-03	Air Vacuum Venturi	\$91.70



Foot Control Activated Air Vacuum System

- Air powered oral evacuation
- Requires 50 PSI minimum to operate
- 1 quart plastic collection bottle
- Mounts to under counter or on the wall, hardware included
- Air connection uses 1/4" OD supply tubing
- Deluxe HVE with 6' of smooth vacuum tubing included
- Requires Scaler Button Style Foot Control to operate - not included

Catalog No.	Description	Retail Price
VCF-372-01	Foot Control Activated Air Vacuum System	\$370.00
VCF-372-01B	Replacement 1 quart collection bottle	\$30.00



Holder Activated Air Vacuum System

- Air powered oral evacuation
- Requires 50 PSI minimum air pressure to operate
- 1 quart plastic collection bottle
- Mounts to underside of dental unit or cabinet
- Deluxe HVE with 6' smooth vacuum tubing included
- Activated by holder & air switch


Catalog No.	Description	Retail Price
VCF-372-02	Holder Activated Air Vacuum System	\$370.00



HVE & SE Cleaning Brushes

- Valve cleaning brushes for HVE (3-3/4" bristles) and for SE (3" bristles)

Catalog No.	Description	Retail Price
VCF-374-01	HVE Cleaning Brushes	\$10.00/pkg of 5
VCF-374-02	SE Cleaning Brushes	\$11.00/pkg of 5

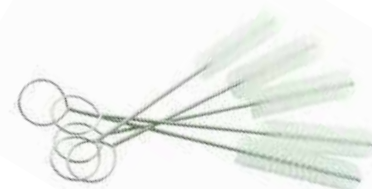

Automatic Handpiece Holder

- Two-piece HP holder assembly with miniature valve and air operated actuator
- Provides air signal to automatic control block for holdback when handpiece is placed into the holder
- Available with On/Off lockout toggle with yellow indicator
- Non-marring, plastic holder
- Aluminum mounting adapter for 1/2" dia. holder bar
- Gray or Black, Normally Closed

Catalog No.	Description	Retail Price
HOF-381-00	Auto HP Holder without Toggle, Normally Closed, gray only	\$37.60
HOF-381-01	Auto Holder with Lockout Toggle, Normally Closed gray or black	\$40.00



Automatic Accessory Holder

- Two-piece Accessory holder assembly with miniature valve and air operated actuator
- Provides air signal when accessory is out of holder
- Available with On/Off lockout toggle with yellow indicator
- Non-marring, plastic holder
- Aluminum mounting adapter for 1/2" dia. holder bar
- Gray or Black, Normally Open
- * * Accessory holders are not compatible with handpiece use.

Catalog No.	Description	Retail Price
HOF-381-10	Auto Accessory Holder without Toggle, Normally Open, gray only	\$37.60
HOF-381-11	Auto Accessory Holder with Lockout Toggle, Normally Open, gray or black	\$42.00



Seamless (1 piece) Automatic Handpiece Or Accessory Holder

- Miniature signal valve and actuator
- Durable, one-piece plastic holder
- Designed for 1/2" bar-mounting
- Includes 10-32 mounting screws
- 5/8 diameter bell
- Normally Open or Normally Closed Valve

Catalog No.	Description	Retail Price
HOF-381-02	Auto Holder, Normally Open, gray only	\$40.00
HOF-381-03	Auto Holder, Normally Closed, gray only	\$40.00



Asepsis Style Oval Holder

- Molded Nylon Construction
- Resistant to Discoloration
- Includes Two 10-32 Mounting Screws
- Normally Closed Valve or Normally Open Valve
- Other colors available upon request

Catalog No.	Description	Retail Price
HOF-381-04	Asepsis Oval Auto Holder, Normally Open, white or gray	\$33.00
HOF-381-05	Asepsis Oval Auto Holder, Normally Closed, white or gray	\$33.00



Electric Automatic Holder with Toggle

Holder assembly has a low voltage switch for direct activation of an electrical accessory while its handpiece is out of the holder. Normal position is when handpiece, vacuum valve or accessory is in holder.

- U.L. listed low voltage switch with 24" pigtail
- On/off lockout toggle with indicator
- Non-marring, plastic holder
- Aluminum mounting adapter for 1/2" dia. holder bar
- Normally Open or Normally Closed Valve



Catalog No.	Description	Retail Price
HOF-391-01	Electric Auto Holder w/Lockout Toggle, Normally Open, gray or black	\$50.00
HOF-391-05	Electric Auto Holder, w/Lockout Toggle, Normally Closed, gray or black	\$50.00
HOF-392-05	Electric Auto Vacuum Holder, w/Lockout Toggle, Normally Closed, gray only	\$55.00

Electric Automatic Holder without Toggle

Holder assembly has a low voltage switch for direct activation of an electrical accessory while its handpiece is out of the holder. Normal position is when handpiece, vacuum valve or accessory is in holder.

- U.L. listed low voltage switch with 24" pigtail
- Non-marring, plastic holder
- Aluminum mounting adapter for 1/2" dia. holder bar
- Normally Open or Normally Closed



Catalog No.	Description	Retail Price
HOF-391-03	Electric Auto Holder, Normally Open, gray	\$40.20
HOF-391-07	Electric Auto Holder, Normally Closed, gray	\$40.20

Asepsis Electric Automatic Holder without Toggle

Asepsis style electric auto holder assembly has a low voltage switch for direct activation of an electrical accessory while its handpiece is out of the holder. Normal position is when handpiece, vacuum valve or accessory is in holder.

- Resists chemical discoloration
- Without lockout, normally open or closed
- Activates electrical accessories
- Available in white (W) or gray (G)



Catalog No.	Description	Retail Price
HOF-391-04	Asepsis Style Electric Auto Holder, Normally Open	\$34.00
HOF-391-08	Asepsis Style Electric Auto Holder, Normally Closed	\$34.00
HOF-391-04-SII	Asepsis Style Electric Auto Holder, Normally Open, style II	\$50.00

Standard Two-Piece Holder Assembly

- Plastic with aluminum mount back
- Adjustable tilt angle
- For 1/2" dia. holder bars
- Available in gray or black

Catalog No.	Description	Retail Price
HOF-393-01	Standard Holder Assembly for HP, Syringe, SE's, & Universal HVE, gray or black	\$18.00
HOF-394-01	Standard HVE Holder, holds vacuum valves from base, gray or black	\$18.00



Universal Holder One-Piece

- Holder assembly for use with handpiece, syringe, saliva ejector or universal central vacuum valve.
- One-piece, non-marring plastic holder
- Adapts to a 1/2" dia. holder bar
- Available in gray or black

Catalog No.	Description	Retail Price
HOF-401-01	Universal holder with stainless steel screws	\$14.50



Universal Holder, Asepsis Style

- Asepsis Oval Holder assembly for use with handpiece, syringe, saliva ejector or universal central vacuum valve
- One-piece, molded plastic holder
- Adapts to 1/2" dia. holder bar
- Resistant to discoloration
- Choice of colors: White or Gray

Catalog No.	Description	Retail Price
HOF-401-03	Universal Asepsis style holder with stainless steel screws	\$20.00



Asepsis HVE Holder

- Smooth, seamless aseptic design
- Exclusive non-yellowing material
- Easy to wipe down
- Fits 1/2" diameter holder bars

Catalog No.	Description	Retail Price
HOF-401-04	Asepsis HVE Holder	\$17.00



Asepsis Standard Holder

- Asepsis Oval Holder assembly for use with handpiece, syringe, saliva ejector or universal central vacuum valve with tip extension
- One-piece, molded plastic holder
- Adapts to 1/2" dia. holder bar; adjustable tilt angle
- Resistant to discoloration
- Choice of colors: White, Gray, Lt. Sand

Catalog No.	Description	Retail Price
HOF-401-05	Asepsis Standard Holder	\$30.00



Syringe Holder

- Standard flat back, without step, drilled

Catalog No.	Description	Retail Price
HOF-402-01	Syringe holder, flat, gray or black	\$6.00



Saliva Ejector Holder

- Standard flat back, with step, not drilled

Catalog No.	Description	Retail Price
HOF-403-01	SE holder, flat, gray or black	\$7.00



HVE Holder

- Standard flat back, drilled

Catalog No.	Description	Retail Price
HOF-404-06	HVE Holder, flat, gray or black	\$7.00



Economy HP Holder

- Use for handpieces, saliva ejectors, syringes
- Standard round back with step
- Mounting bracket included
- Available in black or gray

Catalog No.	Description	Retail Price
HOF-404-01	Economy HP holder	\$6.00



Economy HVE Holder

- Mounting bracket included
- Use for holding HVE by the base
- Standard round back with step
- Available in black or gray

Catalog No.	Description	Retail Price
HOF-404-07	Economy HVE holder	\$7.00



Flat Mount Holder Adapter

- Connects holders to a flat plate
- Use on 1/8" or thinner material
- Gray plastic
- Fits most standard holders

Catalog No.	Description	Retail Price
HOF-405-01	Flat mount holder adapter	\$5.00



Automatic Holder Valve, Side Ported

- For use with automatic handpiece or automatic accessory holders

Catalog No.	Description	Retail Price
HOF-411-01	HP Holder Valve, side ported, normally closed	\$10.90
HOF-411-02	Accessory Holder Valve, side ported, normally open	\$10.90



Automatic Holder Valve, Rear Ported

- For use with automatic handpiece or automatic accessory holders

Catalog No.	Description	Retail Price
HOF-411-03	HP Holder Valve, rear ported, normally closed	\$11.50
HOF-411-04	Accessory Holder Valve, rear ported, normally open	\$11.50



Wall or Cabinet Mount Holders

- Mounts to any vertical surface
- Holds HVE by tip
- Adjustable tilt angle
- Durable plastic holder & machined aluminum mount

Catalog No.	Description	Retail Price
HOF-415-01	3 position holder bar	\$73.00
HOF-415-02	4 position holder bar	\$93.00



Holder Bar with Holders

- 3 or 4 position
- 1 HVE for 3 position / 2 HVE holders for 4 position
- 1 SE holder
- 1 Syringe holder
- 1/2" dia. pin mounts to horizontal folding arm

Catalog No.	Description	Retail Price
HOF-416-01	3 Position Holder Bar with holders & 1/2" dia. pin	\$68.00
HOF-416-02	4 Position Holder Bar with holders & 1/2" dia. pin	\$73.00



Holder Bar with 1/2" Dia. Pin

- 3 or 4 position
- Mounting holes for holders
- 1/2" dia. pin mounts to horizontal folding arms

Catalog No.	Description	Retail Price
HOF-417-01	3 position holder bar with 1/2" dia. pin	\$33.00
HOF-417-02	4 position holder bar with 1/2" dia. pin	\$35.00



Under Counter Holder Bar

- For use with molded asepsis holders
- Cabinet mounted holder bar
- Available in white or gray powder coat
- Bar is 10½" x 3½"

Catalog No.	Description	Retail Price
HOF-433-03	Under Counter Holder Bar	\$140.00



Swing Mount Assistant's Holders

- Mounts to any vertical or horizontal surface
- Holds HVE by tip
- 3 position – 11-3/4" long, 5-1/2" tall
- 4 position – 13-3/4" long, 5-1/2" tall
- Includes mount for vacuum trap

Catalog No.	Description	Retail Price
HOF-433-04	3 Position Swing Holder Bar	\$210.00
HOF-433-05	4 Position Swing Holder Bar	\$240.00



Holder Bars

- For use with automatic, universal, and vacuum handpiece holders.
- Vertical surface mount
- Precision milled aluminum
- Available in Gray - Painted (G) or Black - Anodized (B)

Catalog No.	Description	Retail Price
HOF-434-01	Single, 2"	\$14.50
HOF-434-02	Duo, 4"	\$17.00
HOF-434-03	Tri, 6"	\$19.00
HOF-434-04	Quad, 8"	\$21.00
HOF-434-05	Small Universal, 10"	\$25.00
HOF-434-06	Universal, 14"	\$28.00



Single Holder Bracket

- Durable plastic
- Mounting hardware included

Catalog No.	Description	Retail Price
HOF-436-01-x	Single Holder Bracket, gray * Also available in White (W) or Black (B), subject to custom color charge.	\$3.00



Economy Aluminum Holder Brackets

- Aluminum Holder Brackets for use with Economy Holders
- Precision milled aluminum
- Mounting hardware included

Catalog No.	Description	Retail Price
HOF-435-01	2-position, 2 ⁵ / ₈ "	\$25.00
HOF-435-02	3-position, 4 ¹ / ₄ "	\$30.00
HOF-435-03	4-position, 5 ⁷ / ₈ "	\$64.00



Wrap Around Holder Bars

- ½" dia. holder bar for custom applications.

Catalog No.	Description	Retail Price
HOF-441-01	12" Holder Bar	\$40.00
HOF-441-02	16" Holder Bar	\$40.00



Holder Bar Brackets

- For use with wrap around holder bars. Includes set screws.

Catalog No.	Description	Retail Price
HOF-442-01	Holder bracket with 1/2" hole	\$22.00
HOF-442-02	Standard holder bracket	\$22.00



Curing Light Holder

- For use within standard holders.

Catalog No.	Description	Retail Price
HOF-443-01	Curing Light Holder	\$36.00



Cavitron Handpiece Sleeve Assembly

- Adapts handpiece assembly to fit standard holders
- Molded plastic with set screws

Catalog No.	Description	Retail Price
HOF-444-01	Sleeve Assembly	\$9.00



Low Voltage Holder Switch

- Installs into most auto holders
- Normally open or Normally closed
- 6" leads

Catalog No.	Description	Retail Price
HOF-450-01	Low Voltage Holder Switch, N/O	\$12.40
HOF-450-02	Low Voltage Holder Switch, N/C	\$12.40



Asepsis Assistant's Holder Bar

- For use with Folding or Telescoping Arm System and accepts standard vacuum and universal holders
- Includes Holder Bar
- 1/2" O.D. pin

Catalog No.	Description	Retail Price
HOF-451-01	Holder Bar	\$57.50



Asepsis Assistant's Holder Bar

- Universal central vacuum holder for use with Folding or Telescoping Arm Systems
- Includes 1/2" O.D. pin

Catalog No.	Description	Retail Price
HOF-452-01	Tri (3 Universal)	\$136.00
HOF-452-02	Quad (4 Universal)	\$136.00
HOF-453-01	Tri (2 Universal, 1 Standard)	\$136.00
HOF-453-02	Quad (3 Universal, 1 Standard)	\$136.00
HOF-454-01	Tri (1 Universal, 2 Standard)	\$136.00
HOF-454-02	Quad (2 Universal, 2 Standard)	\$136.00



Surface/Side Mount Swing Holders

- Smooth, easy to wipe down
- Accepts barrier protection
- 180° rotation, full fore/aft pivot
- Choice of cavities hold SE, HVE by tip or HVEs by base
- Holds many combinations of instruments
- Universal bracket allows top, bottom or side mounting

Catalog No.	Description	Retail Price
HOF-454-10	Tri (3 Universal)	\$150.00
HOF-454-11	Quad (4 Universal)	\$150.00
HOF-454-12	Tri (2 Universal, 1 HVE base)	\$150.00
HOF-454-13	Quad (3 Universal, 1 HVE base)	\$150.00
HOF-454-14	Quad (2 Universal, 2 HVE bases)	\$150.00



Surface/Side Mount Swing Holder Bar

- Bar accepts up to 4 hangers of choice (not included)
- Rotates 180°
- Pivots for optimum instrument access
- Tough powder-coat finish
- Universal bracket allows top, bottom or side mounting

Catalog No.	Description	Retail Price
HOF-454-20	Surface Mount Swing Holder Bar	\$72.00



Pivoting Bent Holder Bar

- Mounts under cabinet or work surface
- ½" dia. mounting bar
- powder-coat
- sterling color

Catalog No.	Description	Retail Price
HOF-455-01	7", with 3 asepsis universal holders	\$160.00
HOF-455-02	7", Bar and Mount only	\$120.00
HOF-455-03	9", Bar and Mount only	\$135.00



Pivoting Holder Bar

- Mounts under cabinet or work surface
- ½" dia. mounting bar
- powder-coat
- sterling color

Catalog No.	Description	Retail Price
HOF-456-01	8½", Bar and Mount only	\$76.00
HOF-456-02	8½", with 3 asepsis universal holders	\$105.00
HOF-456-03	10½", Bar and Mount only	\$95.00



Cam Lock Holder Bar

- Mounts under cabinet or work surface
- ½" dia. mounting bar
- Converts left or right
- Rotational pivot

Catalog No.	Description	Retail Price
HOF-457-01	8½" with 3 Automatic asepsis HP Holders, N/C	\$265.00
HOF-457-02	8½" with 3 Universal asepsis holders	\$215.00
HOF-457-03	8½" with 4 Universal asepsis holders	\$230.00
HOF-457-04	8-1/2", without holders	\$175.00



4-Position Holder Bar

- One piece design, easy to mount
- Holds four universal instruments
- Gray, non-marring plastic

Catalog No.	Description	Retail Price
HOF-458-01	4-position holder bar	\$22.20



Toggle Valve, On/Off

- Positive on/off control
- Brass construction
- Nickel plated
- Port size 10-32
- Mounting hardware, 15/32-32

Catalog No.	Description	Retail Price
VAF-451-01	2-Way w/o Exhaust, gray or black	\$13.40
VAF-451-03	3-Way w/Exhaust, gray or black	\$13.40
VAF-451-05	Toggle only, gray or black	\$1.25



Toggle Valve, On/Off, Momentary

- Positive on/off control
- Brass construction
- Nickel plated
- Port size 10-32
- Mounting hardware, 15/32-32

Catalog No.	Description	Retail Price
VAF-452-01	2-Way w/o Exhaust, gray or black	\$13.40
VAF-452-03	3-Way w/Exhaust, gray or black	\$13.40
VAF-451-05	Toggle only, gray or black	\$1.25



Push Button, On/Off, Momentary

- Positive on/off control
- Brass construction
- Nickel plated
- Port size 10-32
- Mounting hardware, 15/32-32

Catalog No.	Description	Retail Price
VAF-453-01	2-Way w/o Exhaust N/C, gray or black	\$17.00
VAF-453-03	3-Way w/Exhaust, N/C, gray or black	\$17.00
VAF-453-04	3-Way w/Exhaust, N/O, gray	\$17.50



Push Button, On/Off, Momentary - Side Ports

Catalog No.	Description	Retail Price
VAF-454-01	2-way w/o Exhaust, N/O	\$17.35
VAF-454-02	3-way w/ Exhaust N/O	\$17.35
VAF-454-03	2-way w/o Exhaust, N/C	\$17.35
VAF-454-04	3-way w/ Exhaust N/C	\$17.35


Available push button colors: Gray (G), Black (B), Light Gray (LG), Plated (PT), White (W), Blue (BL), Red (R), Black Ball (BB), White Ball (WB)



Push Button, On/Off, Momentary - with Soldered Side Ports

Catalog No.	Description	Retail Price
VAF-455-01	2-way w/o Exhaust, N/C, with 2 soldered $\frac{1}{16}$ " barbs	\$18.00
VAF-455-02	3-way w/ Exhaust N/C, with 2 soldered $\frac{1}{16}$ " barbs	\$18.00

Available push button colors: Gray (G), Black (B), Light Gray (LG), Plated (PT), White (W), Blue (BL), Red (R), Black Ball (BB), White Ball (WB)

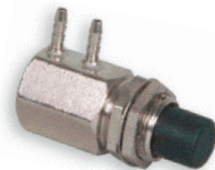

Needle Control Valve

- Variable control for air and water
- Brass construction
- Nickel plated
- Port size 10-32
- Mounting hardware, 15/32-32

Catalog No.	Description	Retail Price
VAF-461-01	Needle Control Valve with knob	\$15.50

Available knobs: Smooth Gray (SG), Smooth Black (SB), Ridged Gray (RG), Ridged Black (RB)



Pinch Flow Control Valve

Catalog No.	Description	Retail Price
VAF-462-01	Pinch Valve with knob, $\frac{9}{16}$ " hex valve, accepts $\frac{1}{8}$ " and $\frac{1}{4}$ " OD tubing	\$13.00


Available knobs: Smooth Gray (SG), Smooth Black (SB), Ridged Gray (RG), Ridged Black (RB)



Plastic Knobs

- Replacement knobs for control valves with $\frac{1}{4}$ " O.D. shaft.

Catalog No.	Description	Retail Price
VAF-464-01	Style I, gray or black	\$2.00
VAF-464-03	Style II, gray or black	\$2.00

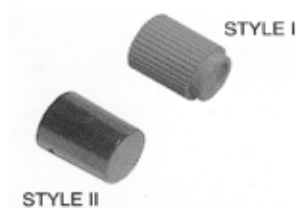

Flow Restrictor

- Variable restriction for air or water (pinch type)
- Accommodates tubing to $\frac{1}{4}$ " O.D.
- Brass construction, nickel plated
- Mounts directly on tubing

Catalog No.	Description	Retail Price
VAF-463-01	Flow restrictor	\$3.50



Momentary Cup Filler Valve

Catalog No.	Description	Retail Price
VAF-466-01	With cup-fill outlet fitting	\$22.70



Tracer Washers

- Use on panel mount valves
- Yellow for air, Blue for water
- Fits toggles, flow controls etc.
- Ideal for female quick-disconnects
- Non-corroding plastic

Catalog No.	Description	Retail Price
FIF-687-02	Blue Water washer for $\frac{1}{4}$ " QD, pkg of 10	\$3.10
FIF-687-03	Blue Water washer for $\frac{3}{8}$ " QD, pkg of 10	\$3.10
FIF-687-04	Yellow Air washer for $\frac{1}{4}$ " QD, pkg of 10	\$3.10
FIF-687-05	Yellow Air washer for $\frac{3}{8}$ " QD, pkg of 10	\$3.10



Toggle Valve Panel Extension

- Extends Valve Neck to 1"
- Fits Most Valves with 15/32-32" Thread

Catalog No.	Description	Retail Price
VAF-468-01	Toggle valve panel extension	\$12.50



Air Pilot, Valve, On/Off

- Positive on/off control
- Brass construction, nickel plated
- 1/16" Barb fitting
- Port size 10-32
- Panel mount, 15/32-32 thread with nut

Catalog No.	Description	Retail Price
VAF-471-01	2-way w/o exhaust, N/C	\$21.00
VAF-471-02	3-way w/exhaust, N/C	\$21.00
VAF-471-03	3-way w/exhaust, N/O	\$21.00 NEW



Mini Water Valve

- Positive on/off control
- Non-retracting
- Nickel plated, Brass construction
- 1/16" Barb fitting
- Port size 10-32

Catalog No.	Description	Retail Price
VAF-472-01	2-way, w/o exhaust, N/C	\$20.00



Check Valve

- May be installed in water coolant system of most handpiece controls to eliminate retraction
- Nickel plated, brass construction
- 1/16" or 1/8" Barb fitting
- Port size 10-32
- Flow rate, .5 SCFM @ 50 PSI
- Crackling pressure, 2.5 PSI
- Operating range, 0-150 PSI

Catalog No.	Description	Retail Price
VAF-473-01	Check valve, 1/16"	\$9.80
VAF-473-02	Check valve, 1/8"	\$9.80
VAF-473-03	Check Valve, 1/8" FPT	\$20.50



Shuttle Valve

- Designed with two inlets and one outlet
- Valve permits free flow from the high pressure side, while blocking the low pressure side
- For use with air or water
- Inport (shuffled), 2
- Outport, 1
- Nickel plated, brass construction
- Port size 10-32

Catalog No.	Description	Retail Price
VAF-474-01	Shuttle Valve, 1/16"	\$10.40
VAF-474-02	Shuttle Valve, 1/8"	\$10.40



Water Relay And Flow Control Combo Valve

- Non-retraction, Combination pilot valve & needle valve
- 10-32 side ports
- Metal construction
- Low activation pressure of 10-15 PSI
- Standard panel mount

Catalog No.	Description	Retail Price
VAF-475-01	Water relay and flow control	\$50.90
VAF-475-02	Water relay without flow control	\$42.50
VAF-475-03	Repair kit	\$9.00
FIF-629-01	Replacement Duck Bill Barb, each	\$5.50



Water Relay (Plastic)

- Most reliable water relay made, tested over 10 million cycles
- Air-activated On/Off for handpiece water coolant
- Non-retraction
- Plastic construction
- 1/16" Barb fitting

Catalog No.	Description	Retail Price
VAF-476-01	Water relay valve, plastic	\$34.00



Water Relay

- Designed for on/off control of handpiece water coolant
- Air activated water retraction
- Brass construction
- 1/16" Barb fitting
- Port size 10-32

Catalog No.	Description	Retail Price
VAF-482-02	Water Relay, 2-way without exhaust, N/C	\$19.90
VAF-482-02E	Quick Exhaust Valve 10-32 ports to adapt to water relay Releases air pressure immediately when not activated	\$20.60



Non-Retraction Water Relay

- Designed for on/off control of handpiece water coolant
- Air activated water relay, Non-retraction
- Positive on/off control
- Brass construction
- 1/16" Barb fitting
- Port size 10-32

Catalog No.	Description	Retail Price
VAF-482-03	2-way Water Relay Valve with Actuator , N/C	\$35.30
VAF-482-04	2-way Miniature Poppet valve (relay), N/C	\$18.20
VAF-482-05	Miniature air pilot actuator	\$19.00



Pilot Actuated Relay Valve

- Requires actuator item# VAF-482-05
- For air or water
- Brass construction
- Port size 10-32, 2-way (non relieving)

Catalog No.	Description	Retail Price
VAF-482-06	Pilot Actuated Relay Valve, N/O	\$23.60
VAF-482-08	Pilot Actuated Relay Valve, N/C	\$19.30



Pilot Actuated Valve

- Requires actuator item# VAF-482-05
- For air or water
- Brass construction
- 1/8" FPT ports, 2-way (non relieving)

Catalog No.	Description	Retail Price
VAF-482-07	Pilot Actuated Relay Valve, N/C	\$54.00



Routing Valve

- Routing control for air or water
- Brass construction, nickel plated
- Port size 10-32
- Mounting hardware 15/32-32

Catalog No.	Description	Retail Price
VAF-483-01	Routing Valve, gray or black	\$14.50

** Available in momentary, as push button, and type "C"



Routing Valve, Push Button Style

Catalog No.	Description	Retail Price
VAF-483-02	2-way w/o Exhaust, N/O	\$14.00
VAF-483-03	3-way w/ Exhaust, N/O	\$14.00

Available push button colors: Gray (G), Black (B), Light Gray (LG), Plated (PT), White (W), Blue (BL), Red (R), Black Ball (BB), White Ball (WB)



Air Activated Routing Valve

- Air-signal actuated—no toggle
- All plated brass body and valve
- 2-Way valve for air or water service

Catalog No.	Description	Retail Price
VAF-483-04	Air Activated Routing Valve	\$20.10



Cuspidor/Cup Filler Time Delay Valve Assembly

- Controls length of time bowl or cup fill is activated
- 1/4" poly inlet & outlet

Catalog No.	Description	Retail Price
VAF-483-05	Time Delay Valve	\$74.70
VAF-483-06	Time Delay, relay assembly	\$38.00



Valve Mounting Bracket

- Used for mounting toggles, indicators, QD's & air electric switches.
- 1/2" I.D. holes
- White powder-coat finish

Catalog No.	Description	Retail Price
VAF-484-01	Single	\$6.90
VAF-484-02	Dual	\$9.20
VAF-484-03	Triple	\$14.70



Toggle Mounting Bracket

- Used for mounting toggles & QD's
- Great for Cavitron supply
- One 1/2" I.D. mounting hole
- One 1 1/16" I.D. mounting hole
- White powder-coat finish

Catalog No.	Description	Retail Price
VAF-485-01	Toggle Mounting Bracket	\$12.25



Mini Automatic Handpiece Control Block

- Mini control blocks of 2 or 4 handpieces
- Solid brass construction, nickel plated
- Extra mounting holes allow a variety of placements
- 10-32 threaded barbs
- Longest lasting inter-woven diaphragm
- High flow design

Catalog No.	Description	Retail Price
VAF-490-00	Mini 2-HP Control Block	\$80.00
VAF-490-01	Mini 3-HP Control Block	\$90.00
VAF-490-02	Mini 4-HP Control Block	\$125.00
VAF-490-03	Mini 5-HP Control Block	\$210.00
VAF-490-04	Drive Air Adjusting Screw	\$5.70
VAF-490-05	Mounting bracket	\$9.50



Replacement Diaphragms

Catalog No.	Description	Retail Price
VAF-490-00D	2-HP Block Diaphragm	\$2.60
VAF-490-01D	3-HP Diaphragm	\$3.00
VAF-490-02D	4-HP Diaphragm	\$3.50


Automatic Handpiece Control Block

Forest Style miniature assembly for selection and operation of drive air and air/water coolants for 2 to 4 air driven handpieces. Designed for use with automatic handpiece holders.

Catalog No.	Description	Retail Price
VAF-491-01	2-HP Control Block	\$150.00
VAF-491-02	3-HP Control Block	\$130.00
VAF-491-03	4-HP Control Block	\$180.00

Replacement Diaphragms

Catalog No.	Description	Retail Price
SIF-844-01	2-HP Block Diaphragm	\$3.00
SIF-844-02	3-HP Block Diaphragm	\$3.00
SIF-844-03	4-HP Block Diaphragm	\$7.00

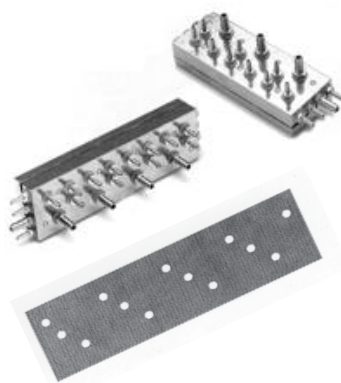

Piston/Pinch Water Valve

- Designed for use with any non-toggle style automatic handpiece holder
- No metal to water contact, Dependable, simple to maintain
- Replaces OEM unit water valves, 1/16" barbed signal fitting

Catalog No.	Description	Retail Price
VAF-492-02	Pinch Style (1 per HP)	\$13.70



Flex Valve

- Unique handpiece control block
- Activates drive air
- Air coolant & water coolant
- Excellent for asepsis applications
- Black nylon construction

Catalog No.	Description	Retail Price
FVD-101-01	Flex Valve, side port	\$13.70
FVD-101-02	Flex Valve, bottom port	\$13.70
FVD-101-03	Flex Valve, assembly	\$16.90



Manual Handpiece Control Block

- For use with manual handpiece control systems
- Includes all barb fittings
- Handpiece pressure adjustment
- Brass construction, nickel plated

Catalog No.	Description	Retail Price
VAF-492-05	Manual Control Block	\$24.30



Scaler Manifold with Flow Adjustments

Features: 6 x 10-32 outlet ports, Dual flow controls, Plated Brass, Multiple mounting configurations, Perfect for scaler or other uses

Catalog No.	Description	Retail Price
VAF-492-06	Scaler manifold with flow adjustment	\$40.00



Drive Pattern Block

Features: Mini distribution manifold, Mounts securely into chassis, Air or water service

Catalog No.	Description	Retail Price
VAF-492-07	Drive pattern block	\$5.20



Distribution Manifold

Features: Separate water and air banks, Internal water + air flow controls, 10-32 ports, metal

Catalog No.	Description	Retail Price
VAF-492-08-M	Distribution manifold	\$32.00



1" Mini Round Pressure Gauge

Features: 1" round face, 1/8" NPT / 10-32 Female
For use with air or water

Catalog No.	Description	Retail Price
VAF-493-01	1" dia. Round Pressure Gauge (0-60 psi)	\$15.30
VAF-493-02	1" dia. Round Pressure Gauge (0-100 psi)	\$15.30



1 1/2" Round Pressure Gauge

Features: 1 1/2" round face, 1/8" MPT thread / 10-32 Female
For use with air or water
Center back-mounted

Catalog No.	Description	Retail Price
VAF-493-04	1 1/2" dia. Round Pressure Gauge (0-160 psi)	\$13.70



COMPONENTS

Valves / Filters / Regulators

2" Round Gauge

- 2" round face, 1/4" MPT center back connection
- 0-160 PSI

Catalog No.	Description	Retail Price
VAF-493-05	2" dia. Round Pressure Gauge	\$19.00



Rectangular Gauge

- Two pressure ranges
- 1/2" high x 1 3/4" wide
- Panel mount
- 10-32 port
- For use with air only

Catalog No.	Description	Retail Price
VAF-494-01	Rectangular Gauge, 0-60 psi	\$30.00



1 1/2" Pressure Gauge, Bottom Mount

- Round 0-100 PSI pressure reading
- 1-1/2" diameter
- 1/8" MPT thread
- For air or water

Catalog No.	Description	Retail Price
VAF-494-03	1 1/2" Pressure Gauge, 0-100 psi, bottom mount	\$15.00
VAF-494-04	1 1/2" Pressure Gauge, 0-160 psi, bottom mount	\$14.00



Vacuum Gauge

- Vacuum Gauge
- 0-3 inches of mercury
- 1/4" NPT bottom-mount
- Fits wet type vacuum pumps
- Vibration and shock resistant

Catalog No.	Description	Retail Price
VAF-495-01	Vacuum Gauge	\$21.00



Panel Mount Gauge

- 0-100 PSI pressure range
- 1 3/4" diameter, use 1 5/8" hole
- 1/8" MPT & 10-32 female port
- For use with air or water

Catalog No.	Description	Retail Price
VAF-501-01	Panel Mount Gauge	\$17.00



Pneumatic Indicator

- 25-120 PSI operating range
- 10-32 port
- Mounting hardware, 15/32-32

Catalog No.	Description	Retail Price
VAF-502-01	Green Indicator, gray or black	\$12.80
VAF-502-03	Red Indicator, gray or black	\$12.80



Standard Utility Center

Utility Center, Standard with Housing & Cover (Black or Gray Cover)

- Mounts on floor or wall (plumbing template included)
- Slotted holes allow installation without disconnecting tubing
- Accommodates up to 1 $\frac{3}{4}$ " umbilicals
- Two master shut-off valves with $\frac{1}{2}$ " MPT to $\frac{3}{8}$ " compression fittings
- Stainless steel housing 14" x 9" x 4 $\frac{1}{2}$ "
- Two non-relieving regulators
- Two pressure gauges
- Two in-line filters
- Two manual shut-off valves



Catalog No.	Description	Retail Price
VAF-511-01N	Complete Utility Center	\$384.00
VAF-511-03N	Utility Center without Housing and Cover	\$283.50
VAF-511-04N	Housing and Cover only, standard	\$82.70
VAF-511-06N	Cover Only	\$18.40

Economy Utility Center

- Utility Center, Economy with Housing & Cover (Black or Gray)
- Economy utility center with housing and cover
- Provides filtered, regulated air and water for assistant's unit
- Mounts on floor or wall (plumbing template included)
- Slotted holes allow installation without disconnecting tubing
- Accommodates up to 1 $\frac{3}{4}$ " umbilicals
- Two integrated master shut-off valves with regulator and filter
- Two manual shut-off valves
- 3-way on/off toggle valve
- Stainless steel housing 14" x 9" x 4 $\frac{1}{2}$ "



Catalog No.	Description	Retail Price
VAF-512-01N	Complete Utility Center, Economy	\$304.00
VAF-512-02N	Utility Center without Housing and Cover	\$218.25
VAF-512-03N	Housing and Cover only, standard	\$82.70
VAF-512-04N	Cover Only	\$18.40

Utility Center and Master Valves

- Master valve assembly with filter and regulator
- Manual air and water shut-off valves
- Enameled metal housing
- Fitting and hole plug kit
- One 1 $\frac{3}{4}$ ", two $\frac{3}{4}$ " and one 2 $\frac{1}{2}$ " holes
- Utility size - 14" x 9" x 5 $\frac{1}{8}$ "

Catalog No.	Description	Retail Price
VAF-144-017	Complete Utility Center Assembly	\$400.00
VAF-144-019	Utility Center without housing	\$280.00
VAF-144-120	Utility Center housing only	\$120.00
VAF-144-121A	Master Valve Assembly with filter & regulator, for air	\$130.00
VAF-144-121W	Master Valve Assembly with filter & regulator, for water	\$130.00
VAF-144-122	Regulator Kit, includes VAF-144-121A & VAF-144-121W and toggle valve with tubing & plastic tee	\$280.00



Mini Utility Assembly

- Used for assistant's accessory installation
- Air only
- includes shut-off valve, regulator, and in-line filter

Catalog No.	Description	Retail Price
VAF-514-01	Mini Utility Assembly	\$63.30



COMPONENTS

Valves / Filters / Regulators

Regulator Kit

- Non-relieving regulator for air or water service
- Includes Regulator, brass body, In-line filter, Round gauge, Mounting bracket, 1/4" poly x 1/8" MPT connector, 1/8" MPT close nipple

Catalog No.	Description	Retail Price
VAF-522-01	Regulator Kit	\$66.30
VAF-521-01	Regulator only	\$42.00
VAF-521-02	Mounting Bracket & Nut	\$18.00



Regulator, NORGREN Brand

Catalog No.	Description	Retail Price
VAF-521-01N	Regulator Kit	\$35.00

- Non-relieving regulator for air or water service
- 1/8" NPT port size, brass body, plastic bonnet
- 5-100 psi, (0.3 to 7 bar),
- without gauge



Master Shut-Off Valves

- Master Shut-Off Valve w/Single Barb Cap
- For remote or utility center installation
- Master shut-off valves with
 - Pilot-actuated shut-off
 - 1/8" FPT inlet, 1/4" poly outlet
 - 10-32 barbs with sleeves to connect 1/8" O.D. tubing

Catalog No.	Description	Retail Price
VAF-521-10	Master Shut-off Valve with single-barb cap	\$42.00
VAF-521-11	Master Shut-off Valve with double-barb cap	\$42.00
VAF-521-13	Replacement Diaphragm	\$2.60



Dual Master Shut-Off Valve Kit

- For remote or utility center installation
- Two master shut-off valves
- Two 1/8" MPT to 3/8" copper compression fittings
- 3-way on/off toggle valve

Catalog No.	Description	Retail Price
VAF-521-12	Dual Master Shut-off Valve Kit	\$112.30
VAF-521-13	Replacement Diaphragm	\$2.60



Dual Master Shut-Off Valve Assembly

- Dual Master Shut-Off Valve w/Filters
- Provides shut-off for air and water
- One 1/4" barb fitting for air
- One 1/8" barb fitting for master control valve

- One 1/8" barb fitting for water
- 40 micron water and air filter
- 3/8" compression tube
- connection for air and water

Catalog No.	Description	Retail Price
VAF-521-15	Dual Master Shut-off Valve Assembly	\$115.60



Master Shut-Off with Regulator & Filter

- Complete pack for air or water service
- Includes filter and replacement element
- Adjustable mini-regulator
- 10-32 outlet, connects directly to angle stop

Catalog No.	Description	Retail Price
VAF-523-02	Master shut-off assembly	\$80.00
VAF-523-03	Air/Water Filter only	\$47.30
VAF-523-04	Cap Filter with o-ring, gray	\$8.00



Economy Master Shut-Off Valve Assembly, for Air or Water

- Select for air or water service
- Includes filter, auto shut-off and gauge
- Adjustable mini-regulator
- 10-32 outlet, connects directly to angle stop

Catalog No.	Description	Retail Price
VAF-523-00-A	Economy Master shut-off assembly with gauge, for air	\$105.00
VAF-523-00-W	Economy Master shut-off assembly with gauge, for water	\$105.00



In-line Filter Assembly

- Nickel plated brass construction, 1" dia. x 1" long
- 35 micron filter disc
- Port size, 1/8 FPT
- Fittings, 1/4 poly x 1/8 MPT included

Catalog No.	Description	Retail Price
VAF-524-00	Filter Assembly with 10-32 x 1/16" barbs	\$25.00
VAF-524-01	Filter Assembly with 1/4 poly x 1/8 MPT fittings	\$30.00
VAF-524-02	Filter Only	\$18.00
VAF-524-03	Replacement Bronze Filter Disc, pkg of 2	\$3.00
VAF-524-03S	Replacement Stainless Steel Filter Disc, each	\$24.00



Filter Element

- Fits #VAF-523-03 or any other mini filter
- 3/4" tall x 1 1/16" OD
- 40 Micron filtration

Catalog No.	Description	Retail Price
VAF-524-04	Filter Element only, each	\$2.10



In-line Air Filter (Plastic)

- Supplied with removable 1/8" barbs, as pictured
- 10-32 female ports accept alternate fittings
- Non-corrosive plastic housing
- 43 micron or 5 micron stainless steel filter screen
- 1/2" diameter x 5/8" long with barbs removed

Catalog No.	Description	Retail Price
VAF-524-05	In-line Air Filter, 43 micron	\$12.40
VAF-524-15	In-line Air Filter, 5 micron	\$12.40



Mini Fixed In-line Regulator

Catalog No.	Description	Retail Price
VAF-524-06	Mini Fixed 30 PSI In-line Regulator	\$20.70

- Constant 30 PSI output
- Rated input from 0-150 PSI
- Non-relieving
- In-line mounting
- 1/16" barbed fittings



Air Filter / Moisture Separator

Catalog No.	Description	Retail Price
VAF-525-01	Air Filter Moisture Separator, 5 micron filter	\$60.00
VAF-525-02	Moisture Separator Bowl only	\$25.50
VAF-525-03	Filter Element, 5 micron, black threads	\$28.00
VAF-525-05	Bowl O-ring, pkg of 12	\$2.95

- 1/8" FPT ports
- Quick removal threaded bowl



Air Filter Regulator (Brass)

Catalog No.	Description	Retail Price
VAF-527-01	Air Filter Regulator	\$120.00
VAF-527-02	Filter Element, 5 Micron with bowl seal	\$28.00
VAF-527-03	Moisture Separator Bowl	\$42.00
VAF-527-04	Mounting Bracket & Nut	\$18.00

- Wilkerson Air Filter Regulator
- 1/8" FPT ports
- Quick removal threaded bowl
- Self-relieving
- Filter with 5 micron element
- Corrosion resistant construction
- Diaphragm operated
- Must be mounted in vertical position



Air Filter Regulator (Coated Aluminum)

Catalog No.	Description	Retail Price
VAF-528-01	Air Filter Regulator- Coated Aluminum	\$100.50
VAF-528-02	Moisture Separator Bowl	\$25.50
VAF-528-03	Filter Element, 40 Micron green threads	\$20.00
VAF-528-04	Bowl Seal	\$2.90
VAF-528-05	Mounting Bracket & Nut	\$18.00

- 1/8" FPT ports
- Quick removal threaded bowl
- Self-relieving
- Filter with 40 micron element
- Diaphragm operated
- Must be mounted in vertical position



COMPONENTS

Valves / Filters / Regulators

Air Regulator

- Installs on air compressor
- Rated to 125 PSI
- 1/4" gauge port
- Increases life of utility center regulator
- 1/2" MPT
- Relieving type

Catalog No.	Description	Retail Price
VAF-529-01	Air Filter Regulator	\$170.00



Miniature Regulator

- For air or water service
- O-150 PSI operating range
- Nickel plated brass construction
- Port size, 10-32, Mounting hardware, 5/32-32

Catalog No.	Description	Retail Price
VAF-532-01	Miniature Regulator, without exhaust	\$25.00
VAF-532-02	Miniature Regulator, with exhaust	\$26.50



Miniature Regulator - Piston Type

- O-150 PSI operating range
- Nickel plated brass construction
- Piston type
- Port size 10-32
- Non-or self-relieving
- Mounting hardware 15/32-32 thread with nut

Catalog No.	Description	Retail Price
VAF-536-01	Miniature Regulator, Piston type, without exhaust	\$26.50
VAF-536-02	Miniature Regulator, Piston type, with exhaust	\$26.50



Non-Adjustable Regulator

- Supplies a maximum of 36 to 42 psi
- 15/32" -32 panel mount with nut
- Limits air pressure in water bottle systems
- Tamper proof

Catalog No.	Description	Retail Price
VAF-536-03	Non-adjustable regulator, 40 psi	\$37.00



Miniature Pressure Regulator

- To fit the following units: Cavitron 3000 and 2002, Prophy-Jet 30, Cavi-Jet 30, Cavi-Endo 25, Cavi-Med

Catalog No.	Description	Retail Price
VAF-536-04	Miniature Pressure Regulator, adjustable, non-relieving, brass, each (Dentsply #60407)	\$51.80
FIF-665-01	1/4" Poly x 10-32 Connector	\$3.20/ea



NEW

Manual Shut-Off Valve

- Used at utility center to turn off air or water during installation or service.
- Chrome finish
- Inlet, 1/2 FPT, Outlet, 3/8 O.D. tube compression fitting

Catalog No.	Description	Retail Price
VAF-535-01	Manual Shut-off Valve 1/2" FPT Inlet	\$9.50
VAF-535-02	Manual Shut-off Valve 1/2" Nominal Compression	\$9.50



Wilkerson Coalescing Filters, 5oz or 10oz

- Fits: Apollo, Matrix, Midmark, Tech West & others
- Complete assembly with indicator
- Includes auto drain
- Rated to 0.001 microns
- 3/8" NPT

Catalog No.	Description	Retail Price
VAF-537-01	5 oz Filter	\$302.80
VAF-537-02	10 oz Filter	\$366.70



Wilkerson Coalescing Filter Elements, 5oz or 10oz

- Rated to 0.01 microns
- 5 oz or 10 oz sizes
- Includes O-rings

Catalog No.	Description	Retail Price
VAF-537-03	5 oz Filter Element, 2" long	\$84.00
VAF-537-04	10 oz Filter Element, 2-1/2" long	\$90.00



Wilkerson Filter Indicator

- Informs user of clogged element
- Changes from green to red
- Includes O-rings

Catalog No.	Description	Retail Price
VAF-537-05	Filter Indicator	\$40.00



Wilkerson Indicator Blanking Plate Kit

- Replaces existing indicator
- Rated to 175 PSI
- Includes O-rings

Catalog No.	Description	Retail Price
VAF-537-06	Blanking Kit	\$16.00



Wilkerson Auto Drain

- Auto Drain, Coalescing Filter, DCI & Wilkerson
- Automatically expels condensation
- Includes O-rings

Catalog No.	Description	Retail Price
VAF-537-07	Auto Drain	\$96.00



Desiccant Chamber Purge Solenoid Valves

- Electronically purges moisture from tank
- Normally Open
- 115 Volt
- 1/4" NPT

- Used on Midmark Matrix oil-less & lubricated compressors
- Used on Apollo oil-less & lubricated compressors
- Used on Air Techniques lubricated compressors

Catalog No.	Description	Retail Price
VAF-538-01	115 Volt	\$170.00
VAF-538-02	230 Volt	\$170.00
VAF-538-03	Debris Strainer for Purge Valves	\$20.00



Water Security Systems

- Filters incoming water supply and shuts off water from control panel
- Remote shut-off of all water supply to dental equipment to prevent accidental water damage
- Allows servicing without disruption to the dental practice
- 20 Micron Filter

Catalog No.	Description	Retail Price
VAF-539-01	3/4" System with Filter and 24 Volt Relay	\$980.00
VAF-539-02	1" System with Filter and 24 Volt Relay	\$1,300.00



Water Filter Replacement Cartridge

- Fits most standard housings
- Pleated design captures more debris
- 20 Micron fiber

Catalog No.	Description	Retail Price
VAF-540-01	Water Filter 2 1/2" x 10" fits 3/4" systems	\$41.20/pkg of 2
VAF-540-02	Water Filter 4 1/2" x 10" fits 1", 1 1/4", 1 1/2" systems	\$82.40/pkg of 2



Economy Self-Contained Water Systems, Standard Or Deluxe

- Deluxe models have a selector toggle allowing to switch between city & bottled water
- 750 ml, 1 liter, & 2 liter heavy duty bottles
- Water pressure is controlled by adjustable regulator
- Standard 28mm thread for water bottle
- Gauge shows bottle pressure
- Bracket mounts to any flat vertical surface
- Optional post-mount available
- On/off toggle releases bottle pressure

Catalog No.	Description	Retail Price
ACF-542-01	Standard System (includes extra bottle), 750ml, surface mount	\$155.00
ACF-542-01D	Deluxe System (includes extra bottle), 750ml, surface mount	\$185.00
ACF-542-02	Standard System, 1-liter, surface mount	\$160.00
ACF-542-02D	Deluxe System, 1-liter, surface mount	\$190.00
ACF-542-03	Standard System, 2-liter, surface mount	\$160.00
ACF-542-03D	Deluxe System, 2-liter, surface mount	\$190.00
ACF-542-AA	Adapter to mount water system to 2" post, anodized aluminum	\$64.00
ACF-542-AW	Adapter to mount water system to 2" post, white powder coat	\$64.00



Self-Contained Water System with Gauge - Single Bottle

- Controlled dispensing system for antiseptic coolant or sterile water. May be used with handpiece controls and 3-way syringes
- Asepsis Design
- High density bottle, 1 liter
- Mini air pressure regulator
- 3-way toggle valve
- Interconnect tubing
- Mounting hardware

Catalog No.	Description	Retail Price
ACF-544-01	2" dia. Post Mount System	\$225.00
ACF-544-02	Wall Mount System	\$215.00
ACF-541-06	Replacement Pressure Head	\$28.50
ACF-541-07	Replacement Washer	\$0.60



Self-Contained Water System with Gauge - Dual Bottle

Controlled dispensing system for antiseptic coolant or sterile water. May be used with handpiece controls and 3-way syringes.

- Asepsis Design
- High density bottle, 1 liter (2)
- Mini air pressure regulator
- 3-way toggle valve
- Interconnect tubing
- Mounting hardware

Catalog No.	Description	Retail Price
ACF-545-01	2" dia. Post Mount System	\$320.00
ACF-545-02	Wall Mount System	\$280.00
ACF-545-03	Beneath Shelf Mount System	\$325.00
ACF-541-06	Replacement Pressure Head	\$28.00
ACF-541-07	Replacement Washer	\$0.60



High Density Bottle

- For use with clean water systems
- Heavy wall construction, safely operated at normal water system pressure
- High Density Polyethylene (HDPE) is compatible with coolant liquids used in dentistry
- Can be used on most dental systems with the Universal Manifold

Catalog No.	Description	Retail Price
ACF-543-01	High Density 750ml Bottle	\$7.00
ACF-543-05	High Density 1 Liter Bottle	\$7.40
ACF-543-02	High Density 2 Liter Bottle	\$12.00
ACF-543-03	Universal Manifold with barbs and washer, plastic	\$30.00
ACF-543-03L	Universal Manifold, long self centering neck, with barbs	\$25.50
ACF-543-07	Universal Manifold with barbs and washer, aluminum	\$46.00
ACF-543-06	Manifold Washer	\$3.00

Bottle dimensions:
750 ml - 3" W x 9 1/2" H
1 liter - 3 1/2" W x 11" H
2 liter - 4" W x 13 7/8" H



Quick-Switch Water Bottle System

- Easy installation requires no tools
- Saves time in changing bottles
- Self-locking feature ensures proper positioning when pressurized
- Automatic pressure release when removing bottle

Catalog No.	Description	Retail Price
ACF-543-04	Quick-Switch Water Bottle System	\$70.00
ACF-543-08	Kit includes converter and adapter	
	O-ring for Converter Manifold	\$2.90



Universal Pressure Head

- Use to build your custom fluid delivery system

Catalog No.	Description	Retail Price
ACF-543-10	Universal Pressure Head without 1/8" FPT port for gauge	\$21.70
ACF-543-11	Universal Pressure head with 1/8" FPT port for gauge	\$24.20



Slim-Line X-Ray View Panel

- Flat screen only 3/4" deep
- 12" x 8" viewing area (30 cm x 20 cm x 2 cm)
- Available with or without AC adapter
- Cold cathode lamp
- Lies flat, stands on desk or wall hung
- 110V or 220V

Catalog No.	Description	Retail Price
SLS-101-01	Slim-Line View Panel, with AC Adapter	\$240.00
SLS-101-02	Slim-Line View Panel, without AC Adapter	\$218.00
SLS-101-03	Cold cathode lamp, replacement	\$60.00



Cuspidor Bowl

- Quality porcelain.
- Bowl 8" OD x 6 3/4" ID
- Stem 2 1/2" OD x 1 1/2" ID

Catalog No.	Description	Retail Price
ACF-574-01	Cuspidor Bowl, with 2 o-rings and screen	\$96.50
ACF-574-02	Cuspidor Bowl only	\$91.50
ACF-574-03	Screen only, A-dec® style	\$5.00



Cuspidor Time Delay Kit

- Convert manual bowl rinse or cup filler to time delay
- Replace expensive and difficult to repair diaphragm systems
- Use to add a time delay to HP flush systems

Catalog No.	Description	Retail Price
ACF-575-01	Cuspidor Time Delay Kit	\$83.00



Cuspidor Cup/Bowl Spouts

Catalog No.	Description	Retail Price
ACF-574-05	Cup Filler Spout, stainless steel	\$24.00
ACF-574-06	Cuspidor Bowl Rinse Spout, stainless steel	\$28.00



Handpiece Lab System

- Requires only air supply, no water connections
- HP holder and tubing included
- Foot control, regulator and gauge pre-connected

Catalog No.	Description	Retail Price
ACF-375-01	Handpiece Lab System	\$140.00
ACF-375-02	Handpiece Lab System with gauge	\$170.00

See page E-18 for alternative style HP Lab System.



Light Reflector, A-dec Cascade 6300

Fits Cascade "6300" Lights • Reflector is NEW, not a re-coated • 4 year warranty on "optic" coating

Catalog No.	Description	Retail Price
LIT-100-01	A-dec Light Reflector, Cascade 6300	\$401.00



Light Reflector, Marus - Knight - Schein

New Glass, not a re-coat • Fits Ritter K Super Star & J. Starlite • Fits Marus, Knight, Henry Schein without Handles • No Core Charge

Catalog No.	Description	Retail Price
LIT-100-02	Light Reflector, Marus - Knight - Schein	\$550.00



Light Reflector, Pelton LF-II

New Glass, not a re-coat • Fits Pelton LF-II • No Core Charge

Catalog No.	Description	Retail Price
LIT-100-03	Light Reflector, Pelton LF-II	\$370.00



Light Reflector, Pelton LF-I

New Glass, not a re-coat • Fits Pelton LF-I • No Core Charge


Catalog No.	Description	Retail Price
LIT-100-04	Light Reflector, Pelton LF-I	\$370.00



Light Reflector, Ritter

New Glass, not a re-coat • Fits Ritter K Super Star & J. Starlite
• Fits Marus, Henry Schein with handles • No Core Charge

Catalog No.	Description	Retail Price
LIT-100-05	Light Reflector, Ritter	\$550.00

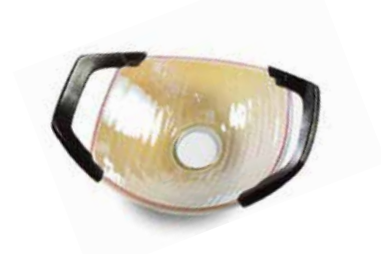

A-dec Cascade Light Shield Bracket Assembly

(A-dec #28.1012.00)

Catalog No.	Description	Retail Price
LIT-101-00	A-dec Cascade Light Shield Bracket Assembly	\$13.90



Light Shield to fit A-dec 6300

(A-dec #28.0503.01)

Fits all A-dec 6300 and 500 halogen dental lights, • Optical coating for easy cleaning and scratch resistance

Catalog No.	Description	Retail Price
LIT-101-01	Light Shield to fit A-dec 6300	\$21.00



Replacement Light Shield for Belmont/Healthcotm Lights

Fits HLC, HLU, HLT, HLW, BLC, BLU, BLT, BLW model lights with the square head (HLU2040)

• OEM part #HLU 2005 • Also fits Healthco Aristocrat™ lights

Catalog No.	Description	Retail Price
LIT-101-02	Rectangular Light Shield for Belmont/Healthco™ Lights	\$8.00



Belmont Light Shield

Belmont #8000 2005 • Fits DCI Equipment and Beaverstate lights

Catalog No.	Description	Retail Price
LIT-101-03	Belmont Light Shield	\$7.90



Belmont Clesta Light Shield

Catalog No.	Description	Retail Price
LIT-101-04	Belmont Clesta Light Shield (Ref: 501701)	\$58.00



Marus & Ritter Light Shield

Catalog No.	Description	Retail Price
LIT-101-05	Marus & Ritter Light Shield (Ritter #122845) Redesigned to include notches for the bulb housing	\$14.00



Pelton & Crane LF-II Light Shield

2 piece Set, OEM #3321598

Catalog No.	Description	Retail Price
LIT-101-06	Pelton & Crane LF-II Light Shield	\$11.40



Pelton & Crane Light Shield, LF-I

Catalog No.	Description	Retail Price
LIT-101-07	Pelton & Crane LF-I Light Shield (OEM #006777)	\$7.50



Light Bulb, 17 Volt / 95 Watt

Fits Cascade "6300" lights (A-dec #041.709.00) • Replaces 24V/150 Watt bulbs in older 6300 lights • Fits Midmark dental light #002-1154-00 • Fits Forest 9000 Series

Catalog No.	Description	Retail Price
LIT-102-00	A-dec Light bulb	\$24.00



Light Bulb, 12 VAC 55 Watt

H-3 Bulb Type, fits:

- A-dec Performer (A-dec #041.513.00)
- Ampco Contender
- Fiad Healthco
- SDS Fantastic Belmont
- Celux, Chayes Virginia, DCI
- Equipment Pro 72
- Marus Pro 70
- Belmont X-Calibur Bulb (Belmont #8000 1006)

Catalog No.	Description	Retail Price
LIT-102-01	H-3 bulb type, 12 Volt AC / 55 watt	\$12.50



Belmont Clesta Light Bulb 24V 60W Black Tip

24 Volt / 60 Watt Black Tip • Fits all Belmont Clesta Lights (Belmont 501074)

Catalog No.	Description	Retail Price
LIT-102-02	24 Volt / 60 Watt Black Tip	\$28.00



Belmont ECO Bulb with Black Tip, 12V 55 Watt

Fits: All Belmont ECO Lights • AL-711M, AL-712M, AL-715M • Belmont #G04CADN29A0

Catalog No.	Description	Retail Price
LIT-102-03	12 Volt / 55 Watt Black Tip	\$48.00



Light Bulb, Ceramic, 24 VAC 150 Watt

24 Volt AC / 150 watt

Fits: Dental Ez Distinction, Forest Faro Cool Beam, Henry Schein Knight Marus, Old Forest, Porter Daystar, Ritter K, Vacudent

Catalog No.	Description	Retail Price
LIT-102-04	Light Bulb, Ceramic, 24 VAC 150 Watt	\$18.50



Light Bulb, 24 VAC 150 Watt

Fits: Dentech, Eurodent, Faro Sunlight 70, J. Morita, Knight L LSM, Planmeca
Delight, SDS 1350, Siemens A1, E1 Sirona

Catalog No.	Description	Retail Price
LIT-102-05	Light Bulb, 24 VAC 150 Watt	\$11.50



Pelton & Crane Light Bulb, 24 VAC 150 Watt

24 Volt AC / 150 watt • Pelton & Crane LF+ / LF I / LF II / LF III

Catalog No.	Description	Retail Price
LIT-102-06	Light Bulb, 24 VAC 150 Watt	\$23.00



COMPONENTS

Lighting / Replacement Light Bulbs

Ritter Light Bulb, J Starlite, 28 VAC 75 Watt

28 Volt AC / 75 watt

Catalog No.	Description	Retail Price
LIT-102-07	Ritter Light Bulb, J Starlite, 28 VAC 75 Watt	\$32.00



MDT Light Bulb, 24 VAC 75 Watt

Catalog No.	Description	Retail Price
LIT-102-08	MDT Light Bulb, 24 VAC 75 Watt	\$51.00



DCI Equipment Light Bulb, 24 VAC 100 Watt

24 Volt AC / 100 watt • DCI Equipment light

Catalog No.	Description	Retail Price
LIT-102-09	DCI Equipment Light Bulb, 24 VAC 100 Watt	\$23.00



Replacement LED Handpiece Bulbs



Catalog No.	Description	Retail Price
HRB-110-01	Fits Kavo Turbines & Micro Motors, requires 3.5 Volt DC	\$68.20
HRB-110-02	Fits Bien-Air Micro Motors, Requires 3.5 Volt DC (Not a replacement for MX Motor Handpiece)	\$68.20
HRB-110-03	Fits Sirona Turbines & Micro Motors, Requires 3.5 Volt DC	\$68.20

Replacement Handpiece Bulbs



Catalog No.	Description	Retail Price
HRB-220-01	3.5 Volt / 750 MA, Fits Midwest XGT, Stylus & Stylus ATC HP's	\$36.00
HRB-220-02	3.5 Volt / 750 MA, Fits A-dec / W&H HP's	\$36.00
HRB-220-03	A-dec 3.5V 750 MA (A-dec #041.317.00), HPs after year 1988	\$70.00
HRB-220-04	3.3 Volt / 4 Watt, Fits Kavo Lux coupler previous PN 553-3881, current Kavo #1002.2928, Dabi Atlante	\$28.00

Light Bulb, 120 VAC 250 Watt

Fits: Chapman Huffman 60-575-00, 60-580-00, 60-585-00, 60-590-00 Dentsply Triad 2000, Kavo Vicon DLS, LS Illuminator, Polar II Lares Quartz Fire SX, Midwest Insight I and II, Tradition I Sci Can Dual Lux 2 Curing

Catalog No.	Description	Retail Price
LIT-103-00	Light Bulb, 120 VAC 250 Watt	\$30.50



Light Bulb, 12 VAC 75 Watt

Fits: 3M Visilux I Visilux II, XL3000, Almore OP II, Coltene Coltolux CuringUnits II, III and IV, Darby Dental Superdent 2001, 2005, Degussa Degulux, Demetron Optilux 101, Optilux 400, Optilux 403, Demetron VCL-100, VCL-200, VCL-300, First Medica Hilux 200, Lares Apollo CP, Lasac Tech LTI-5000, Megadental Megalux, Motloid Compolite 660, Penn Med 300L and 310 Curing Guns, Phasealloy Witelite Gun, Unitek Ortholux, Ortholux XT, Viva Dent Heliolux GT, COE Lunar TA

Catalog No.	Description	Retail Price
LIT-103-01	Light Bulb, 12 VAC 75 Watt	\$43.00



Light Bulb, 12 VAC 100 Watt

Fits: EFOS lite, Proden Executor, Henry Schein, Kinetic Instrument Sunspot, S.S. White Profile Curing Lite, Teledyne Getz Focus, L.D. Caulk Prismetics Lite

Catalog No.	Description	Retail Price
LIT-103-02	Light Bulb, 12 VAC 100 Watt	\$43.00



Light Bulb, 14 VAC 35 Watt

Fits: Healthco Aristocrat II, L.D. Caulk Maxlight, Demetron Optilux 150

Catalog No.	Description	Retail Price
LIT-103-03	Light Bulb, 14 VAC 35 Watt	\$36.00



Light Bulb, 15 VAC 150 Watt

Fits: AirTech VistaCam & VistaCam Elite, Bosworth Superlight, Fuji DnetaCam 2 & 50, Healthco VLC Illumination Wand, Deller Fiberlight Deadlight

Catalog No.	Description	Retail Price
LIT-103-04	Light Bulb, 15 VAC 150 Watt	\$22.00



Light Bulb, 12 VAC 75 Watt

Fits: Dentsply/Healthco Celebrity Perspective 1 & 2, Integra ViperCam, Medtronics Dent03, QuickCam DDS, Welsh Allyn Perspective 2

Catalog No.	Description	Retail Price
LIT-103-05	Light Bulb, 12 VAC 75 Watt	\$41.00



Light Bulb, 12 VAC 100 Watt

Fits: Acu-Cam Polo, Dentech Ultima II, Kavo Vicon LS & LSII, Micro Motors LHS 710 Fiberoptic

Catalog No.	Description	Retail Price
LIT-103-06	Light Bulb, 12 VAC 100 Watt	\$19.00



Plastic Fittings Kit

Contains 120 plastic fittings packaged in a sturdy 12 compartment clear plastic box, size 1 1/4" x 3 1/2" x 7". Kit Includes:

Qty 10 - FIF-626-01 - 10-32 x 1/16" Barb Elbows	Qty 10 - FIF-627-04 - 1/16" x 1/8" x 1/16" Tees
Qty 10 - FIF-626-02 - 10-32 x 1/8" Barb Elbows	Qty 10 - FIF-631-01 - 1/16" x 1/8" Reducing Barbs
Qty 10 - FIF-626-03 - 10-32 x 1/16" Threaded Barb Tees	Qty 10 - FIF-632-01 - 1/16" "Y" Barbs
Qty 10 - FIF-627-01 - 1/16" Barb Tees	Qty 10 - FIF-632-02 - 1/8" "Y" Barbs
Qty 10 - FIF-627-02 - 1/8" Barb Tees	Qty 10 - FIF-632-03 - 1/16" Cross Barbs
Qty 10 - FIF-627-03 - 1/8" x 1/16" x 1/8" Tees	Qty 10 - FIF-632-04 - 1/8" Cross Barbs



Catalog No.	Description	Retail Price
SIF-700-01	Plastic Fittings Kit	\$80.00

Barb Fitting Kit

Contains 8 different fittings packaged in a sturdy 8 compartment clear plastic box, 1 1/8" x 2 3/4" x 4 1/2". Kit Includes:

Qty 5 - 10-32 x 4-way cross female connector	Qty 10 - 10-32 x 1/16" barbs and washers
Qty 5 - 10-32 female connector	Qty 10 - 10-32 hex plugs and washers
Qty 5 - 10-32 male coupler	Qty 10 - 1/8" sleeve clamps
Qty 10 - 10-32 x 1/8" barbs and washers	Qty 10 - 1/16" sleeve clamps



Catalog No.	Description	Retail Price
SIF-832-01	Barb Fitting Kit	\$50.00

Technician's Fitting Kit

Kit includes:

Qty 1 - 12 compartment storage box	Qty 10 - 10-32 female tee
Qty 100 - 10-32 plastic washers	Qty 50 - 10-32 plug
Qty 10 - 1/16" in-line metal barb	Qty 50 - 10-32 1/16" metal barb
Qty 10 - 1/8" in-line metal barb	Qty 50 - 10-32 1/8" metal barb
Qty 100 - 1/8" plastic clamp sleeve	Qty 20 - 1/16" metal tee
Qty 25 - 1/4" plastic clamp sleeve	Qty 20 - 1/16" metal cross
Qty 10 - 10-32 female cross	



Catalog No.	Description	Retail Price
SIF-834-01	Technician Fitting Kit	\$221.50

Sleeves

Catalog No.	Description	Retail Price
FIF-605-01	1/8" I.D.	\$1.50/pkg of 10
FIF-605-02	1/4" I.D.	\$1.95/pkg of 10



Uni-Clamp, 1/8" I.D.

Catalog No.	Description	Retail Price
FIF-606-01	Uni-Clamp, 1/8" I.D.	\$3.20/pkg of 10
FIF-606-02	Uni-Clamp, 1/8" I.D.	\$24.00/pkg of 100



COMPONENTS

Fittings

Plastic Washer

Catalog No.	Description	Retail Price
FIF-607-01	For 10-32 Fittings	\$4.00/pkg of 100



10-32 Hex Plug Plated

Catalog No.	Description	Retail Price
FIF-611-01	10-32 Hex with washer	\$4.50/ea
FIF-611-02	10-32 Hex with washer	\$38.00/pkg of 100



MPT Socket Head

Catalog No.	Description	Retail Price
FIF-612-01	1/8" MPT Socket Head	\$1.10/ea
FIF-612-02	1/4" MPT Socket Head	\$1.50/ea
FIF-612-03	3/8" MPT Socket Head	\$2.00/ea
FIF-612-04	1/2" MPT Socket Head	\$2.00/ea



Barb Plugs

Catalog No.	Description	Retail Price
FIF-614-01	1/8" Barb Plug	\$5.70/pkg of 10
FIF-614-02	1/4" Barb Plug	\$6.00/pkg of 10



10-32 Male Coupler

Catalog No.	Description	Retail Price
FIF-615-01	10-32 Male Coupler	\$4.90/pkg of 5



8-32 x 1/16" Barbs

Catalog No.	Description	Retail Price
FIF-616-01	6-32 x 1/16" Barb	\$7.50/pkg of 10
FIF-616-02	8-32 x 1/16" Barb	\$6.50/pkg of 10



10-32 x 1/16" Barbs and 1/8" Barbs

Catalog No.	Description	Retail Price
FIF-617-01	10-32 x 1/16" Barb	\$3.60/pkg of 10
FIF-618-01	10-32 x 1/8" Barb	\$3.80/pkg of 10



Barb x 1/8" MPT Adapters

Catalog No.	Description	Retail Price
FIF-621-01	1/8" MPT x 1/16" Barb	\$2.30/ea
FIF-622-01	1/8" MPT x 1/8" Barb	\$2.30/ea
FIF-621-02	1/8" MPT x 1/16" Barb, 90 degree elbow	\$5.00/ea



Stepped Barb

Catalog No.	Description	Retail Price
FIF-623-01	10-32 x $\frac{1}{8}$ " x $\frac{1}{16}$ " Stepped Barb	\$2.40/ea



A-dec[®] $\frac{5}{16}$ " Barb Fitting

Catalog No.	Description	Retail Price
FIF-623-02	10-32 x $\frac{5}{16}$ " Barb	\$5.00/pkg of 10



In-line Tubing QD Swivel Fittings

- Facilitates installation and service
- Fits $\frac{1}{8}$ " and $\frac{1}{16}$ " tubing
- Swivel QD connection
- Mix or match tubing sizes

Catalog No.	Description	Retail Price
FIF-624-02	$\frac{1}{16}$ " Male tubing quick disconnect-swivel	\$6.00/pkg of 5
FIF-624-03	$\frac{1}{8}$ " Male tubing quick disconnect-swivel	\$7.00/pkg of 5



$\frac{1}{8}$ " & $\frac{1}{16}$ " Swivel with O-Ring

- Facilitates installation and service
- Fits $\frac{1}{8}$ " and $\frac{1}{16}$ " tubing
- Swivel QD connection
- Mix or match tubing sizes

Catalog No.	Description	Retail Price
FIF-624-04	$\frac{1}{16}$ " Swivel with O-ring	\$7.00/pkg of 5
FIF-624-05	$\frac{1}{8}$ " Swivel with O-ring	\$7.50/pkg of 5



In-line Tubing QD Collar

- Facilitates installation and service
- Fits $\frac{1}{8}$ " and $\frac{1}{16}$ " tubing
- Swivel QD connection
- Mix or match tubing sizes

Catalog No.	Description	Retail Price
FIF-624-06	Collar	\$5.40/pkg of 5



O-Ring for In-line Tubing QD Swivel Fittings

- O-Ring, Buna-n, .144 I.D. x .025 Width

Catalog No.	Description	Retail Price
FIF-624-07	O-ring	\$7.00/pkg of 12



COMPONENTS

Fittings

$\frac{1}{16}$ " In-line Tubing QD Assembly

Catalog No.	Description	Retail Price
FIF-624-10	$\frac{1}{16}$ " In-line Tubing QD Assembly	\$21.00/pkg of 5



$\frac{1}{8}$ " In-line Tubing QD Assembly

Catalog No.	Description	Retail Price
FIF-624-11	$\frac{1}{8}$ " In-line Tubing QD Assembly	\$21.00/pkg of 5



$\frac{1}{16}$ " & $\frac{1}{8}$ " In-line Tubing QD Assembly

Catalog No.	Description	Retail Price
FIF-624-12	$\frac{1}{16}$ " x $\frac{1}{8}$ " In-line Tubing QD Assembly	\$21.00/pkg of 5



In-line Needle Valve

- Poly to MPT
- 150 PSI working pressure
- Positive seal and flow adjustment

Catalog No.	Description	Retail Price
FIF-624-20	In-line Needle Valve, $\frac{1}{4}$ " x $\frac{1}{4}$ "	\$26.00



Right Angle Needle Valves

- Poly to MPT
- 150 PSI working pressure
- Positive seal and flow adjustment

Catalog No.	Description	Retail Price
FIF-624-21	Right Angle Needle Valve, $\frac{1}{4}$ " x $\frac{1}{8}$ "	\$26.00
FIF-624-22	Right Angle Needle Valve, $\frac{1}{4}$ " x $\frac{1}{4}$ "	\$26.00



Self Drilling Saddle Valve

- Used to tap into copper lines
- Self drills into $\frac{3}{8}$ " to $1\frac{3}{8}$ " copper tubing
- Outlet is $\frac{1}{4}$ " O.D. compression tube fitting

Catalog No.	Description	Retail Price
FIF-624-30	Self Drilling Saddle Valve, $\frac{1}{4}$ " x $\frac{1}{8}$ "	\$14.00



In-line Miniature QD with Shut-Off

- Valve stem for shut-off

Catalog No.	Description	Retail Price
FIF-624-40	In-line Miniature QD with shut-off, $\frac{1}{16}$ " x 10-32	\$19.50
FIF-624-41	In-line Miniature QD with shut-off, $\frac{1}{8}$ " x 10-32	\$17.00



In-line Barb Connectors

Catalog No.	Description	Retail Price
FIF-625-01	$\frac{1}{16}$ " Barb Connector	\$5.30/pkg of 10
FIF-625-02	$\frac{1}{8}$ " Barb Connector	\$8.20/pkg of 10
FIF-625-03	$\frac{1}{8}$ " x $\frac{1}{16}$ " Barb Connector	\$10.00/pkg of 10



Manifold Block w/10-32 Ports

- 9 each 10-32 ports for barbs
- Reduces tees and crosses
- Helps manage tubings

Catalog No.	Description	Retail Price
FIF-625-04	Distribution Manifold, 9 female 10-32 ports	\$30.00



Barb Elbow Connectors (Plastic)

Catalog No.	Description	Retail Price
FIF-626-01	10-32 x $\frac{1}{16}$ " Barb Elbow	\$9.30/pkg of 10
FIF-626-02	10-32 x $\frac{1}{8}$ " Barb Elbow	\$10.30/pkg of 10



Threaded Plastic Tees (Plastic)

Catalog No.	Description	Retail Price
FIF-626-03	10-32 x $\frac{1}{16}$ " Barb Tee	\$12.40/pkg of 10
FIF-626-04	10-32 x $\frac{1}{8}$ " Barb Tee	\$13.40/pkg of 10



Barb Tee Connectors (Plastic)

Catalog No.	Description	Retail Price
FIF-627-01	$\frac{1}{16}$ " Barb Tee	\$8.25/pkg of 10
FIF-627-02	$\frac{1}{8}$ " Barb Tee	\$9.30/pkg of 10
FIF-627-03	$\frac{1}{8}$ " x $\frac{1}{16}$ " x $\frac{1}{8}$ " Tee	\$9.30/pkg of 10
FIF-627-04	$\frac{1}{16}$ " x $\frac{1}{8}$ " x $\frac{1}{16}$ " Tee	\$9.30/pkg of 10



Barb Connectors (Metal)

Catalog No.	Description	Retail Price
FIF-628-01	$\frac{1}{16}$ " Barb Tee	\$4.50/ea
FIF-628-02	$\frac{1}{8}$ " Barb Tee	\$4.60/ea
FIF-628-03	$\frac{1}{16}$ " Cross	\$7.00/ea



Duck Bill Barb "Check Valve"

- For handpiece and syringe water lines. Used to provide water line suck back control
- Different appearance from standard barb
- Will fit any port where standard barb is used

Catalog No.	Description	Retail Price
FIF-629-01	Duck Bill Barb	\$5.50/ea



Barb x 10-32 "T" Fitting

- For 1/4" tubing

Catalog No.	Description	Retail Price
FIF-630-01	1/8" Barb x 10-32 Port	\$14.00/pkg of 5



Straight Barb Connectors (Plastic)

Catalog No.	Description	Retail Price
FIF-631-01	1/16" x 1/8" Reducing Barb	\$8.25/pkg of 10
FIF-631-02	1/16" Straight Through Barb	\$8.25/pkg of 10
FIF-631-03	1/8" Straight Through Barb	\$8.25/pkg of 10



"Y" Barb Connectors (Plastic)

Catalog No.	Description	Retail Price
FIF-632-01	1/16", "Y" Barb	\$10.30/pkg of 10
FIF-632-02	1/8", "Y" Barb	\$10.30/pkg of 10



Cross Barb Connectors (Plastic)

Catalog No.	Description	Retail Price
FIF-632-03	1/16" Barb Cross	\$12.40/pkg of 10
FIF-632-04	1/8" Barb Cross	\$12.40/pkg of 10



Press Barbs (Brass)

- Use where thread clearance is low.


Catalog No.	Description	Retail Price
FIF-633-01	1/16" Brass Barb	\$1.50/ea
FIF-633-02	1/8" Brass Barb	\$1.00/ea



Hookup Tubing Quick-Connects (Plastic)

- Simple twist-disconnect fittings

Catalog No.	Description	Retail Price
FIF-634-01	Male luer for 1/8" tubing	\$5.20/pkg of 10
FIF-634-02	Male luer for 1/4" tubing	\$5.20/pkg of 10
FIF-634-03	Female luer for 1/8" tubing	\$5.20/pkg of 10
FIF-634-04	Female luer for 1/4" tubing	\$5.20/pkg of 10

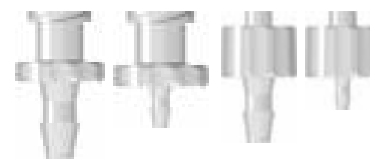

10-32 Female Coupler (Plated)

Catalog No.	Description	Retail Price
FIF-635-01	10-32 Female Coupler	\$6.00/pkg of 5



10-32 Cross / Tee Fitting (Plated)

Catalog No.	Description	Retail Price
FIF-636-01	10-32, 4-Way Female Cross	\$1.60/ea
FIF-636-02	10-32, Female Tee	\$1.60/ea



MPT x 10-32 Cross (Plated)

- One x 1/8" MPT Port
- Three #10-32 Female Ports
- Convenient Manifold for Utility Center

Catalog No.	Description	Retail Price
FIF-636-05	MPT x 10-32 Cross	\$3.20



10-32 1/2" Cross / Tee Fitting (Plated)

Catalog No.	Description	Retail Price
FIF-638-01	10-32, 4-Way Female Cross, 1/2" body	\$1.60/ea
FIF-638-02	10-32, Female Tee, 1/2" body	\$1.60/ea



Female Cross (Brass)

Catalog No.	Description	Retail Price
FIF-637-01	1/8 FPT, 4-Way Cross	\$8.00/ea



Pipe Tee Female (Brass)

Catalog No.	Description	Retail Price	Catalog No.	Description	Retail Price
FIF-639-01	1/8" FPT	\$4.70/ea	FIF-639-03	3/8" FPT	\$8.50/ea
FIF-639-02	1/4" FPT	\$4.70/ea	FIF-639-04	1/2" FPT	11.00/ea



Elbow 90 Degree Female (Brass)

Catalog No.	Description	Retail Price	Catalog No.	Description	Retail Price
FIF-641-01	1/8" FPT Elbow	\$5.00/ea	FIF-641-03	3/8" FPT Elbow	\$6.50/ea
FIF-641-02	1/4" FPT Elbow	\$7.50/ea	FIF-641-04	1/2" FPT Elbow	\$8.00/ea



Street Elbow Male x Female (Brass)

Catalog No.	Description	Retail Price	Catalog No.	Description	Retail Price
FIF-642-01	1/8" Elbow	\$3.20/ea	FIF-642-03	3/8" Elbow	\$4.40/ea
FIF-642-02	1/4" Elbow	\$3.90/ea	FIF-642-04	1/2" Elbow	\$5.20/ea



Cross Connector (Brass)

Catalog No.	Description	Retail Price
FIF-643-01	1/8" MPT x 10-32 Cross	\$3.30/ea



Pipe Tee / Elbow Connector (Brass)

Catalog No.	Description	Retail Price
FIF-644-01	1/8" MPT x 10-32 Elbow	\$6.00/ea
FIF-644-02	1/8" MPT x 10-32 Tee	\$6.00/ea



COMPONENTS

Fittings

10-32 Swivel Elbow (Plated)

Catalog No.	Description	Retail Price
FIF-645-01	10-32 Male Swivel x 10-32 Female	\$6.80/ea



1/8" FPT x 10-32 Adapter (Plated)

Catalog No.	Description	Retail Price
FIF-646-01	1/8" FPT x 10-32 Adapter	\$8.30/ea



Cross Connectors (Plated)

Catalog No.	Description	Retail Price
FIF-646-10	10-32 Male, 10-32 Female (3)	\$4.50/ea



Compression Barb (Plated)

• Fits 3/8" Compression fitting,
Use to connect 1/4" tube

Catalog No.	Description	Retail Price
FIF-646-15	Compression Barb	\$3.60/ea



Reducing Coupler (Brass)

Catalog No.	Description	Retail Price	Catalog No.	Description	Retail Price
FIF-647-01	1/4" x 1/8" Reducer	\$3.00/ea	FIF-647-03	1/2" x 1/4"	\$7.60/ea
FIF-647-02	3/8" x 1/4"	\$4.00/ea	FIF-647-04	1/2" x 3/8"	\$7.60/ea



Reducer Bushing (Brass)

Catalog No.	Description	Retail Price
FIF-648-01	1/8" MPT x 10-32	\$2.60/ea



NEW

Female Coupler (Brass)

Catalog No.	Description	Retail Price	Catalog No.	Description	Retail Price
FIF-649-01	1/8" x 1/8" Coupler	\$1.35/ea	FIF-649-03	3/8" x 3/8" Coupler	\$2.50/ea
FIF-649-02	1/4" x 1/4" Coupler	\$2.00/ea	FIF-649-04	1/2" x 1/2" Coupler	\$4.10/ea



NEW

MPT Hex Plug (Brass)

Catalog No.	Description	Retail Price	Catalog No.	Description	Retail Price
FIF-650-01	1/8" MPT Hex Plug	\$1.00	FIF-650-03	3/8" MPT Hex Plug	\$1.90
FIF-650-02	1/4" MPT Hex Plug	\$1.30			



NEW

Bushing MPT x FPT (Brass)

Catalog No.	Description	Retail Price	Catalog No.	Description	Retail Price
FIF-651-01	1/4" MPT x 1/8" Bushing	\$1.50/ea	FIF-651-02	1/2" MPT x 1/4" Bushing	\$3.20/ea
FIF-651-03	3/8" MPT x 1/8" Bushing	\$2.00/ea	FIF-651-06	1/2" MPT x 3/8" Bushing	\$2.70/ea
FIF-651-04	3/8" MPT x 1/4" Bushing	\$2.00/ea	FIF-651-07	3/4" MPT x 3/8" Bushing	\$8.00/ea
FIF-651-05	1/2" MPT x 1/8" Bushing	\$3.80/ea	FIF-651-08	3/4" MPT x 1/4" Bushing	\$12.00/ea
			FIF-651-09	3/4" MPT x 1/2" Bushing	\$4.80/ea



Female x Compression Tube Adapters (Plated)

• With barbs and washers.

Catalog No.	Description	Retail Price
FIF-652-01	10-32 to $\frac{3}{8}$ " Compression Nut	\$5.00/ea
FIF-652-02	10-32 to $\frac{1}{4}$ " Compression Nut	\$4.50/ea



10-32 Female x Compression Tube (Brass)

Catalog No.	Description	Retail Price
FIF-652-03	10-32 Female x $\frac{3}{8}$ "	\$3.00/ea
FIF-652-04	10-32 Female x $\frac{1}{4}$ "	\$3.40/ea



Compression Tube x MPT Adapters (Brass)

Catalog No.	Description	Retail Price
FIF-653-01	$\frac{1}{2}$ " MPT to $\frac{3}{8}$ " Compression	\$8.20/ea
FIF-653-02	$\frac{1}{8}$ " MPT to $\frac{3}{8}$ " Compression	\$3.80/ea
FIF-653-03	$\frac{1}{4}$ " MPT to $\frac{3}{8}$ " Compression	\$8.20/ea



MPT Close Nipple (Brass)

Catalog No.	Description	Retail Price
FIF-654-01	$\frac{1}{8}$ " MPT Close Nipple	\$0.80/ea
FIF-654-02	$\frac{1}{4}$ " MPT Close Nipple	\$1.20/ea
FIF-654-03	$\frac{3}{8}$ " MPT Close Nipple	\$1.70/ea
FIF-654-04	$\frac{1}{2}$ " MPT Close Nipple	\$2.00/ea



MPT Long Nipple (Brass)

Catalog No.	Description	Retail Price
FIF-655-01	$\frac{1}{8}$ " MPT x $1\frac{1}{2}$ " Nipple	\$2.00/ea
FIF-655-02	$\frac{1}{4}$ " MPT x $1\frac{1}{2}$ " Nipple	\$1.70/ea
FIF-655-03	$\frac{3}{8}$ " MPT x $1\frac{1}{2}$ " Nipple	\$3.20/ea
FIF-655-04	$\frac{1}{2}$ " MPT x 2" Nipple	\$6.88/ea



Hex Nipple MPT (Brass)

Catalog No.	Description	Retail Price
FIF-657-01	$\frac{1}{8}$ " MPT	\$1.35
FIF-657-02	$\frac{1}{4}$ " MPT	\$2.05
FIF-657-03	$\frac{3}{8}$ " MPT	\$2.90



MPT Reducing Nipple (Brass)

Catalog No.	Description	Retail Price
FIF-658-01	$\frac{1}{4}$ " MPT x $\frac{1}{8}$ " MPT	\$2.10
FIF-658-02	$\frac{3}{8}$ " MPT x $\frac{1}{4}$ " MPT	\$3.50



Poly Union (Brass)

Catalog No.	Description	Retail Price
FIF-662-01	$\frac{1}{4}$ " x $\frac{1}{4}$ " Poly Union	\$2.50/ea
FIF-662-02	$\frac{3}{8}$ " x $\frac{3}{8}$ " Poly Union	\$3.00/ea
FIF-662-04	$\frac{5}{16}$ " x $\frac{5}{16}$ " Poly Union	\$4.10/ea



COMPONENTS

Fittings

Poly x MPT Fittings (Brass)

Catalog No.	Description	Retail Price
FIF-663-01	1/4" Poly x 1/8" MPT	\$2.00/ea
FIF-663-02	1/4" Poly x 1/4" MPT	\$2.70/ea
FIF-663-03	3/8" Poly x 1/8" MPT	\$2.55/ea
FIF-663-04	3/8" Poly x 1/4" MPT	\$3.00/ea
FIF-663-05	5/16" Poly x 1/8" MPT	\$2.80/ea
FIF-663-06	5/16" Poly x 1/4" MPT	\$3.80/ea
FIF-663-07	3/8" Poly x 3/8" MPT	\$6.70/ea
FIF-663-08	1/4" poly x 3/8" MPT	\$7.40/ea



Poly Reducing Union (Brass)

Catalog No.	Description	Retail Price
FIF-664-01	3/8" x 1/4" Union	\$5.10



NEW

Poly x Compression Tube Unions (Brass)

Catalog No.	Description	Retail Price
FIF-664-02	1/4" x 1/4"	\$7.90
FIF-664-03	3/8" x 3/8"	\$8.00



NEW

Poly x 10-32 Male Fitting (Brass)

Catalog No.	Description	Retail Price
FIF-665-01	1/4" Poly x 10-32 Connector	\$3.20/ea



Poly Tees (Brass)

Catalog No.	Description	Retail Price
FIF-666-01	1/4" Poly Tee	\$6.00/ea
FIF-666-02	3/8" Poly Tee	\$6.80/ea
FIF-666-03	3/8" x 3/8" x 1/4" Poly Tee	\$10.00/ea
FIF-666-04	5/16" Poly Tee	\$9.00/ea



Poly x MPT Elbow Fittings (Brass)

Catalog No.	Description	Retail Price
FIF-667-01	1/8" Poly x 1/8" MPT Elbow	\$2.70/ea
FIF-667-02	1/4" Poly x 1/8" MPT Elbow	\$3.00/ea
FIF-667-03	1/4" Poly x 1/4" MPT Elbow	\$3.10/ea
FIF-667-04	1/4" Poly x 3/8" MPT Elbow	\$7.50/ea
FIF-667-05	5/16" Poly x 1/8" MPT Elbow	\$4.90/ea
FIF-667-06	3/8" Poly x 1/8" MPT Elbow	\$5.00/ea
FIF-667-07	3/8" Poly x 1/4" MPT Elbow	\$6.00/ea



Poly x Long MPT Elbow Fitting (Brass)

- For water supply to vacuum pump impeller housing
- When using on vacuum pump, replace poly nut with metal sleeve

Catalog No.	Description	Retail Price
FIF-667-10	1/4" Poly X 1/4" MPT	\$20.00/ea



NEW

Poly Run Tee x MPT (Brass)

Catalog No.	Description	Retail Price
FIF-668-01	1/4" Poly x 1/8" MPT	\$6.20/ea
FIF-668-02	1/4" Poly x 1/4" MPT	\$6.20/ea



Male Poly Tee (Brass)

Catalog No.	Description	Retail Price
FIF-668-03	1/4" x 1/4" Poly Run x 1/8" MPT	\$6.75/ea
FIF-668-04	5/16" x 5/16" Poly Run x 1/8" MPT	\$15.50/ea



Union Tee Compression (Brass)

Catalog No.	Description	Retail Price
FIF-669-01	1/4" Union Tee	\$5.80/ea
FIF-669-02	5/16" Union Tee	\$5.80/ea
FIF-669-03	3/8" Union Tee	\$6.00/ea
FIF-669-04	1/2" Union Tee	\$18.00/ea
FIF-669-05	5/8" Union Tee	\$17.00/ea



Union Compression (Brass)

Catalog No.	Description	Retail Price
FIF-670-01	1/8" Union	\$2.50/ea
FIF-670-02	1/4" Union	\$3.50/ea
FIF-670-03	5/16" Union	\$3.00/ea
FIF-670-04	3/8" Union	\$3.80/ea



Compression Tube x MPT Straight (Brass)

Catalog No.	Description	Retail Price
FIF-670-10	1/4" tube x 1/8" MPT straight	\$2.20/ea
FIF-670-11	3/8" tube x 1/8" MPT straight	\$2.90/ea
FIF-670-12	1/4" tube x 1/4" MPT straight	\$2.70/ea
FIF-670-13	3/8" tube x 1/4" MPT straight	\$2.80/ea
FIF-670-14	3/8" tube x 3/8" MPT straight	\$3.20/ea
FIF-670-15	3/8" tube x 1/2" MPT straight	\$4.80/ea



Compression Tube x MPT Elbow (Brass)

Catalog No.	Description	Retail Price
FIF-670-20	1/4" tube x 1/8" MPT elbow	\$3.50/ea
FIF-670-21	1/4" tube x 1/4" MPT elbow	\$4.00/ea
FIF-670-22	3/8" tube x 1/8" MPT elbow	\$5.50/ea
FIF-670-23	3/8" tube x 1/4" MPT elbow	\$5.50/ea



Compression Nut (Brass)

Catalog No.	Description	Retail Price
FIF-670-30	1/8" Compression nut, pkg of 10	\$3.60/pkg of 10
FIF-670-31	1/4" Compression nut, pkg of 10	\$4.00/pkg of 10
FIF-670-32	5/16" Compression nut, pkg of 10	\$9.00/pkg of 10
FIF-670-33	3/8" Compression nut, pkg of 10	\$9.00/pkg of 10



Poly x FPT Straight Connectors (Brass)

Catalog No.	Description	Retail Price
FIF-671-01	1/4" Poly x 1/8" FPT	\$1.80/ea
FIF-671-02	1/4" Poly x 1/4" FPT	\$3.30/ea
FIF-671-03	3/8" Poly x 1/4" FPT	\$3.30/ea
FIF-671-04	5/16" Poly x 1/8" FPT	\$7.50/ea



COMPONENTS

Fittings

Poly x Compression Tube (Brass)

Catalog No.	Description	Retail Price
FIF-672-01	1/4" Poly to 3/8" Compression	\$5.40/ea
FIF-672-02	3/8" Poly to 3/8" Compression	\$5.50/ea



Push To Connect Fittings, Swivel

Catalog No.	Description	Retail Price
FIF-672-10	Swivel Elbow 1/8" NPT x 1/4" Quick Connect	\$8.50/ea
FIF-672-12	Swivel Elbow 1/4" NPT x 1/4" Quick Connect	\$19.50/ea



Push To Connect Fitting, Straight

Catalog No.	Description	Retail Price
FIF-672-11	Straight Fitting 1/8" NPT x 1/4" Quick Connect	\$6.00/ea



Poly Nut & Sleeve (Brass)

Catalog No.	Description	Retail Price
FIF-673-01	1/4" Poly Nut with Poly Sleeves	\$4.50/pkg of 5
FIF-673-02	3/8" Poly Nut with Poly Sleeves	\$5.00/pkg of 5
FIF-673-03	5/16" Poly Nut with Poly Sleeves	\$4.50/pkg of 5



Poly Nut with Brass Sleeve (Brass)

- Recommended for long term water connection

Catalog No.	Description	Retail Price
FIF-673-04	1/4" Poly Nut with Brass Sleeve	\$1.30/ea



NEW

Poly Sleeves (Plastic)

Catalog No.	Description	Retail Price
FIF-674-01	1/4" Poly Sleeve	\$4.55/pkg of 10
FIF-674-02	3/8" Poly Sleeve	\$4.55/pkg of 10
FIF-674-03	5/16" Poly Sleeve	\$4.55/pkg of 10



Compression Sleeve, (Delrin)

Catalog No.	Description	Retail Price
FIF-674-04	1/4" delrin compression sleeve	\$0.50/ea
FIF-674-05	5/16" delrin compression sleeve	\$0.60/ea
FIF-674-06	3/8" delrin compression sleeve	\$0.60/ea
FIF-674-07	1/2" delrin compression sleeve	\$0.60/ea



NEW

Poly Plugs (Plastic)

Catalog No.	Description	Retail Price
FIF-675-01	1/4" Poly Plug	\$3.00/pkg of 5
FIF-675-02	3/8" Poly Plug	\$3.00/pkg of 5
FIF-675-03	5/16" Poly Plug	\$5.50/pkg of 5



Brass Sleeves (Brass)

Catalog No.	Description	Retail Price	FIF-676-03	5/16" Brass Sleeve	\$2.00/pkg of 10
FIF-676-01	1/8" Brass Sleeve	\$2.30/pkg of 10	FIF-676-04	3/8" Brass Sleeve	\$4.10/pkg of 10
FIF-676-02	1/4" Brass Sleeve	\$2.10/pkg of 10	FIF-676-05	1/2" Brass Sleeve	\$4.20/pkg of 10



Male Connector Flare x MPT (Brass)

Catalog No.	Description	Retail Price	FIF-677-03	3/8" x 1/4"	\$2.65/ea
FIF-677-01	1/4" x 1/4"	\$1.50/ea	FIF-677-04	1/2" x 1/4"	\$4.35/ea
FIF-677-02	5/16" x 1/4"	\$3.50/ea			



Union Tee Flare (Brass)

Catalog No.	Description	Retail Price
FIF-678-01	3/8" Union Tee	\$4.20



Male Elbow Flare x MPT (Brass)

Catalog No.	Description	Retail Price
FIF-679-01	1/4" x 1/4"	\$2.30
FIF-679-02	1/4" x 1/8"	\$3.10
FIF-679-03	3/8" x 1/4"	\$2.30



Flare Nut (Brass)

Catalog No.	Description	Retail Price
FIF-679-05	1/4" Flare Nut	\$0.80
FIF-679-06	3/8" Flare Nut	\$1.00



Mini Quick Disconnect (Plated)

- 2-piece threaded assembly

Catalog No.	Description	Retail Price
FIF-680-01	1/16" Mini Quick Disconnect	\$6.20/ea
FIF-680-02	1/8" Mini Quick Disconnect	\$6.20/ea



Poly Male Insert QD without Shut-Off (Plated)

- For use with air only

Catalog No.	Description	Retail Price
FIF-681-01	1/4" OD Male x 1/4" Poly without Shut-Off	\$8.00/ea
FIF-681-02	3/8" OD Male x 3/8" Poly without Shut-Off	\$10.50/ea



Poly Male Insert QD with Shut-Off (Plated)

- For use with water or air

Catalog No.	Description	Retail Price
FIF-682-01	1/4" OD Male x 1/4" Poly with Shut-Off	\$11.00/ea
FIF-682-02	3/8" OD Male x 3/8" Poly with Shut-Off	\$21.00/ea



COMPONENTS

Fittings

Male Insert QD, Rear Port (Plated)

- Fits all 1/4" Female QD's
- 10-32 Female Threaded Rear Port

Catalog No.	Description	Retail Price
FIF-682-03	1/4" QD with 10-32 Rear Port	\$14.40/ea



Male QD with Shut-Off (Plated)



NEW

- 1/4" Male Quick Disconnect with shut-off

Catalog No.	Description	Retail Price
FIF-682-04	Male QD Valve with 10-32F end	\$15.50/ea
FIF-682-05	Male QD Valve with 10-32M end	\$15.50/ea
FIF-682-06	Male QD Valve with 1/16" barb end	\$15.50/ea
FIF-682-07	Male QD Valve with 3/32" barb end	\$15.50/ea
FIF-682-08	Male QD Valve with 1/8" barb end	\$15.50/ea



Male QD with Needle Valve (Plated)

- 1/4" Male QD with Needle Valve for Flow Adjustment
- 10-32 Female

Catalog No.	Description	Retail Price
FIF-682-09	1/4" QD with Needle Valve	\$15.50/ea



NEW

Poly QD Female Panel Mount with Shut-Off (Plated)

- Each, includes both yellow and blue indicators.

Catalog No.	Description	Retail Price
FIF-683-01	1/4" ID Female x 1/4" Poly with Shut-Off	\$13.20/ea
FIF-683-02	3/8" ID Female x 3/8" Poly with Shut-Off	\$24.00/ea



QD Female x MPT with Shut-Off (Plated)

Catalog No.	Description	Retail Price
FIF-684-01	1/4" ID Female x 1/8" MPT with Shut-Off	\$7.00/ea
FIF-684-02	1/4" ID Female x 1/4" MPT with Shut-Off	\$8.50/ea
FIF-684-03	3/8" ID Female x 1/4" MPT with Shut-Off	\$9.70/ea



Panel Mount x 10-32 with Shut-Off (Plated)

- Includes 1/16" and 1/8" barb

Catalog No.	Description	Retail Price
FIF-685-01	1/4" ID Coupler with Shut-Off	\$16.00/ea



CPC Type Panel Mount Female QD with Shut-Off (Plated)

Catalog No.	Description	Retail Price
FIF-686-01	1/4" Coupler with Shut-Off (use with FIF-687-01 Poly Fitting)	\$43.00/ea



CPC Type Panel Mount Male QD with Shut-Off (Plated)

- For use with air only, poly fitting

Catalog No.	Description	Retail Price
FIF-687-01	1/4" Quick-Disconnect Insert, with air only (use with FIF-686-01 Poly Fitting)	\$17.00/ea



Air Hose Quick Connects Male (Plated)

Catalog No.	Description	Retail Price
FIF-687-10	M Type Male Thread 1/4"	\$10.00/ea



Air Hose Quick Connects Female (Plated)

Catalog No.	Description	Retail Price
FIF-687-11	M Type Male Thread 1/4"	\$16.00/ea



Quick Disconnect Barb Reducers

- For use on QD when installing Dentsply CaviJet & Cavitron. Single Piece Design

Catalog No.	Description	Retail Price
FIF-687-20	1/16" x 1/4" Compression QD Barb Reducers	\$5.00/ea
FIF-687-21	1/16" x 3/8" Compression QD Barb Reducers	\$5.00/ea
FIF-687-22	1/8" x 1/4" Compression QD Barb Reducers	\$5.00/ea
FIF-687-23	1/8" x 3/8" Compression QD Barb Reducers	\$4.20/ea



Tracer Washers

- Use on panel mount valves
- Yellow for air, Blue for water
- Fits toggles, flow controls etc.
- Ideal for female quick-disconnects
- Non-corroding plastic

Catalog No.	Description	Retail Price
FIF-687-02	Blue Water Indicator Washer for 1/4" QD	\$3.10/pkg of 10
FIF-687-03	Blue Water Indicator Washer for 3/8" QD	\$3.10/pkg of 10
FIF-687-04	Yellow Air Indicator Washer for 1/4" QD	\$3.10/pkg of 10
FIF-687-05	Yellow Air Indicator Washer for 3/8" QD	\$3.10/pkg of 10



QD Manifold Assembly

Catalog No.	Description	Retail Price
FIF-688-01	6-port Quick Disconnect Manifold Assembly	\$40.00
FIF-688-02	8-port Quick Disconnect Manifold Assembly	\$50.00
FIF-688-03	6-port Gasket	\$2.50
FIF-688-04	8-port Gasket	\$3.50



Plastic Vacuum Adapters (Plastic)

Catalog No.	Description	Retail Price			
FIF-689-01	1/2" x 1/2"	\$2.00/ea	FIF-689-04	5/8" x 5/8"	\$2.00/ea
FIF-689-02	1/2" x 5/8"	\$2.00/ea	FIF-689-05	5/8" x 7/8"	\$8.00/ea
FIF-689-03	1/2" x 7/8"	\$6.30/ea	FIF-689-06	3/4" x 3/4"	\$2.50/ea



Plastic Straight Male Adapters (Plastic)

Catalog No.	Description	Retail Price			
FIF-691-01	1/2" MPT x 1/2" Barb	\$1.40/ea	FIF-691-05	1/2" MPT x 3/4" Barb	\$2.10/ea
FIF-691-02	3/4" MPT x 1/2" Barb	\$2.00/ea	FIF-691-06	3/4" MPT x 3/4" Barb	\$2.10/ea
FIF-691-03	1/2" MPT x 5/8" Barb	\$2.30/ea	FIF-691-07	3/4" MPT x 1" Barb	\$2.70/ea
FIF-691-04	3/4" MPT x 5/8" Barb	\$2.10/ea	FIF-691-08	1" MPT x 1" Barb	\$2.50/ea



Plastic Elbow Male Adapters (Plastic)

Catalog No.	Description	Retail Price			
FIF-692-01	1/2" MPT x 1/2" Barb	\$2.10/ea	FIF-692-05	1/2" MPT x 3/4" Barb	\$2.60/ea
FIF-692-02	3/4" MPT x 1/2" Barb	\$2.60/ea	FIF-692-06	3/4" MPT x 3/4" Barb	\$3.10/ea
FIF-692-03	1/2" MPT x 5/8" Barb	\$1.80/ea	FIF-692-07	3/4" MPT x 1" Barb	\$3.10/ea
FIF-692-04	3/4" MPT x 5/8" Barb	\$2.60/ea	FIF-692-08	1" MPT x 1" Barb	\$3.10/ea



Plastic Elbow Female Adapters (Plastic)

Catalog No.	Description	Retail Price
FIF-693-01	1/2" FPT x 1/2" Barb	\$2.30/ea
FIF-693-02	1/2" FPT x 5/8" Barb	\$3.50/ea



Plastic Barbed "T" Fittings (Plastic)

Catalog No.	Description	Retail Price
FIF-694-01	1/2" x 1/2" x 1/2" Barbed Tee	\$1.60/ea
FIF-694-02	5/8" x 5/8" x 5/8" Barbed Tee	\$1.80/ea
FIF-694-03	5/8" x 5/8" x 1/2" Barbed Tee	\$1.80/ea



Plastic Barb x Barb Elbow (Plastic)

Catalog No.	Description	Retail Price
FIF-695-01	1/2" x 1/2" Barbed Elbow	\$2.30/ea
FIF-695-02	5/8" x 5/8" Barbed Elbow	\$1.80/ea



Plastic Barbed "Y" Fitting (Plastic)

Catalog No.	Description	Retail Price
FIF-696-01	1/2" x 1/2" x 1/2" Barbed Y, each	\$4.10/ea



Adjustable Tie

Catalog No.	Description	Retail Price
FIF-601-01	4" length	\$1.00/pkg of 10
FIF-601-02	5" length	\$1.10/pkg of 10
FIF-601-03	8" length	\$3.20/pkg of 10



Adjustable Tie Mount

Catalog No.	Description	Retail Price
FIF-602-01	Adjustable Tie Mount, adhesive	\$4.60/pkg of 10



Heyco Snap-In Bushings

Catalog No.	Description	Retail Price
FIF-603-08	1/2" Black	\$1.00/pkg of 5
FIF-603-09	3/4" Black	\$1.50/pkg of 5
FIF-603-10	1 1/2" Black	\$3.50/pkg of 5



Cable Clamp

Catalog No.	Description	Retail Price
FIF-604-01	1/4"	\$2.20/pkg of 10
FIF-604-02	1/2"	\$3.30/pkg of 10
FIF-604-03	3/4"	\$4.40/pkg of 10
FIF-604-04	1 1/4"	\$8.60/pkg of 10



Stainless Steel Hose Clamps

Stainless Steel Hose Clamp. 1/2" band width.

Catalog No.	Description	Retail Price
SIF-608-01	3/8" to 7/8"	\$2.15
SIF-608-02	1/2" to 1-1/4"	\$2.30
SIF-608-03	3/4" to 1-3/4"	\$2.40
SIF-608-04	1-1/4" to 2-1/4"	\$2.50
SIF-608-05	1-3/4" to 2-3/4"	\$2.60
SIF-608-06	2-1/4" to 3-1/4"	\$2.70



Post Strap (Plastic)

- Plastic strap fits 2" post
- 2-3/4" hole distance.

Catalog No.	Description	Retail Price
FIF-609-01	Plastic hub strap	\$3.10



Sleeve Installation Tools

Used for installation of tubing sleeve clamps.

Catalog No.	Description	Retail Price
SIF-801-01	Standard 1/8" Sleeve Tool	\$20.00
SIF-801-02	Standard 1/4" Sleeve Tool	\$20.00
SIF-801-03	Standard 1/8" x 1/4" Multi-Sleeve Tool	\$20.00



Service Tool

Adapts to 4-hole Midwest handpiece tubing to provide instant air source for cleaning of plugged tubings or testing of components.

- Nickel plated brass construction
- Universal barb fits 1/8" or 1/4" QD tubing

Catalog No.	Description	Retail Price
SIF-802-01	Service Tool	\$25.00
FIF-623-01	Replacement Barb	\$4.40



Economy Test Gauge

Replaceable barb for use with 1/8" or 1/16" ID tubing. 0-160 PSI gauge.

Catalog No.	Description	Retail Price
SIF-811-01	Test Gauge	\$26.00
FIF-623-01	Replacement Barb	\$4.40



Blow Gun

- 125 PSI Working Pressure
- Lever Operated
- 1/4 FPT Connection
- Available with 10' Coiled Tubing (yellow) with 1/4 MPT Fittings

Catalog No.	Description	Retail Price
SIF-813-01	Blow Gun with 10ft Coil Tubing	\$47.00
SIF-813-02	Blow Gun only	\$19.00
SIF-813-03	10ft Coil Tubing only	\$30.00
SIF-813-04	Blow Gun Tip Kit	\$29.00



O-Ring Lubricant

- For synthetic and rubber based o-rings
- Silicon base lubricant
- Moisture resistant.

Catalog No.	Description	Retail Price
SIF-814-01	O-Ring Lubricant, 2 oz.	\$32.00
SIF-814-02	Silicone Lubricant, net wt. .04 oz, 5/pkg	\$7.00



Loctite Service Items

- For bonding almost any materials
- For sealing threaded connections

Catalog No.	Description	Retail Price
SIF-815-02	Thread Locker 242, 0.34 oz	\$33.00
SIF-815-03	Super Bonder 409, 0.10 oz	\$11.00



All Purpose Lube

- Moisture resistant, Teflon base lubricant.
- Not recommended for O-rings

Catalog No.	Description	Retail Price
SIF-816-01	All Purpose Lube, 1.5 oz.	\$16.50



O-Ring Tools

- Aiding in installing o-rings onto stems
- 3 sizes: .199, .166, .100


Catalog No.	Description	Retail Price
SIF-817-01	O-Ring Tools Kit (A-dec Ref: #009.013.00)	\$24.00



Hex Wrench Sets

- Long arm series
- Heat treated alloy steel, black finish
- Available in standard and metric sizes
- Ball driver end

Catalog No.	Description	Retail Price
SIF-821-01	Hex Wrench Set, standard	\$27.00
SIF-821-02	Hex Wrench Set, metric	\$34.00

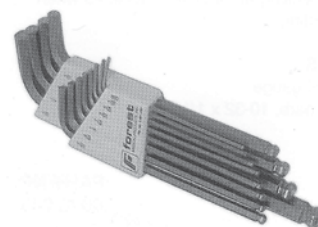

O-Ring Service Kit

- Contains 120 o-rings in 12 popular sizes from 7/64" O.D. to 1 1/2" O.D.
- Packaged in a sturdy 12 compartment clear plastic box, 1 1/4" x 3 1/2" x 7".

Catalog No.	Description	Retail Price
SIF-823-01	O-Ring Kit	\$45.00



Dentist's Emergency Repair Kit

Use for repair of minor problems occurring most frequently in dental units. Kit Includes:

Qty 10 – 1/4" Plastic Sleeve	Qty 2 – HP Gasket, Midwest 5-Hole
Qty 10 – 1/8" Plastic Sleeve	Qty 6 – O-Ring HVE Top
Qty 5 – 1/4" Poly Sleeve	Qty 6 – O-Ring HVE /Barrel
Qty 2 – 1/4" Poly Nut & Sleeve	Qty 6 – O-Ring Saliva Ejector Barrel
Qty 20 – 10-32 Plastic Washer	Qty 3 – O-Ring Adec Style Syringes
Qty 5 – 10-32 x 1/8" Plated Barb	Qty 3 – O-Ring Std Syringe Button
Qty 5 – 10-32 x 1/16" Plated Barb	Qty 3 – O-Ring Syringe Tip
Qty 3 – 1/16" Plated Inline Barb	Qty 3 – O-Ring Syringe Button
Qty 3 – 1/8" Plated Inline Barb	Qty 1 – Tubing Stripper Tool
Qty 2 – 1/16" Plastic Tee	Qty 10 – 4" Cable Tie
Qty 2 – 1/8" Plastic Tee	Qty 1 – Saliva Ejector Tip, Gray

Catalog No.	Description	Retail Price
SIF-830-01	Dentist's Repair Kit	\$65.00



Diaphragm Sheet Stock

11" x 8.5" Sheet, .013" Thick

Catalog No.	Description	Retail Price
SIF-844-04	Diaphragm Sheet, each	\$25.00









Post Leveling Tool

Used for quick and easy leveling of light and unit posts. Vertical & horizontal.

Catalog No.	Description	Retail Price
SIF-860-02	Magnetic, attaches to posts	\$12.00



O-Rings

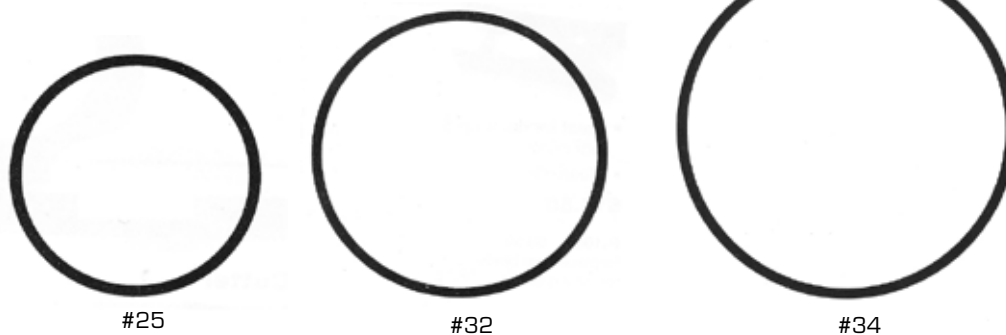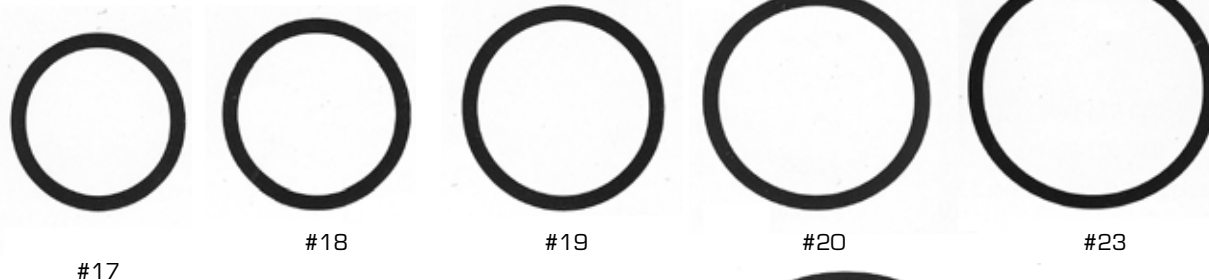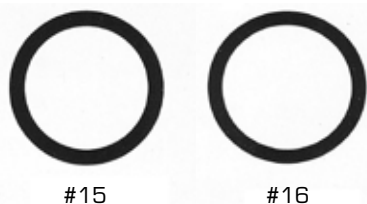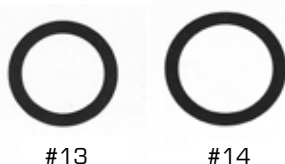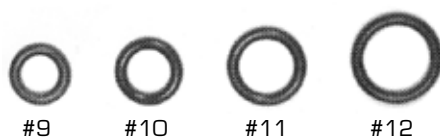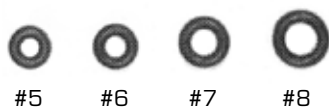

Package of 12, (unless otherwise noted).

Catalog No.	Description	
SIF-831-01	#1	.029 ID, .109 OD, .040 Wall
SIF-831-02	#2	.042 ID, .142 OD, .050 Wall
SIF-831-03	#3	.056 ID, .176 OD, .060 Wall
SIF-831-04	#4	.070 ID, .210 OD, .070 Wall
SIF-831-05	#5	.101 ID, .241 OD, .070 Wall
SIF-831-06	#6	.114 ID, .254 OD, .070 Wall
SIF-831-07	#7	.145 ID, .285 OD, .070 Wall
SIF-831-08	#8	.176 ID, .316 OD, .070 Wall
SIF-831-09	#9	.208 ID, .348 OD, .070 Wall
SIF-831-10	#10	.239 ID, .379 OD, .070 Wall
SIF-831-11	#11	.301 ID, .441 OD, .070 Wall
SIF-831-12	#12	.645 ID, .504 OD, .070 Wall
SIF-831-13	#13	.426 ID, .566 OD, .070 Wall
SIF-831-14	#14	.489 ID, .629 OD, .070 Wall
SIF-831-15	#15	.551 ID, .070 Wall
SIF-831-16	#16	.614 ID, .754 OD, .070 Wall
SIF-831-17	#17	.676 ID, .816 OD, .070 Wall
SIF-831-18	#18	.739 ID, .883 OD, .070 Wall
SIF-831-19	#19	.801 ID, .941 OD, .070 Wall
SIF-831-20	#20	.864 ID, 1.004 OD, .070 Wall
SIF-831-23	#23	1.051 ID, .070 Wall
SIF-831-25	#25	1.1875 ID, 1.3125 OD, .0625 Wall, pkg of 10
SIF-831-32	#32	1.864 ID, 2.004 OD, .070 Wall, pkg of 10
SIF-831-34	#34	2.114 ID, .070 Wall

Retail Price: \$2.90 per package

HANDPIECES

Fiber Optic High Speed Handpiece

• Sable Access and Rotamax Handpieces are manufactured to meet and exceed ISO 7785-1 Handpiece Standard. All Sable Handpieces include Push Button chuck operation and a Quick Connect, 360 degree swivel feature. All products are supplied with a comprehensive manufacturers 12-Month warranty.



**Quality
Handpieces
Made in Canada**



Catalog No.	Description	Retail Price
2000015	Access S4500KL Mini Head Fiber Optic Highspeed Handpiece (No Coupler)	\$670
2000016	Rotamax S6500KL Torque Head Fiber Optic Highspeed (No Coupler)	\$670

Kavo Type Couplers

2000114	4-Hole Kavo MULTIflex Non-Fiber Optic Coupler	\$131
2000105	5-Hole Kavo MULTIflex Lux Type Fiber Optic Coupler	\$165
2000116	6-Hole Kavo MULTIflex Lux Type Fiber Optic Coupler	\$237
2000119	6-Hole Kavo MULTIflex Lux Type LED Fiber Optic Coupler	\$190

NSK Type Couplers

2000107-6	6-Hole NSK NL400 Fiber Optic Coupler	\$167
2000107-4	4-Hole NSK NL400 Non Optic Coupler	\$151
2000107-6L	6-Hole NSK NL400 LED Bulb Fiber Optic Coupler	\$217

Non-Fiber Optic Highspeed Handpieces



Catalog No.	Description	Retail Price
2000005	Sable Access II, PB, Miniature Head, 360/Quick Connect, Coupler Included	\$423
2000003	Sable Access II, PB, Miniature Head, Handpiece only, 360/Quick Connect, Coupler Required	\$372
2000006	Sable Rotamax II, PB, Standard Head, 360/Quick Connect, Coupler Included	\$423
2000004	Sable Rotamax II, PB, Standard Head Handpiece only, 360/Quick Connect, Coupler Required	\$372
2000101	4-Hole 360/Quick Connect Coupler, To Fit Access II & Rotamax II	\$51
2000102	2-Hole 360/Quick Connect Coupler, To Fit Access II & Rotamax II	\$51

COMPONENTS

Handpieces / HP Accessories / HP Connectors

• Sable Premium Style Slowspeed Products carry a 12-Month Warranty on Complete Handpieces from date of install

Sable Premium Slowspeed Handieces



Catalog No.	Description	Retail Price
1600008	Sable NSK Style - Doctor's Package, includes: 20K Air Motor; 1:1 Contra Angle, 1:1 Straight Attachment, Standard Swing Latch Head	\$448
1600009	Sable NSK Style - Doctor's Package, includes: 20K Air Motor; 1:1 Contra Angle, 1:1 Straight Attachment, Ball Bearing PB Latch Head	\$503
1600010	Sable Kavo Style - Doctor's Package, includes: 20K Air Motor; 1:1 Contra Angle, 1:1 Straight Attachment, Ball Bearing Latch Head	\$1,090

Sable Hygentist Handpieces



Catalog No.	Description	Retail Price
1600011	Hygienist E-Type Prophy Package 1, includes 20K Air Motor, 4:1 Reduction Contra Angle, Snap-on Prophy Head	\$358
1600012	Hygienist E-Type Prophy Package 2, includes 20K Air Motor, 4:1 Reduction Contra Angle, Screw-in Prophy Head	\$358
1600013	Hygienist E-Type Prophy Package 3, includes 20K Air Motor, 4:1 Reduction Contra Angle, Standard Latch Head	\$358
1600014	Hygienist E-Type Prophy Package 4, includes 20K Air Motor, 4:1 Reduction Contra Angle, Ball Bearing Latch Head	\$420
1600015	RDH Style Hygiene Prophy Handpiece	\$595

Sable Premium Slowspeed Motor & Attachments



Catalog No.	Description	Retail Price
1600110	20K Air Motor, with Internal Air & Water, 20,000 RPM, Accepts all ISO E-type Attachments	\$409
1600210	1:1 Ball Bearing Contra Angle, with Internal Air & Water, 40,000 RPM, Accepts all Kavo & Sable type Heads	\$200
1600310	1:1 Ball Bearing Straight Nosecone, with Internal Air & Water, 40,000 RPM	\$298
1600810	Ball Bearing Push Button Latch Head	\$184
1601201	Lubrication Nozzle for Premium Slowspeed Handpieces	\$14

6-Pin Fiber Optic Assembly

- Fully assembled with Connector & Stainless Steel Nut
- Highly flexible PVC or Silcryn™ tubing with copper wire, assembled with 6-pin connector, including sleeve, electrical contacts and barbs.

Catalog No.	Description	Retail Price
4E-82347-xx-7'	4-hole PVC tubing w/ 2 wires & 6-pin connector, 7' length	\$80.00
4E-82347-xx-7'S	4-hole Silcryn tubing w/ 2 wires & 6-pin connector, 7' length	\$120.00
4E-82347-xx-10'	4-hole PVC tubing w/2 wires & 6-pin connector, 10' length	\$90.00
4E-82347-xx-10'S	4-hole Silcryn tubing w/ 2 wires & 6-pin connector, 10' length	\$130.00
OP-823-00	Euro / 6-pin connector assembly (SS nut, sleeve, electric contacts & barb)	\$36.00

Color Codes: 01-black, 14-light gray, 54-sterling, 72-sand, 92-surf



8-Pin Micromotor Tubing Assembly

- 8-Pin Tubing (compatible with BienAir MC3® micromotor)
- For micromotors with and without light
- Short, rotating coupling with stainless steel nut & plastic connector
- Silicone jacket tubing

Catalog No.	Description	Retail Price
4V-050-xx	6' (1.8 m), straight - 14 only (light gray)	\$90.00



Silicone Pneumatic Tubing - Euro Style *made in Italy by TeKne*

* Color Codes:
 Black (01)
 Light Gray (14)
 Brown (02)
 White (11)



- Pneumatic straight silicone tubing for turbines, air micro motors and scalers.
- The silicone sheath and the connectors can be easily and safely cleaned with any commercial cleansers.
- Nuts and sleeves are chrome-plated and specially satined to reduce surface glare.
- Sizes of the internal polyurethane tubings have been optimized to get the best performance of the attached instrument and minimize the kinking effect of sharp bends.
- Internal tubings are extruded in controlled environment and meet FDA standards
- On special order, the inner waterline tubing can be extruded with an antimicrobial additive which prevents proliferation of bacteria on the inner surface
- Different types of silicone tubings are available upon request

Catalog No.	Description	Retail Price
4V-060-xx	Silicone Pneumatic Tubing with 3-hole connector & sleeve, 5-1/2" length, price/each	\$90.00
4V-061-xx	Silicone Pneumatic Tubing with 4-hole connector & sleeve, 5-1/2" length, price/each	\$90.00
4V-062-xx	Silicone Pneumatic Tubing, 6-pin connector & sleeve, 5-1/2" length, price/each	\$125.00
4V-063-xx	Silicone Pneumatic Tubing with 5-hole connector & sleeve and integrated bulb, 5-1/2" length, price/each	\$225.00

COMPONENTS

Handpieces / HP Accessories / HP Connectors

Silicone Dental Tubing for Micromotors (NSK or Bien-Air)



Straight silicone tubing for electric micromotors.

- The silicone sheath can be easily and safely cleaned with any commercial cleansers.
- Inner tubings are extruded in controlled environment without any added phthalates, meeting FDA and Directive 2002/72/EC & 2007/47/EC standards.
- Nuts and sleeves are chrome-plated with satin finish to reduce uncomfortable surface glare.
- All electric contacts are gold-plated and precisely crimped to guarantee the best electric connection.
- All electric wires are extra-flexible, adequately

- sized to passing current and insulated with a special oil-resistant material.
- Standard length of silicone tubing: 5 feet (1.60 m).
- Available in 4 standard colors: black, gray, white or brown.

Options:

- On special order, the inner waterline tubing can be extruded with an antimicrobial additive which prevents proliferation of bacteria on the inner surface.
- Different types of silicone hoses are available on special request.



made in Italy by TeKne

Catalog No.	Description	Retail Price
4V-080-xx	Silicone tubing w/4-hole connector & 4-pin electric socket for use w/NSK M40 (N, LED) micromotor	\$110.00
4V-090-xx	Silicone tubing w/4-hole connector & 4-pin electric socket for use w/DEFINITIVE or BienAir MC3 LK/LED micromotors	\$136.00
4V-100-xx	Silicone tubing w/4-hole connector & 4-pin electric socket for use w/BienAir MC2 IR micromotor	\$100.00
4V-110-xx	Silicone tubing w/4-hole connector & 4-pin electric socket for use w/BienAir MC2 ISOLITE 300/LED micromotor	\$125.00

Power-Optic 5-Hole Handpiece Tubing Assembly

Includes brass plated nut and connector, with or without bulb

Catalog No.	Description	Retail Price
4V-070-xx-07	Power-Optic 5-hole HP tubing, 7ft	\$164.00
4V-070-xx-07B	Power-Optic 5-hole HP tubing, 7ft, with 3.5V halogen bulb	\$211.00
4V-070-xx-12	Power-Optic 5-hole HP tubing, 12ft	\$190.00
4V-070-xx-12B	Power-Optic 5-hole HP tubing, 12ft, with 3.5V halogen bulb	\$237.00

Colors (xx): lt. gray (14), sand (72), sterling (54), surf (92)



ISO-C 6-Pin Economy HP Illumination System, 5'

- For system where lamp is in handpiece or swivel adapter
- Transformer, 120V to 3.85 VDC 1 AMP
- Air electric switch, 10 PSI activated
- Installation kit with instructions
- 4-hole asepsis handpiece tubing
- Available in Gray (G), Light Sand (LS), Sterling (S), Black (B), Dark Surf (DS)
- ISO-C 6-pin, 5ft standard length, straight tubing (also available in coiled)

Catalog No.	Description	Retail Price
HPA-100-01	ISO-C 6-pin Economy HP Illumination System	\$210.00



ISO 5-Hole Economy HP Illumination System, 5'

- Lamp in handpiece tubing connector
- Transformer, 120V to 3.85 VDC 1 AMP
- Air electric switch, 10 PSI activated
- Installation kit with instructions
- Lamp module, ISO, 3.5 VAC, 720 MA
- ISO 4-hole asepsis handpiece tubing
- Available in Gray (G), Lt. Sand (LS), Sterling (S), Black (B), Dark Surf (DS)
- 5 ft. standard length straight tubing (available in coiled tubing)

Catalog No.	Description	Retail Price
HPA-101-01	ISO 5-hole Economy HP Illumination System	\$265.00



5-Hole Lamp Module

- For use with ISO fiber optic tubing
- 3.5 VDC, 720 MA
- Red gasket
- Ultra high pressure provides 50% more light & longer bulb life

Catalog No.	Description	Retail Price
HPA-105-02	5-pin Lamp Module	\$68.00



Dual Switch Fiber Optic Light Control Set with Power Transformer

- Dual Switch - LED & Halogen For 1 or 2 Optic Handpieces
- Operates all brands of handpieces
- Power Source: AC/DC 24V
- Adjustable Output: DC 2V-3.5V
- Adjustable Illumination time delay: 0.5 - 4.5 second for each handpiece
- Weight: 4 oz
- Size: W 2.25" x L 3.5" x H 1"
- Universal Power Transformer:
Input: 100 to 240VAC, 50/60 Hertz 0.5 Amps,
Output: 12VDC, 1.0 Amps

Catalog No.	Description	Retail Price
HPA-105-05	Dual Fiber Optic Light Control Set	\$310.00



Scaler Tubing, 25K, Connector

- 25Khz 3-Wire Handpiece Assembly with connector
- Fits 2002 Cavitron, 25Khz Cavijet, & 25 Khz Sonatron
- Handle length is 6"

Catalog No.	Description	Retail Price
HPA-107-01	Scaler Tubing 25k with Connector	\$340.00



Scaler Tubing, 30K, Connector

- 30Khz 3-Wire Handpiece Assembly with connector
- For models without prophylaxis function
- Handle length is 5-1/8"

Catalog No.	Description	Retail Price
HPA-107-02	Scaler Tubing 30k with Connector	\$340.00



Terminal Nut, Short

- Metal construction, threaded
- For use on fiber optic tubing using no lamp module
- Converts ISO to ISO-C tubing

Catalog No.	Description	Retail Price
HPA-106-01	Terminal Nut, Short	\$46.50



Terminal Nut, Long

- Metal construction, threaded
- For use with ISO fiber optic tubing with lamp module #HPA-105-02
- Converts ISO-C to ISO style tubing

Catalog No.	Description	Retail Price
HPA-106-02	Terminal Nut, Long	\$46.50



COMPONENTS

Handpieces / HP Accessories / HP Connectors

Midwest 4-Hole Connector with Stainless Steel Nut

Catalog No.	Description	Retail Price
TB-103	SS Nut & Metal Connector	\$12.80
TB-102	SS Nut & Plastic Connector	\$11.80
TB-100	Metal Connector Only	\$8.20
TB-101	SS Nut Only	\$4.60
TB-104	Plastic Connector Only	\$7.30



Borden 3-Hole Connector with Stainless Steel Nut

Catalog No.	Description	Retail Price
TB-203	SS Nut & Metal Connector	\$12.80
TB-202	SS Nut & Plastic Connector	\$11.80
TB-200	Metal Connector Only	\$8.20
TB-201	SS Nut Only	\$4.60
TB-204	Plastic Connector Only	\$7.30



Midwest 4-Hole Connector with Nickel Plated Nut

Catalog No.	Description	Retail Price
OP-560-00	4-hole Plated Nut with Plastic Connector	\$10.50
OP-570-00	4-hole Plated Nut with Plated Connector	\$11.00
OP-567-00N	4-hole Plated Nut only	\$4.60



Borden 3-Hole Connector with Nickel Plated Nut

Catalog No.	Description	Retail Price
OP-550-00	3-hole Plated Nut with Plastic Connector	\$9.00
OP-550-00N	3-hole Plated Nut only	\$4.60



Autoclavable Handpiece Gasket (325°F)

Catalog No.	Description	Retail Price
SIF-793-01	Borden, 2/3 hole	\$2.00/ea
SIF-793-02	Midwest, 4-hole	\$2.00/ea
SIF-793-03	Midwest, 5-hole	\$2.00/ea
SIF-793-04	6-pin ISO	\$3.50/ea



4-Hole Tubing Connector

Catalog No.	Description	Retail Price
SIF-794-01	4-hole HP Tubing Connector	\$9.00



Handpiece Tubing Adapter

Catalog No.	Description	Retail Price
SIF-795-01	Install 2/3-Hole HP to 4-Hole HP Connector	\$80.00
SIF-795-02	Install 4-Hole HP to 2-Hole HP Connector	\$80.00
SIF-795-03	Install 4-Hole HP to 3-Hole HP Connector	\$90.00



Handpiece Purge Station with Hepa Filter

- Purges excess lubricant and water from handpiece
- Patented Bio-Cleanse HEPA exhaust filter
- Hepa Filter helps protect operator from air-borne pathogens
- Provides for quick and easy handpiece maintenance
- Exhaust reservoir safely collects excess lubricant
- Purges handpiece coolant air and water lines
- Easy view gauge with adjustable air regulator
- Air supply line with 1/4" male quick-disconnect included
- Compact size of a foot control - requires little counter space
- Handpiece adapter required - see below

Catalog No.	Description	Retail Price
ACF-566-01	Handpiece Purge Station with Hepa Filter	\$204.00
ACF-566-02	Single Replacement Hepa Filter	\$8.30

Adapters:

ACF-566-10	Midwest® Adapter	\$39.00	ACF-566-14	Borden Adapter	\$58.50
ACF-566-11	Star® Adapter	\$41.00	ACF-566-15	Midwest XGT Adapter	\$52.00
ACF-566-12	Kavo® Adapter	\$52.00	ACF-566-16	Lares Adapter	\$52.00
ACF-566-13	NSK, Kinetics, Healthco, Viper, Atech Adapter	\$58.50	ACF-566-17	W&H® Adapter	\$52.00



Bottle Upgrade for HP Purge Station

- Use as a replacement for most brands of HP Purge Stations
- Upgrades all units to accept Bio-Cleanse Hepa Filter System

Catalog No.	Description	Retail Price
ACF-566-03	Bottle with milled cap (no filter)	\$12.40
ACF-566-02	Single Replacement Hepa Filter	\$8.30



Handpiece Muffler and Oil Collector

- Compact size to fit inside control head
- Can be mounted thru frame for removal from outside control head
- Will handle up to four handpieces
- Minimum back pressure
- Filter media is standard cotton swab material

Catalog No.	Description	Retail Price
ACF-567-01	Handpiece Muffler and Oil Collector	\$28.00
ACF-567-02	Exhaust Muffler	\$22.00



Air Actuated Electric Switch

- Air activated electric switch, 4 PSI cracking pressure
- 10-32 port - 1/16" exchangeable fitting
- Normally open or closed
- U.L. listed switch

Catalog No.	Description	Retail Price
ACF-564-01	Normally Open	\$20.00
ACF-564-02	Normally Closed	\$20.00



Handpiece Pressure Gauge

- Used to test accuracy of drive air pressure with handpiece running
- Connects in-line for accurate measurements
- Range indicator from 0-60 PSI & 1-4 bar
- Fits 6-pin and 5-hole Midwest connectors or 2/3-hole Borden connectors

Catalog No.	Description	Retail Price
SIF-803-01	Handpiece Pressure Gauge, Midwest	\$26.00
SIF-803-02	Handpiece Pressure Gauge, Borden	\$90.00



DENTAL EQUIPMENT	7 - 37
DENTAL COMPONENTS	41 - 125
DENTAL TUBING	129 - 139
NUMERICAL INDEX PAGES	141



Chart / Material Comparison	129
Silcryn™ Tubing	131
Handpiece Tubing	132
Syringe Tubing	135
Foot Control Tubing	137
Vacuum & Saliva Ejector Tubing	139
Supply Tubing	141
Tubing Supplies	143

We proudly offer premium quality Freelin-Wade dental tubing, a leading American manufacturer of quality dental and industrial tubing since 1981. In addition to the wide range of stock tubing; customizing in size, color and combination is available to meet your specific needs. The same premium quality tubing is used in production of our dental units and handpiece assemblies. The next pages feature in-stock dental tubing readily available. Information on the extensive range of custom dental and industrial tubing is available upon request.

Standard Color Codes

 01 Black	 02 Brown	 03 Silver	 04 Gray
 05 Red	 06 Green	 07 Blue	 08 Yellow
 09 Orange	 10 Natural	 11 White	 12 Purple
 14 Light Gray	 25 Transparent Red	 26 Transparent Green	 27 Transparent Blue
 28 Transparent Yellow	 29 Transparent Orange	 32 Transparent Purple	 37 Ice Blue
 45 Neon Pink	 46 Neon Green	 48 Neon Yellow	 49 Neon Orange
 54 Sterling	 72 Sand	 81 Off White	 92 Surf

The colors in this chart are not realistic representations of actual tubing colors.
For better color evaluation, please contact us for a sample.

Authorized Distributor For




Copyright Notice: Freelin-Wade is a registered trademark of Freelin-Wade, Silcryn™ is a trademark of Freelin-Wade.

Tubing Material Comparison

	PVC (Polyvinyl Chloride)	PUR (Polyurethane)	Silcryn™
Cost	\$	\$\$	\$\$\$
Look & Feel	<ul style="list-style-type: none"> · Softer (combined with flexible material) · Smooth & Shiny 	<ul style="list-style-type: none"> · Less Flexible than Silcryn™ and PVC · Smooth & Shiny 	<ul style="list-style-type: none"> · Matte-like, Soft & Silky – similar to silicone · Hangs flat & straight
Qualities	<ul style="list-style-type: none"> · Easy to work with · Tear & kink resistant · Cost effective 	<ul style="list-style-type: none"> · Long lasting · Abrasion Resistant (ex. dentist roll over it with chairs) · Works better (does not stiffen) with use of lubricants · More heat resistant than PVC 	<ul style="list-style-type: none"> · High tear resistance · Better color uniformity than silicone
Disadvantages	<ul style="list-style-type: none"> · Does not work well with lubricants 	<ul style="list-style-type: none"> · Costs more than PVC 	<ul style="list-style-type: none"> · Costs more than PUR & PVC · Does not work well with lubricants
Ideal Uses	<ul style="list-style-type: none"> · Handpieces, syringe & foot control tubing, supply tubing · Most used tubing 	<ul style="list-style-type: none"> · Handpieces, syringe & foot control tubing 	<ul style="list-style-type: none"> · Scalers, handpieces, umbilicals, vacuums & syringe tubing. · Jacket for multi-tube & wire assemblies

Factory Recommendations:

- Avoid exposure to direct sunlight or heat. Ideal working room temperature is 75°F or 20°C.
- Avoid working pressure over 80 psi
- Avoid constant stool roll-overs which cause weak points.
- Check dental unit regulators periodically for proper function.
- Ideally all tubing should be replaced every 3 years.

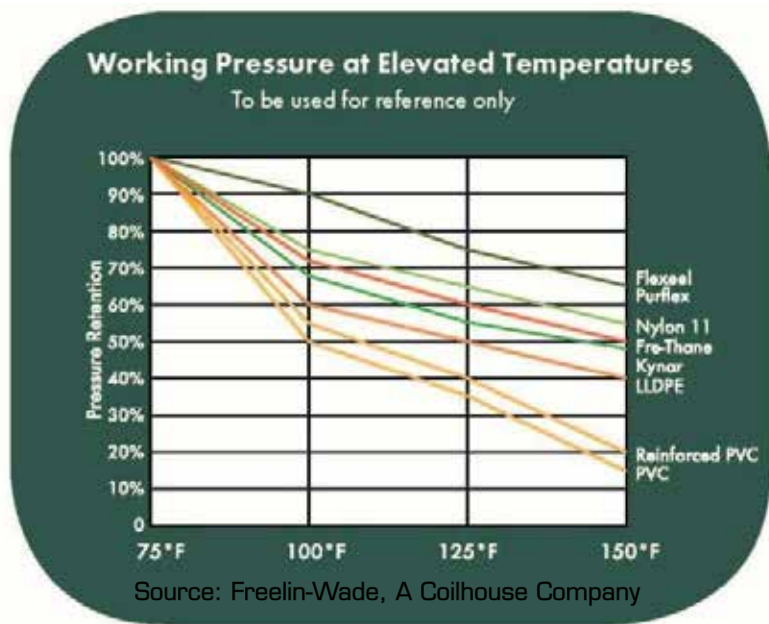

Note: Fre-Thane = Polyurethane (PUR)

*Silcryn is a trademark of Freelin-Wade, A Coilhouse Company

Need
Customized
Tubing?
Send Us
Specs

4-Hole Asepsis Handpiece Tubing, 70A

Catalog No.	Package/Unit Size	Color	Retail Price
4G-549-__	100' Box	01, 14, 54, 72, 92	\$0.899/ft
4A-549-__	500' Reel		\$0.846/ft
4S-549-__	7' length, assembled with metal adapter and nut		\$22.88/ea



2-Hole Syringe Tubing, 70A

Catalog No.	Package/Unit Size	Color	Retail Price
2B-549-__	250' Reel	14, 54	\$0.676/ft
2A-549-__	500' Reel		\$0.621/ft



Saliva Ejector Tubing, 70A

Catalog No.	Package/Unit Size	I.D.	Color	Retail Price
1G-502-__	100' box	3/16"		\$0.468/ft
1A-502-__	500' Reel	3/16"	14, 54, 72, 92	\$0.462/ft
1G-585-__	100' box	7/32"		\$0.468/ft
1A-585-__	500' Reel	7/32"	14, 54, 92	\$0.462/ft



6-Sided Vacuum Tubing, 70A

Catalog No.	Package/Unit Size	Color	Retail Price
1A-516-__	200' Reel	14, 54, 72, 92	\$1.961/ft



Scaler 70A Jacket with Tube & 3 Wires

Catalog No.	Package/Unit Size	Color	Retail Price
4E-533-__	7' Straight	01, 14	\$20.06/ea



Euro-HP Jacket, .420" x .350", 70A

Catalog No.	Package/Unit Size	Color	Retail Price
1A-534-__	500' Reel	14, 54 Custom Colors Available	\$1.142/ft



2-Hole Standard (Polyurethane)

Catalog No.	Package/Unit Size	Color	Retail Price
2G-001-__	100' Reel	01, 14, 54, 72, 92	\$0.469/ft
2A-001-__	1000' Reel		\$0.445/ft
2V-001-__	7' Coil		\$8.29/ea



2-Large Hole Standard (Polyurethane)

Catalog No.	Package/Unit Size	Color	Retail Price
2A-005-01	500' Reel	01	\$0.832/ft
2V-005-01	7' Coil		\$13.74/ea



2-Hole Asepsis (Polyurethane)

Catalog No.	Package/Unit Size	Color	Retail Price
2G-041-__	100' Reel	01, 14, 54, 72	\$0.664/ft
2A-041-__	1000' Reel		\$0.613/ft
2V-041-__	7' Coil		\$8.86/ea



2-Hole Asepsis (PVC)

Catalog No.	Package/Unit Size	Color	Retail Price
2G-047-__	100' Reel	14	\$0.559/ft
2A-047-__	1000' Reel	14, 54	\$0.492/ft



3-Hole Asepsis (Polyurethane)

Catalog No.	Package/Unit Size	Color	Retail Price
3G-042-__	100' Reel	01, 14, 54, 72, 92	\$0.716/ft
3A-042-__	1000' Reel		\$0.665/ft
3V-042-__	7' Coil		\$10.15/ea
3S-042-__	7' Coiled, assembled with metal adapter & nut		\$27.04/ea



3-Hole Asepsis (PVC)

Catalog No.	Package/Unit Size	Color	Retail Price
3G-047-__	100' Reel	01, 14, 54, 72, 92	\$0.519/ft
3A-047-__	1000' Reel		\$0.464/ft
3S-047-__	7' Straight, assembled with metal adapter & nut		\$22.88/ea



3-Hole Standard (PVC)

Catalog No.	Package/Unit Size	Color	Retail Price
3G-001-__	100' Reel	01, 14	\$0.519/ft
3A-001-__	1000' Reel		\$0.463/ft
3S-001-__	7' Straight, assembled with metal adapter & nut		\$21.84/ea



4-Hole Asepsis (Polyurethane)

Catalog No.	Package/Unit Size	Color	Retail Price
4G-042-__	100' Reel	01, 14, 54, 72, 92	\$0.870/ft
4A-042-__	500' Reel		\$0.817/ft
4V-042-__	7' Coil		\$11.05/ea
4S-042-__	7' Coiled, assembled with metal adapter & nut		\$24.44/ea



4-Hole Asepsis (PVC)

Catalog No.	Package/Unit Size	Color	Retail Price
4G-047-__	100' Reel	01, 14, 54, 72, 92	\$0.573/ft
4A-047-__	500' Reel		\$0.519/ft

*available with wires and/or Euro Connector Assembly

4-Hole Standard (PVC)

Catalog No.	Package/Unit Size	Color	Retail Price
4G-001-__	100' Reel	01, 14	\$0.582/ft
4A-001-__	500' Reel		\$0.523/ft
4S-001-__	7' Straight, assembled with metal adapter & nut		\$19.34/ea



2-Hole Sheathed Syringe Tubing

Catalog No.	Package/Unit Size	Color	Retail Price
2S-005-__	7' straight tubing with 5' sheath	14, 54, 72	\$13.42/ea



2-Hole (Polyurethane)

Catalog No.	Package/Unit Size	Color	Retail Price
2G-006-__	100' Reel	01, 14, 54, 72, 92	\$0.305/ft
2B-006-__	1000' Reel		\$0.295/ft
2A-006-__	2000' Reel		\$0.269/ft
2V-006-__	7' Coil		\$9.08/ea




2-Hole Asepsis (Polyurethane)

Catalog No.	Package/Unit Size	Color	Retail Price
2G-046-__	100' Reel	01, 14, 54, 72, 92	\$0.320/ft
2B-046-__	1000' Reel		\$0.302/ft
2A-046-__	2000' Reel		\$0.295/ft
2V-046-__	7' Coil		\$7.30/ea



3-Hole (Polyurethane)

Catalog No.	Package/Unit Size	Color	Retail Price
3G-013-__	100' Reel	14	\$0.550/ft
3B-013-__	500' Reel		\$0.518/ft
3A-013-__	1000' Reel		\$0.498/ft
3V-013-__	7' Coil		\$9.22/ea

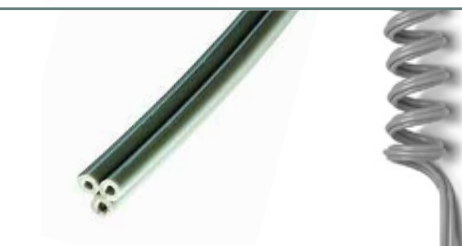

3-Hole Asepsis (Polyurethane)

Catalog No.	Package/Unit Size	Color	Retail Price
3G-413-__	100' Reel	14	\$0.639/ft
3B-413-__	500' Reel		\$0.617/ft
3A-413-__	1000' Reel		\$0.589/ft
3V-413-__	7' Coil		\$8.80/ea



3-Hole Sheathed Syringe Tubing

Catalog No.	Description	Color	Retail Price
3S-005-xx	3-Hole Sheathed Syringe Tubing, price per 7 ft length	14, 92	\$8.74



5-Hole (Polyurethane)

- Jacket tube - .030 PVC light gray color only
- Two tubes - .240 OD x .150 ID PUR 85A, colors: orange and yellow
- Three tubes - .125 OD x .066 ID PUR 85A, colors: brown, clear, green

Catalog No.	Package/Unit Size	Color	Retail Price
5G-200-__	100' Reel	14	\$1.713/ft
5A-200-__	500' Reel		\$1.663/ft



4-Hole Asepsis (Polyurethane)

Catalog No.	Package/Unit Size	Color	Retail Price
4G-415-__	100' Reel	14, 54, 72	\$1.584/ft
4A-415-__	500' Reel		\$1.534/ft



4-Hole (Polyurethane)

Catalog No.	Package/Unit Size	Color	Retail Price
4G-015-__	100' Reel	14, 54, 72	\$1.181/ft
4A-015-__	500' Reel		\$1.128/ft



4-Hole (PVC)

Catalog No.	Package/Unit Size	Color	Retail Price
4G-013-__	100' Reel	14, 54, 72	\$0.683/ft
4A-013-__	500' Reel		\$0.628/ft



3-Hole (Polyurethane)

Catalog No.	Package/Unit Size	Color	Retail Price
3G-017-__	100' Reel	14, 54, 72	\$0.972/ft
3A-017-__	500' Reel		\$0.920/ft



2-Hole Asepsis (Polyurethane)

Catalog No.	Package/Unit Size	Color	Retail Price
2G-413-__	100' Reel	14, 72	\$1.045/ft
2A-413-__	500' Reel		\$1.036/ft



2-Hole (Polyurethane)

Catalog No.	Package/Unit Size	Color	Retail Price
2G-013-__	100' Reel	14, 54, 72	\$0.831/ft
2A-013-__	500' Reel		\$0.780/ft



2-Hole Asepsis (PVC)

Catalog No.	Package/Unit Size	Color	Retail Price
2G-410-__	100' Reel	14, 54, 72	\$0.547/ft
2A-410-__	500' Reel		\$0.492/ft



2-Hole (PVC)

Catalog No.	Package/Unit Size	Color	Retail Price
2G-010-__	100' Reel	14, 54, 72	\$0.437/ft
2A-010-__	500' Reel		\$0.383/ft



Tubing

Vacuum & Saliva Ejector

Vacuum, Saliva Ejector (PVC)

Catalog No.	Package/Unit Size	I.D.	Color	Retail Price
1G-002-__	100' Reel	3/16"	01, 14, 54, 72, 92	\$0.416/ft
1A-002-__	500' Reel	3/16"		\$0.308/ft
1G-085-__	100' Reel	7/32"		\$0.355/ft
1A-085-__	500' Reel	7/32"		\$0.300/ft



Vacuum and Drain Tubing (PVC)

Catalog No.	Package/Unit Size	I.D.	Color	Retail Price
1A-004-__	200' Reel	7/16"	14, 72	\$0.938/ft
1A-164-__	100' Reel	5/8"	14, 54	\$1.776/ft



6-Sided Vacuum Drain (PVC)

Catalog No.	Package/Unit Size	I.D.	Color	Retail Price
1A-016-__	200' Reel	7/16"	01, 14, 54, 72, 92	\$1.311/ft
1A-008-__	100' Reel	5/8"	54 only	\$2.262/ft



Asepsis-Flex (PVC)

Catalog No.	Package/Unit Size	I.D.	Color	Retail Price
1E-501-__	25' Length	1/2"	14, 54, 72, 92	\$1.963/ft
1F-501-__	100' Length	1/2"		\$1.617/ft
1E-625-__	25' Length	5/8"		\$2.887/ft
1F-625-__	100' Length	5/8"		\$2.469/ft



Flexi-Corr® Hose & Cuffs

Catalog No.	Package/Unit Size	I.D.	Color	Cuff No.	Retail Price/ft
1E-047-__	25' Length	1/2"	04, 54, 72, 92	N/A	\$1.694
1F-047-__	250' Box (5 ea. 50' Length)				\$1.425
1E-053-__	25' Length	5/8"	"	OC-058	\$1.842
1F-053-__	150' Box (3 ea. 50' Length)				\$1.599
1E-054-__	25' Length	3/4"	"	OC-034	\$2.258
1F-054-__	50' Length				\$1.855
1E-055-__	25' Length	1"	"	OC-100	\$2.420
1F-055-__	50' Length				\$2.017
1E-056-__	25' Length	1 1/4"	"	OC-114	\$2.715
1F-056-__	50' Length				\$2.312
1E-057-__	25' Length	1 1/2"	"	OCH-112	\$3.307
1F-057-__	50' Length				\$2.903
1E-058-__	25' Length	1 3/4"	"	N/A	\$3.763
1F-058-__	50' Length				\$3.360
1E-093-__	25' Length	2"	"	OC-200	\$5.080
1F-093-__	50' Length				\$4.543



Tubing Cuffs

Catalog No.	I.D.	Color	Retail Price	Cuff No.	I.D.	Retail Price
OC-012	1/2"	cuff	\$3.76/ea	OC-114	1 1/4"	\$3.02/ea
OC-058	5/8"	4 (gray)	\$2.52/ea	OCH-112	1 1/2"	\$4.02/ea
OC-034	3/4"		\$3.95/ea	OC-200	2"	\$6.87/ea
OC-100	1"		\$2.85/ea	Custom designs and sizes are available.		



85A (Softer Polyurethane)

Catalog No.	Package/Unit Size	O.D.	Color	Retail Price
1L-013-__	250' Box	1/8"	01, 02, 04, 05,	\$0.158/ft
1G-013-__	250' Reel		06, 07, 08, 09,	\$0.156/ft
1C-013-__	500' Reel		10, 11, 12, 14,	\$0.121/ft
1A-013-__	2500' Reel		25, 26, 27, 28, 28, 29, 32	\$0.107/ft

85A (Softer Polyurethane)

Catalog No.	Package/Unit Size	O.D.	Color	Retail Price
1G-006-__	100' Reel	.240"	01, 02, 04, 06,	\$0.400/ft
1B-006-__	500' Reel		07, 08, 09, 10,	\$0.345/ft
1A-006-__	1000' Reel		11, 14, 25, 26, 27, 54	\$0.320/ft

90A (Harder Polyurethane)

Catalog No.	Package/Unit Size	O.D.	Color	Retail Price
1D-132-__	100' Box	1/4"	01, 04, 05, 06, 07,	\$0.411/ft
1B-132-__	500' Reel		08, 09, 10, 11, 14,	\$0.345/ft
1A-132-__	1000' Reel		25, 27, 28	\$0.320/ft

90A (Harder Polyurethane)

Catalog No.	Package/Unit Size	O.D.	Color	Retail Price
1G-005-__	100' Reel	3/8"	01, 04, 05, 06, 07,	\$0.808/ft
1B-005-__	250' Reel		08, 09, 10, 11, 14,	\$0.742/ft
1A-005-__	500' Reel		25, 26, 27, 28, 48, 49	\$0.703/ft

Braided Supply Tubing (PVC)

Catalog No.	Package/Unit Size	O.D.	Color	Retail Price
1G-024-__	100' Reel	1/4"	04	\$1.229/ft
1A-024-__	1000' Reel			\$1.092/ft

Water Tube

Catalog No.	Package/Unit Size	O.D.	Color	Retail Price
1G-009-__	100' Reel	.090"	07, 10	\$0.167/ft
1A-009-__	2500' Reel			\$0.114/ft

Umbilical & Ultra-Sonic Scaler Bundle Variations

We offer any combination of wire and tubing into a single, compact jacketed bundle. It's cost-effective and a safe solution for routing water, air and electricity together.

Contact us for a custom quote.



Tubing

Air/Water Supply

Prophy Tubing, 7ft length

Catalog No.	Description	Color	Retail Price
4F-048-__	For Dentsply/Cavitron Prophy Jet Model C-100 & C-200	01, 14	\$26.83/ea
5F-048-__	For Dentstply/Cavitron Prophy Jet Model C300 & Cavi-Jet Model G105		\$31.20/ea



Supply Tubing, 3/16" OD x 1/8" ID, 85A (PUR)

Catalog No.	Package/Unit Size	O.D.	Color	Retail Price
1J-194-__	100' Bag	3/16"	01, 07, 10	\$0.305
1C-194-__	500' Reel	3/16"		\$0.289
1A-194-__	2000' Reel	3/16"		\$0.256



Hose Cutter

Catalog No.	Description	Retail Price
TUB-057-05	Hose Cutter, red	\$31.36
TUB-057-B	Replacement blade	\$8.32

Heavy duty Hose Cutter will cut plastic tubing or hose up to 1" O.D. reinforced hose can be cut perfectly with ease.



Tube Cutter

Catalog No.	Description	Retail Price
TUB-058-__	Tubing Cutter	\$9.26
TUB-058-B	Replacement blade	\$3.43

The most efficient and accurate tubing cutter available. It's durable enough to cut hard plastic tubing yet lightweight and small enough for convenient storage. Suitable for 1/8" to 3/4" O.D. tubing. Available in black (O1) or blue (O7).



"H" Clips

Our "H" Clips are an easy and inexpensive way to route dual tubing or hose runs. Conveniently attached via the countersunk center screw hole (3/16"). Tubing "H" Clips neatly and safely secure tubing for a clean professional look. Manufactured from black Polyacetal plastic with a useful temperature range of -40° to + 212°F.

Catalog No.	Description	Retail Price
HTC-156-200	.156" - .200"	\$0.74
HTC-201-325	.201" - .325"	\$0.92
HTC-326-400	.326" - .400"	\$0.93
HTC-401-485	.401" - .485"	\$0.95
HTC-486-575	.486" - .575"	\$1.10
HTC-576-675	.576" - .675"	\$1.57
HTC-676-825	.676" - .825"	\$1.84
HTC-826-1000	.826" - 1.00"	\$2.53



Tube Channels

Catalog No.	Description	Retail Price
TC3-10	3/16"	\$2.22
TC4-10	1/4"	\$2.26
TC5-10	5/32"	\$2.31
PTC-6-7	3/8" OD - 7 channels	\$5.43
PTC-6-10	3/8" OD - 10 channels	\$7.28
PTC-6-6	1/2" OD - 6 channels	\$7.61

Tube Channel mounts designed to secure multiple tubing lines. The small 10-tube channel can be mounted with its pressure sensitive adhesive backing or with a #4 screw. Available in 3 sizes to fit 5/32", 3/16", and 1/4" tubing. The large style channel holds 3/8" tubing (7 or 10 channels) and 1/2" tubing (6 channels) and mounts with 2 screws.



Tube Stripper

Catalog No.	Description	Retail Price
TUB-059	Tubing Stripper	\$1.80



Umbilical Tubing Stringer

• Used to pull 1/8" and 1/4" tubing through existing umbilical

Catalog No.	Description	Retail Price
TUB-060	Umbilical Tubing Stringer, 12' long	\$34.00



Tube Gauge

• Tube Gauge - a handy tool that provides a quick way to measure tubing O.D. and I.D. It measures in standard and metric sizes.

Catalog No.	Description	Retail Price
TUB-065	Tube Gauge	\$2.25



Coiled Bib Holder

• Professional, functional bib holders are easy to use and comfortable to wear
• Made from high quality polyurethane

Catalog No.	Package/Unit Size	Color	Retail Price
OV-115-__	each	07, 10, 11, 25, 26 27, 32, 45, 46, 48, 49	\$3.83



Catalog No.	Page	Retail Price
1000	10	\$4,827
1007	12	\$66
1014	11	\$116
2000	10	\$5,240
3000	10	\$6,234
4064	13	\$154
4091	13	\$344
2000-060	11	\$3,378
2000-060E	11	\$3,986
2000-060S	11	\$3,986
2000-060SA	11	\$3,492
2000-060T	11	\$3,986
2000-226	11	\$286
2000-243	11	\$228
2000-244	11	\$404
4000-055	13	\$32
4003-123	13	\$10
4003-193-1	13	\$186
4003-240	13	\$8
4003-275	13	\$54
4003-276	13	\$36
4091-7	13	\$212
4095-4	13	\$370
5025-1	13	\$144
7000-021	11	\$116
7000-96-1.1	11	\$108
A-3100	21	\$2,820
A-3101	21	\$2,670
A-3150	21	\$3,400
A-3151	21	\$2,790
AT-96B	9	\$466
CLX-100	34	\$1,250
DR-96B	9	\$342
DR-96T	9	\$739
FPM103G	12	\$408
HG5-001	36	\$860
HG5-002C	36	\$197
HG5-002CR	36	\$338
HG5HV-001	36	\$2,860
HG5-Mini	36	\$860
INT-2001	10	\$7,520
LED-105-CL	34	\$2,790
LED-105-CM	34	\$3,230
LED-105-P	34	\$2,200
LED-105-W	34	\$3,230
N-96S	9	\$337
US-INS-01	35	\$96
US-INS-02	35	\$96
USM- 050A	28	\$340
USM- 100A	28	\$510
USM- 101A	28	\$320
USM- 200	16	\$735
USM- 200A	28	\$610
USM- 200P	25	\$750

Catalog No.	Page	Retail Price
USM- 202A	28	\$650
USM- 220A	28	\$690
USM- 220P	25	\$860
USM- 222A	28	\$740
USM- 250	16	\$815
USM- 300	24	\$1,465
USM- 3200	33	\$1,015
USM- 3250	24	\$1,530
USM- 330	24	\$1,530
USM- 3300	15	\$3,295
USM- 3305	17	\$2,090
USM- 3310	15	\$3,295
USM- 3350	17	\$2,250
USM- 3360	19	\$3,400
USM- 400A	28	\$1,020
USM- 400P	25	\$1,155
USM- 4100	22	\$2,365
USM- 4100A	23	\$2,730
USM- 4150	22	\$2,630
USM- 4150A	23	\$3,265
USM- 4250	22	\$4,000
USM- 4350	18	\$2,600
USM- 4500A	23	\$1,260
USM- 4500P	26	\$1,550
USM- 4501P	22	\$1,520
USM- 4550A	23	\$1,610
USM- 4550P	26	\$1,730
USM- 4551P	27	\$1,720
USM- 4700	20	\$1,435
USM- 4710	20	\$1,630
USM- 5000A	31	\$740
USM- 5035A	31	\$765
USM- 5100	18	\$790
USM- 5125A	31	\$720
USM- 5130A	31	\$750
USM- 5200A	30	\$950
USM- 5220A	30	\$800
USM- 5222A	30	\$870
USM- 5230A	32	\$1,015
USM- 5300	16	\$1,910
USM- 5525A	32	\$1,000
USM- 6000	14	\$3,800
USM- 6010	14	\$3,800
USM- 4062	13	\$972
USM- 4063	13	\$1,678
USM-CA-350	8	\$3,500
USM-CA-450	8	\$3,900
USM-CA-550	8	\$4,200
USM-OP1	10	\$460
USM-OP2	10	\$440
USM-OP3	10	\$425
USM-PC-00	10	\$116
USM-PC-01	10	\$116
USM-PC-02	10	\$602

Catalog No.	Page	Retail Price
USM-PC-03C	10	\$426
USM-PC-03D	10	\$626
USM-PC-04	10	\$284
USM-PC-06	10	\$108
USM-PC-07	10	\$460
USM-PC-08	10	\$246
USM-PC-09	10	\$456
USM-PC-10	10	\$278
USM-PC-11	10	\$116
USM-SP200	12	\$3,120
USM-SP400	12	\$3,660
US-SEI-01	35	\$665
US-SEI-02	35	\$665

DENTAL COMPONENTS & TUBING

Index

Catalog No.	Page	Retail Price	Catalog No.	Page	Retail Price	Catalog No.	Page	Retail Price
1600008	120	\$448.00	1E-055-xx	136	\$2.420	2V-005-xx	132	\$13.74
1600009	120	\$503.00	1E-056-xx	136	\$2.715	2V-006-xx	134	\$9.08
1600010	120	\$1,090.00	1E-057-xx	136	\$3.307	2V-041-xx	132	\$8.86
1600011	120	\$358.00	1E-058-xx	136	\$3.763	2V-046-xx	134	\$7.30
1600012	120	\$358.00	1E-093-xx	136	\$5.080	3A-001-xx	132	\$0.463
1600013	120	\$358.00	1E-501-xx	136	\$1.963	3A-013-xx	134	\$0.498
1600014	120	\$420.00	1E-625-xx	136	\$2.887	3A-017-xx	135	\$0.920
1600015	120	\$595.00	1F-047-xx	136	\$1.425	3A-042-xx	132	\$0.665
1600110	120	\$409.00	1F-053-xx	136	\$1.599	3A-047-xx	132	\$0.464
1600210	120	\$200.00	1F-054-xx	136	\$1.855	3A-413-xx	134	\$0.589
1600310	120	\$298.00	1F-055-xx	136	\$2.017	3B-013-xx	134	\$0.518
1610810	120	\$184.00	1F-056-xx	136	\$2.312	3B-413-xx	134	\$0.617
1611201	120	\$14.00	1F-057-xx	136	\$2.903	3G-001-xx	132	\$0.519
2000003	119	\$372.00	1F-058-xx	136	\$3.360	3G-013-xx	134	\$0.550
2000004	119	\$372.00	1F-093-xx	136	\$4.543	3G-017-xx	135	\$0.972
2000005	119	\$423.00	1F-501-xx	136	\$1.617	3G-042-xx	132	\$0.716
2000006	119	\$423.00	1F-625-xx	136	\$2.469	3G-047-xx	132	\$0.519
2000015	119	\$670.00	1G-002-xx	136	\$0.416	3G-413-xx	134	\$0.639
2000016	119	\$670.00	1G-005-xx	137	\$0.808	3S-001-xx	132	\$21.84
2000101	119	\$51.00	1G-006-xx	137	\$0.400	3S-005-xx	54, 134	\$8.74
2000102	119	\$51.00	1G-009-xx	137	\$0.167	3S-042-xx	132	\$27.04
2000105	119	\$237.00	1G-013-xx	137	\$0.156	3S-047-xx	132	\$22.88
2000107	119	\$167.00	1G-024-04	137	\$1.229	3V-013-xx	134	\$9.22
2000107-4	119	\$151.00	1G-085-xx	136	\$0.355	3V-042-xx	132	\$10.15
2000107-6	119	\$167.00	1G-132-xx	137	\$0.411	3V-413-xx	134	\$8.80
2000107-6L	119	\$217.00	1G-502-xx	131	\$0.468	4A-001-xx	133	\$0.523
2000114	119	\$131.00	1G-585-xx	131	\$0.468	4A-013-xx	135	\$0.628
2000116	119	\$165.00	1J-194-xx	138	\$0.305	4A-015-xx	135	\$1.128
2000119	119	\$190.00	1L-013-xx	137	\$0.158	4A-042-xx	132	\$0.817
1A-002-xx	136	\$0.308	2A-001-xx	132	\$0.445	4A-047-xx	133	\$0.519
1A-004-xx	136	\$0.938	2A-005-xx	132	\$0.832	4A-415-xx	135	\$1.534
1A-005-xx	137	\$0.703	2A-006-xx	134	\$0.269	4A-549-xx	131	\$0.846
1A-006-xx	137	\$0.320	2A-010-xx	135	\$0.383	4E-533-xx	131	\$20.06
1A-008-xx	136	\$2.262	2A-013-xx	135	\$0.780	4E-82347-xx-10	121	\$90.00
1A-009-xx	137	\$0.114	2A-041-xx	132	\$0.613	4E-82347-xx-10S	121	\$130.00
1A-013-xx	137	\$0.107	2A-046-xx	134	\$0.295	4E-82347-xx-7	121	\$80.00
1A-016-xx	136	\$1.311	2A-047-xx	132	\$0.492	4E-82347-xx-7S	121	\$120.00
1A-024-04	137	\$1.092	2A-410-xx	135	\$0.492	4F-048-xx	138	\$26.83
1A-085-xx	136	\$0.300	2A-413-xx	135	\$1.036	4G-001-xx	133	\$0.582
1A-132-xx	137	\$0.320	2A-549-xx	131	\$0.621	4G-013-xx	135	\$0.683
1A-164-xx	136	\$1.776	2B-006-xx	134	\$0.295	4G-015-xx	135	\$1.181
1A-194-xx	138	\$0.256	2B-046-xx	134	\$0.302	4G-042-xx	132	\$0.870
1A-502-xx	131	\$0.462	2B-549-xx	131	\$0.676	4G-047-xx	133	\$0.573
1A-516-xx	131	\$1.961	2G-001-xx	132	\$0.469	4G-415-xx	135	\$1.584
1A-534-xx	131	\$1.142	2G-006-xx	134	\$0.305	4G-549-xx	131	\$0.899
1A-585-xx	131	\$0.462	2G-010-xx	135	\$0.437	4S-001-xx	133	\$19.34
1B-005-xx	137	\$0.742	2G-013-xx	135	\$0.831	4S-042-xx	133	\$24.44
1B-006-xx	137	\$0.345	2G-041-xx	132	\$0.664	4S-549-xx	131	\$22.88
1B-132-xx	137	\$0.345	2G-046-xx	134	\$0.320	4V-042-xx	132	\$11.05
1C-013-xx	137	\$0.121	2G-047-xx	132	\$0.559	4V-050-14	121	\$90.00
1C-194-xx	138	\$0.289	2G-410-xx	135	\$0.547	4V-060-xx	121	\$90.00
1E-047-xx	136	\$1.694	2G-413-xx	135	\$1.045	4V-061-xx	121	\$90.00
1E-053-xx	136	\$1.842	2S-005-xx	54, 134	\$13.42	4V-062-xx	121	\$125.00
1E-054-xx	136	\$2.258	2V-001-xx	132	\$8.29	4V-063-xx	121	\$225.00

DENTAL COMPONENTS & TUBING

Index

Catalog No.	Page	Retail Price	Catalog No.	Page	Retail Price	Catalog No.	Page	Retail Price
4V-070-xx-07	122	\$164.00	ACF-574-01	67, 92	\$96.50	ASF-103-04	43	\$185.00
4V-070-xx-07B	122	\$211.00	ACF-574-02	67, 92	\$91.50	ASF-103-11	42	\$58.80
4V-070-xx-12	122	\$190.00	ACF-574-03	67, 92	\$5.00	ASF-103-12	42	\$58.80
4V-070-xx-12B	122	\$237.00	ACF-574-05	67, 92	\$24.00	ASF-104-01	43	\$40.00
4V-080-xx	122	\$110.00	ACF-574-06	67, 92	\$28.00	ASF-104-02	43	\$35.00
4V-090-xx	122	\$136.00	ACF-575-01	67, 92	\$83.00	ASF-105-01	43	\$150.00
4V-100-xx	122	\$100.00	ACF-581-01	48	\$370.00	ASF-105-02	43	\$190.00
4V-110-xx	122	\$125.00	ACF-582-01	49	\$370.00	ASF-106-01	44	\$20.00
5A-200-xx	135	\$1.663	ACF-582-03	49	\$8.00	ASF-106-350	41	\$160.00
5G-200-xx	135	\$1.713	ACF-582-04	49	\$40.00	ASF-107-01	44	\$60.00
5F-048-xx	138	\$31.20	ACF-582-05	49	\$12.00	ASF-107-02	44	\$80.00
ACF-375-01	92	\$140.00	ACF-582-06	49	\$8.00	ASF-107-03	44	\$64.00
ACF-375-02	92	\$170.00	ACF-584-01	49	\$3.70	ASF-107-04	44	\$70.00
ACF-541-06	91	\$28.50	ACF-584-03	49	\$19.50	ASF-107-11	44	\$70.00
ACF-541-07	91	\$0.60	ACF-584-04	49	\$19.50	ASF-107-12	44	\$80.00
ACF-542-01	91	\$155.00	ACF-584-05	49	\$19.50	ASF-108-01	44	\$90.00
ACF-542-01D	91	\$185.00	ACF-592-01	49	\$94.00	ASF-108-02	44	\$136.50
ACF-542-02	91	\$160.00	ACF-592-02	49	\$180.00	ASF-109-01	44	\$360.00
ACF-542-02D	91	\$190.00	ACF-592-03	49	\$185.00	ASF-109-02	45	\$320.00
ACF-542-03	91	\$160.00	ACF-593-01	50	\$140.00	ASF-109-02K	45	\$595.00
ACF-542-03D	91	\$190.00	ACF-593-01M	50	\$136.00	ASF-109-04	45	\$340.00
ACF-542-AA	91	\$64.00	ACF-593-02	50	\$170.00	ASF-109-04K	45	\$880.00
ACF-542-AW	91	\$64.00	ACF-593-03M	50	\$145.00	ASF-110-01	45	\$154.00
ACF-543-01	91	\$7.00	ACF-593-01MC	50	\$28.00	ASF-111-01	45	\$800.00
ACF-543-02	91	\$12.00	ACF-593-04M	50	\$190.00	ASF-111-02	45	\$800.00
ACF-543-03	91	\$30.00	ACF-593-06	50	\$150.00	ASF-111-03	45	\$800.00
ACF-543-03L	91	\$25.50	ACF-593-06M	50	\$210.00	ASF-111-04	45	\$800.00
ACF-543-04	92	\$70.00	ACF-593-07	50	\$171.00	ASF-112-01	45	\$800.00
ACF-543-05	91	\$7.40	ACF-593-07M	50	\$219.00	ASF-112-02	45	\$800.00
ACF-543-06	91	\$3.00	ACF-594-01	50	\$25.00	ASF-112-03	45	\$800.00
ACF-543-07	91	\$46.00	ACF-594-02	50	\$27.50	ASF-112-04	45	\$800.00
ACF-543-08	92	\$2.90	ACF-594-03	50	\$30.00	ASF-112-05	45	\$800.00
ACF-543-10	92	\$21.70	ACF-600-01	50	\$300.00	ASF-112-06	45	\$800.00
ACF-543-11	92	\$24.20	ACF-601-01	50	\$230.00	ASF-121-02	46	\$520.00
ACF-544-01	91	\$225.00	ACF-601-02	50	\$240.00	ASF-121-03	46	\$520.00
ACF-544-02	91	\$215.00	ASF-100-01	41	\$940.00	ASF-122-01	46	\$280.00
ACF-545-01	91	\$320.00	ASF-100-02	41	\$940.00	ASF-122-02	46	\$280.00
ACF-545-02	91	\$280.00	ASF-100-03	41	\$940.00	ASF-123-01	46	\$280.00
ACF-545-03	91	\$325.00	ASF-100-04	41	\$950.00	ASF-123-02	46	\$280.00
ACF-564-01	125	\$20.00	ASF-100-04TP	41	\$950.00	ASF-123-03	46	\$280.00
ACF-564-02	125	\$20.00	ASF-100-05	41	\$240.00	ASF-133-01	46	\$175.00
ACF-566-01	125	\$204.00	ASF-100-06	41	\$240.00	ASF-134-01	46	\$185.00
ACF-566-02	125	\$8.30	ASF-100-07	41	\$250.00	ASF-155-01	46	\$170.00
ACF-566-03	125	\$12.40	ASF-100-08	41	\$250.00	ASF-155-03	46	\$170.00
ACF-566-10	125	\$39.00	ASF-101-01	42	\$310.00	ASF-156-01	47	\$330.00
ACF-566-11	125	\$41.00	ASF-101-02	42	\$230.00	ASF-160-01	47	\$102.00
ACF-566-12	125	\$52.00	ASF-101-03	42	\$459.00	ASF-161-01	47	\$150.00
ACF-566-13	125	\$58.50	ASF-102-01	42	\$62.00	ASF-161-02	47	\$165.00
ACF-566-14	125	\$58.50	ASF-102-02	42	\$64.00	ASF-162-01	47	\$198.00
ACF-566-15	125	\$52.00	ASF-102-03	42	\$70.00	ASF-163-01	47	\$164.00
ACF-566-16	125	\$52.00	ASF-102-04	42	\$76.00	ASF-164-00	48	\$63.30
ACF-566-17	125	\$52.00	ASF-103-01	42	\$58.80	ASF-164-01	48	\$39.50
ACF-567-01	125	\$28.00	ASF-103-02	41, 43	\$135.00	ASF-164-02	48	\$5.00
ACF-567-02	125	\$22.00	ASF-103-03	43	\$50.00	ASF-165-01	48	\$18.00

DENTAL COMPONENTS & TUBING

Index

Catalog No.	Page	Retail Price	Catalog No.	Page	Retail Price	Catalog No.	Page	Retail Price
ASF-165-02	48	\$9.00	FIF-603-08	115	\$1.00	FIF-627-02	103	\$9.30
ASF-171-01	48	\$120.00	FIF-603-09	115	\$1.50	FIF-627-03	103	\$9.30
ASF-171-02	48	\$4.00	FIF-603-10	115	\$3.50	FIF-627-04	103	\$9.30
ASF-171-03	48	\$33.00	FIF-604-01	115	\$2.20	FIF-628-01	103	\$4.50
ASF-171-04	48	\$33.00	FIF-604-02	115	\$3.30	FIF-628-02	103	\$4.50
ASF-172-01	48	\$102.00	FIF-604-03	115	\$4.40	FIF-628-03	103	\$7.00
ASF-172-02	48	\$127.00	FIF-604-04	115	\$8.60	FIF-629-01	80 , 103	\$5.50
FCF-251-01E-x-SII	55	\$68.00	FIF-605-01	99	\$1.50	FIF-630-01	104	\$14.00
FCF-251-01S-x-SII	55	\$77.00	FIF-605-02	99	\$1.95	FIF-631-01	104	\$8.25
FCF-251-01-x-SI	55	\$75.00	FIF-606-01	99	\$3.20	FIF-631-02	104	\$8.25
FCF-251-02S-x-SII	55	\$107.00	FIF-606-02	99	\$24.00	FIF-631-03	104	\$8.25
FCF-251-02-x-SI	55	\$110.00	FIF-607-01	100	\$4.00	FIF-632-01	104	\$10.30
FCF-252-01-x-SI	55	\$170.00	FIF-611-01	100	\$4.50	FIF-632-02	104	\$10.30
FCF-252-01-x-SII	55	\$130.00	FIF-611-02	100	\$38.00	FIF-632-03	104	\$12.40
FCF-252-02-x-SI	55	\$187.00	FIF-612-01	100	\$1.10	FIF-632-04	104	\$12.40
FCF-252-02-x-SII	55	\$160.00	FIF-612-02	100	\$1.50	FIF-633-01	104	\$1.50
FCF-252-03-x-SII	55	\$130.00	FIF-612-03	100	\$2.00	FIF-633-02	104	\$1.00
FCF-252-04-x-SII	55	\$160.00	FIF-612-04	100	\$2.00	FIF-634-01	104	\$5.20
FCF-253-01-x-SI	55	\$170.00	FIF-614-01	100	\$5.70	FIF-634-02	104	\$5.20
FCF-253-01-x-SII	55	\$130.00	FIF-614-02	100	\$6.00	FIF-634-03	104	\$5.20
FCF-253-02-x-SI	55	\$187.00	FIF-615-01	100	\$4.90	FIF-634-04	104	\$5.20
FCF-253-02-x-SII	55	\$160.00	FIF-616-01	100	\$7.50	FIF-635-01	104	\$6.00
FCF-261-01-x-SI	55	\$180.00	FIF-616-02	100	\$6.50	FIF-636-01	105	\$1.60
FCF-261-01-x-SII	55	\$165.00	FIF-617-01	100	\$3.60	FIF-636-02	105	\$1.60
FCF-261-02-x-SI	55	\$197.00	FIF-618-01	100	\$3.80	FIF-636-05	105	\$3.20
FCF-261-02-x-SII	55	\$195.00	FIF-621-01	100	\$2.30	FIF-637-01	105	\$8.00
FCF-263-01-SI	56	\$70.00	FIF-621-02	100	\$5.00	FIF-638-01	105	\$1.60
FCF-263-01-SII	56	\$24.00	FIF-622-01	100	\$2.30	FIF-638-02	105	\$1.60
FCF-263-02-SI	56	\$50.00	FIF-623-01	101, 116	\$2.40	FIF-639-01	105	\$4.70
FCF-263-02-SII	56	\$24.00	FIF-623-02	101	\$5.00	FIF-639-02	105	\$4.80
FCF-264-01-SI	56	\$35.00	FIF-624-02	101	\$6.00	FIF-639-03	105	\$8.00
FCF-264-01-SII	56	\$24.00	FIF-624-03	101	\$7.00	FIF-639-04	105	\$11.00
FCF-265-01-SI	56	\$50.00	FIF-624-04	101	\$7.00	FIF-641-01	105	\$5.00
FCF-265-02-SI	56	\$58.00	FIF-624-05	101	\$7.50	FIF-641-02	105	\$7.50
FCF-265-03K-SII	56	\$4.50	FIF-624-06	101	\$5.40	FIF-641-03	105	\$6.50
FCF-265-03-SII	56	\$3.00	FIF-624-07	101	\$7.00	FIF-641-04	105	\$8.00
FCF-270-01	56	\$195.00	FIF-624-10	102	\$21.00	FIF-642-01	105	\$3.20
FCF-852-01-SI	56	\$15.00	FIF-624-11	102	\$21.00	FIF-642-02	105	\$3.90
FCF-852-02	56	\$15.00	FIF-624-12	102	\$21.00	FIF-642-03	105	\$7.00
FCF-854-01-SI	56	\$21.00	FIF-624-20	102	\$26.00	FIF-642-04	105	\$8.00
FCF-855-01-SII	57	\$7.80	FIF-624-21	102	\$26.00	FIF-643-01	105	\$3.30
FCF-856-01-SII	57	\$8.90	FIF-624-22	102	\$26.00	FIF-644-01	105	\$6.00
FCF-857-01-SII	57	\$11.50	FIF-624-30	102	\$14.00	FIF-644-02	105	\$6.00
FCF-858-01-SII	57	\$9.00	FIF-624-40	102	\$19.50	FIF-645-01	106	\$6.80
FCF-859-01-SII	57	\$6.50	FIF-624-41	102	\$17.00	FIF-645-02	106	\$7.00
FCF-859-02	57	\$44.00	FIF-625-01	103	\$5.30	FIF-646-01	106	\$8.30
FCF-860-01	57	\$24.00	FIF-625-02	103	\$8.20	FIF-646-10	106	\$4.50
FCF-860-02	57	\$16.50	FIF-625-03	103	\$10.00	FIF-646-15	106	\$3.60
FCF-860-03	57	\$17.00	FIF-625-04	103	\$30.00	FIF-647-01	106	\$3.00
FCF-860-04	57	\$22.00	FIF-626-01	103	\$9.20	FIF-647-02	106	\$4.00
FIF-601-01	115	\$1.00	FIF-626-02	103	\$10.30	FIF-647-03	106	\$7.60
FIF-601-02	115	\$1.10	FIF-626-03	103	\$12.40	FIF-647-04	106	\$7.60
FIF-601-03	115	\$3.20	FIF-626-04	103	\$13.40	FIF-648-01	106	\$2.60
FIF-602-01	115	\$4.60	FIF-627-01	103	\$8.25	FIF-649-01	106	\$1.35

DENTAL COMPONENTS & TUBING

Index

Catalog No.	Page	Retail Price	Catalog No.	Page	Retail Price	Catalog No.	Page	Retail Price
FIF-649-02	106	\$2.00	FIF-666-04	108	\$9.00	FIF-674-06	110	\$0.60
FIF-649-03	106	\$2.50	FIF-667-01	108	\$2.70	FIF-674-07	110	\$0.60
FIF-649-04	106	\$4.10	FIF-667-02	108	\$3.00	FIF-675-01	110	\$3.00
FIF-650-01	106	\$1.00	FIF-667-03	108	\$3.10	FIF-675-02	110	\$3.00
FIF-650-02	106	\$1.30	FIF-667-04	108	\$7.50	FIF-675-03	110	\$5.50
FIF-650-03	106	\$1.90	FIF-667-05	108	\$4.90	FIF-676-01	111	\$2.30
FIF-651-01	106	\$1.50	FIF-667-06	108	\$5.00	FIF-676-02	111	\$2.10
FIF-651-02	106	\$3.20	FIF-667-07	108	\$6.00	FIF-676-03	111	\$2.00
FIF-651-03	106	\$2.00	FIF-667-10	108	\$20.00	FIF-676-04	111	\$4.10
FIF-651-04	106	\$2.00	FIF-668-01	109	\$6.20	FIF-676-05	111	\$4.20
FIF-651-05	106	\$3.80	FIF-668-02	109	\$6.20	FIF-677-01	111	\$1.50
FIF-651-06	106	\$2.70	FIF-668-03	109	\$6.75	FIF-677-02	111	\$3.50
FIF-651-07	106	\$8.00	FIF-668-04	109	\$15.50	FIF-677-03	111	\$2.65
FIF-651-08	106	\$12.00	FIF-669-01	109	\$5.80	FIF-677-04	111	\$4.35
FIF-651-09	106	\$4.80	FIF-669-02	109	\$5.80	FIF-678-01	111	\$4.20
FIF-652-01	107	\$5.00	FIF-669-03	109	\$6.00	FIF-679-01	111	\$2.30
FIF-652-02	107	\$4.50	FIF-669-04	109	\$18.00	FIF-679-02	111	\$3.10
FIF-652-03	107	\$3.00	FIF-669-05	109	\$17.00	FIF-679-03	111	\$2.30
FIF-652-04	107	\$3.40	FIF-670-01	109	\$2.50	FIF-679-05	111	\$0.80
FIF-653-01	107	\$8.20	FIF-670-02	109	\$3.50	FIF-679-06	111	\$1.00
FIF-653-02	107	\$3.80	FIF-670-03	109	\$3.00	FIF-680-01	111	\$6.20
FIF-653-03	107	\$8.20	FIF-670-04	109	\$3.80	FIF-680-02	111	\$6.20
FIF-654-01	107	\$0.80	FIF-670-10	109	\$2.20	FIF-681-01	111	\$8.00
FIF-654-02	107	\$1.20	FIF-670-11	109	\$2.90	FIF-681-02	111	\$10.50
FIF-654-03	107	\$1.70	FIF-670-12	109	\$2.70	FIF-682-01	111	\$11.00
FIF-654-04	107	\$2.00	FIF-670-13	109	\$2.80	FIF-682-02	111	\$21.00
FIF-655-01	107	\$2.00	FIF-670-14	109	\$3.20	FIF-682-03	112	\$14.40
FIF-655-02	107	\$1.70	FIF-670-15	109	\$4.80	FIF-682-04	112	\$15.50
FIF-655-03	107	\$3.20	FIF-670-20	109	\$3.50	FIF-682-05	112	\$15.50
FIF-655-04	107	\$6.88	FIF-670-21	109	\$4.00	FIF-682-06	112	\$15.50
FIF-657-01	107	\$1.35	FIF-670-22	109	\$5.50	FIF-682-07	112	\$15.50
FIF-657-02	107	\$2.05	FIF-670-23	109	\$5.50	FIF-682-08	112	\$15.50
FIF-657-03	107	\$2.90	FIF-670-30	109	\$3.60	FIF-682-09	112	\$15.50
FIF-658-01	107	\$2.10	FIF-670-31	109	\$4.00	FIF-683-01	112	\$13.20
FIF-658-02	107	\$3.50	FIF-670-32	109	\$9.00	FIF-683-02	112	\$24.00
FIF-662-01	107	\$2.50	FIF-670-33	109	\$9.00	FIF-684-01	112	\$7.00
FIF-662-02	107	\$3.00	FIF-671-01	109	\$1.80	FIF-684-02	112	\$8.50
FIF-662-04	107	\$4.10	FIF-671-02	109	\$3.30	FIF-684-03	112	\$9.70
FIF-663-01	108	\$2.00	FIF-671-03	109	\$3.30	FIF-685-01	112	\$16.00
FIF-663-02	108	\$2.70	FIF-671-04	109	\$7.50	FIF-686-01	113	\$43.00
FIF-663-03	108	\$2.80	FIF-672-01	110	\$5.40	FIF-687-01	113	\$17.00
FIF-663-04	108	\$3.00	FIF-672-02	110	\$5.50	FIF-687-02	78, 113	\$3.10
FIF-663-05	108	\$2.55	FIF-672-10	110	\$8.50	FIF-687-03	78, 113	\$3.10
FIF-663-06	108	\$3.80	FIF-672-11	110	\$6.00	FIF-687-04	78, 113	\$3.10
FIF-663-07	108	\$6.70	FIF-672-12	110	\$19.50	FIF-687-05	78, 113	\$3.10
FIF-663-08	108	\$7.40	FIF-673-01	110	\$4.50	FIF-687-10	113	\$10.00
FIF-664-01	108	\$5.10	FIF-673-02	110	\$5.00	FIF-687-11	113	\$16.00
FIF-664-02	108	\$7.90	FIF-673-03	110	\$4.50	FIF-687-20	113	\$5.00
FIF-664-03	108	\$8.00	FIF-673-04	110	\$1.30	FIF-687-21	113	\$5.00
FIF-665-01	88, 108	\$3.20	FIF-674-01	110	\$4.55	FIF-687-22	113	\$5.00
FIF-665-01	108	\$3.20	FIF-674-02	110	\$4.55	FIF-687-23	113	\$4.20
FIF-666-01	108	\$6.00	FIF-674-03	110	\$4.55	FIF-688-01	113	\$40.00
FIF-666-02	108	\$6.80	FIF-674-04	110	\$0.50	FIF-688-02	113	\$50.00
FIF-666-03	108	\$10.00	FIF-674-05	110	\$0.60	FIF-688-03	113	\$2.50

DENTAL COMPONENTS & TUBING

Index

Catalog No.	Page	Retail Price	Catalog No.	Page	Retail Price	Catalog No.	Page	Retail Price
FIF-688-04	113	\$3.50	HOF-401-04	71	\$17.00	HOF-455-02	75	\$120.00
FIF-689-01	114	\$2.00	HOF-401-05	71	\$30.00	HOF-455-03	75	\$135.00
FIF-689-02	114	\$2.00	HOF-402-01	71	\$6.00	HOF-456-01	75	\$76.00
FIF-689-03	114	\$6.30	HOF-403-01	71	\$7.00	HOF-456-02	75	\$105.00
FIF-689-04	114	\$2.00	HOF-404-01	71	\$6.00	HOF-456-03	75	\$95.00
FIF-689-05	114	\$8.00	HOF-404-06	71	\$7.00	HOF-457-01	75	\$265.00
FIF-689-06	114	\$2.50	HOF-404-07	71	\$7.00	HOF-457-02	75	\$215.00
FIF-691-01	114	\$1.40	HOF-405-01	72	\$5.00	HOF-457-03	75	\$230.00
FIF-691-02	114	\$2.00	HOF-411-01	72	\$10.90	HOF-457-04	75	\$175.00
FIF-691-03	114	\$2.30	HOF-411-02	72	\$10.90	HOF-458-01	75	\$22.20
FIF-691-04	114	\$2.10	HOF-411-03	72	\$11.50	HPA-100-01	122	\$210.00
FIF-691-05	114	\$2.10	HOF-411-04	72	\$11.50	HPA-101-01	122	\$265.00
FIF-691-06	114	\$2.10	HOF-415-01	72	\$73.00	HPA-105-02	123	\$68.00
FIF-691-07	114	\$2.70	HOF-415-02	72	\$93.00	HPA-105-05	123	\$310.00
FIF-691-08	114	\$2.50	HOF-416-01	72	\$68.00	HPA-106-01	123	\$46.40
FIF-692-01	114	\$2.10	HOF-416-02	72	\$73.00	HPA-106-02	123	\$46.40
FIF-692-02	114	\$2.60	HOF-417-01	72	\$33.00	HPA-107-01	123	\$340.00
FIF-692-03	114	\$1.80	HOF-417-02	72	\$35.00	HPA-107-02	123	\$340.00
FIF-692-04	114	\$2.60	HOF-433-03	73	\$140.00	HRB-110-01	96	\$68.20
FIF-692-05	114	\$2.60	HOF-433-04	73	\$210.00	HRB-110-02	96	\$68.20
FIF-692-06	114	\$3.10	HOF-433-05	73	\$240.00	HRB-110-03	96	\$68.20
FIF-692-07	114	\$3.10	HOF-434-01	73	\$14.50	HRB-220-01	96	\$36.00
FIF-692-08	114	\$3.10	HOF-434-02	73	\$17.00	HRB-220-02	96	\$36.00
FIF-693-01	114	\$2.30	HOF-434-03	73	\$19.00	HRB-220-03	96	\$70.00
FIF-693-02	114	\$3.50	HOF-434-04	73	\$21.00	HRB-220-04	96	\$28.00
FIF-694-01	114	\$1.60	HOF-434-05	73	\$25.00	HTC-156-200	139	\$0.74
FIF-694-02	114	\$1.80	HOF-434-06	73	\$28.00	HTC-201-325	139	\$0.92
FIF-694-03	114	\$1.80	HOF-435-01	73	\$25.00	HTC-326-400	139	\$0.93
FIF-695-01	114	\$2.30	HOF-435-02	73	\$30.00	HTC-401-485	139	\$0.95
FIF-695-02	114	\$1.80	HOF-435-03	73	\$64.00	HTC-486-575	139	\$1.10
FIF-696-01	114	\$4.10	HOF-436-01-G	73	\$3.00	HTC-576-675	139	\$1.57
FVD-101-01	83	\$13.70	HOF-436-01-x	73	\$6.00	HTC-676-825	139	\$1.84
FVD-101-02	83	\$13.70	HOF-441-01	73	\$40.00	HTC-826-1000	139	\$2.53
FVD-101-03	83	\$16.90	HOF-441-02	73	\$40.00	LIT-100-01	93	\$401.00
HOF-381-00	69	\$37.60	HOF-442-01	74	\$22.00	LIT-100-02	93	\$550.00
HOF-381-01	69	\$40.00	HOF-442-02	74	\$22.00	LIT-100-03	93	\$370.00
HOF-381-02	69	\$40.00	HOF-443-01	74	\$36.00	LIT-100-04	93	\$370.00
HOF-381-03	69	\$40.00	HOF-444-01	74	\$9.00	LIT-100-05	93	\$550.00
HOF-381-04	69	\$33.00	HOF-450-01	74	\$12.40	LIT-101-00	94	\$13.90
HOF-381-05	69	\$33.00	HOF-450-02	74	\$12.40	LIT-101-01	94	\$21.00
HOF-381-10	69	\$37.60	HOF-451-01	74	\$57.50	LIT-101-02	94	\$8.00
HOF-381-11	69	\$42.00	HOF-452-01	74	\$136.00	LIT-101-03	94	\$7.90
HOF-391-01	70	\$50.00	HOF-452-02	74	\$136.00	LIT-101-04	94	\$58.00
HOF-391-03	70	\$40.20	HOF-453-01	74	\$136.00	LIT-101-05	94	\$14.00
HOF-391-04	70	\$34.00	HOF-453-02	74	\$136.00	LIT-101-06	94	\$11.40
HOF-391-04-SII	70	\$50.00	HOF-454-01	74	\$136.00	LIT-101-07	94	\$7.50
HOF-391-05	70	\$50.00	HOF-454-02	74	\$136.00	LIT-102-00	95	\$24.00
HOF-391-07	70	\$40.20	HOF-454-10	75	\$150.00	LIT-102-01	95	\$12.50
HOF-391-08	70	\$34.00	HOF-454-11	75	\$150.00	LIT-102-02	95	\$28.00
HOF-392-05	70	\$55.00	HOF-454-12	75	\$150.00	LIT-102-03	95	\$48.00
HOF-393-01	70	\$18.00	HOF-454-13	75	\$150.00	LIT-102-04	95	\$18.50
HOF-394-01	70	\$18.00	HOF-454-14	75	\$150.00	LIT-102-05	95	\$11.50
HOF-401-01	70	\$14.50	HOF-454-20	75	\$72.00	LIT-102-06	95	\$23.00
HOF-401-03	71	\$20.00	HOF-455-01	75	\$160.00	LIT-102-07	96	\$32.00

DENTAL COMPONENTS & TUBING

Index

Catalog No.	Page	Retail Price	Catalog No.	Page	Retail Price	Catalog No.	Page	Retail Price
LIT-102-08	96	\$51.00	SIF-814-02	116	\$7.00	SYR-031-04	52	\$8.00
LIT-102-09	96	\$23.00	SIF-815-02	116	\$33.00	SYR-233-01	51	\$4.00
LIT-103-00	97	\$30.50	SIF-815-03	116	\$11.00	SYR-234-01	51	\$23.50
LIT-103-01	97	\$43.00	SIF-816-01	116	\$16.50	SYR-234-02	51	\$200.00
LIT-103-02	97	\$43.00	SIF-817-01	116	\$24.00	SYR-234-03	51	\$27.50
LIT-103-03	97	\$36.00	SIF-821-01	117	\$27.00	SYR-234-05	51	\$255.00
LIT-103-04	97	\$22.00	SIF-821-02	117	\$34.00	SYR-235-01	52	\$4.90
LIT-103-05	97	\$41.00	SIF-823-01	117	\$45.00	SYR-235-02	52	\$24.00
LIT-103-06	97	\$19.00	SIF-830-01	117	\$65.00	SYR-235-03	52	\$205.00
OC-012	136	\$3.76	SIF-831-01	118	\$2.90	SYR-236-00	52	\$48.00
OC-034	136	\$3.95	SIF-831-02	118	\$2.90	SYR-236-0A	52	\$14.00
OC-058	136	\$2.52	SIF-831-03	118	\$2.90	SYR-236-0D	52	\$14.00
OC-100	136	\$2.85	SIF-831-04	118	\$2.90	SYR-236-0M	52	\$14.00
OC-114	136	\$3.02	SIF-831-05	118	\$2.90	SYR-236-0P	52	\$14.00
OC-200	136	\$6.87	SIF-831-06	118	\$2.90	SYR-241-01	52	\$9.00
OCH-112	136	\$4.02	SIF-831-07	118	\$2.90	SYR-241-03	52	\$9.20
OP-550-00	124	\$9.00	SIF-831-08	118	\$2.90	SYR-241-04	53	\$9.20
OP-550-00N	124	\$4.60	SIF-831-09	118	\$2.90	SYR-241-05	53	\$9.20
OP-560-00	124	\$10.50	SIF-831-10	118	\$2.90	SYR-241-06	53	\$11.00
OP-567-00N	124	\$4.60	SIF-831-11	118	\$2.90	SYR-242-01	53	\$2.50
OP-570-00	124	\$11.00	SIF-831-12	118	\$2.90	SYR-242-02	53	\$8.00
OP-823-00	121	\$36.80	SIF-831-13	118	\$2.90	SYR-242-03	53	\$12.00
OV-115-xx	139	\$3.83	SIF-831-14	118	\$2.90	SYR-243-01	53	\$18.90
PTC6-10	139	\$7.28	SIF-831-15	118	\$2.90	SYR-245-02	53	\$15.50
PTC6-6	139	\$7.61	SIF-831-16	118	\$2.90	SYR-246-01	54	\$16.00
PTC6-7	139	\$5.43	SIF-831-17	118	\$2.90	SYR-246-02	54	\$26.80
SIF-608-01	115	\$2.15	SIF-831-18	118	\$2.90	SYR-246-03	54	\$3.80
SIF-608-02	115	\$2.30	SIF-831-19	118	\$2.90	SYR-246-04	54	\$6.00
SIF-608-03	115	\$2.40	SIF-831-20	118	\$2.90	SYR-246-05	54	\$6.90
SIF-608-04	115	\$2.50	SIF-831-23	118	\$2.90	SYR-246-06	54	\$14.00
SIF-608-05	115	\$2.60	SIF-831-25	118	\$2.90	SYR-246-07	54	\$26.00
SIF-608-06	115	\$2.70	SIF-831-32	118	\$3.00	SYR-246-08	54	\$22.00
SIF-609-01	115	\$3.10	SIF-831-34	118	\$3.90	SYR-CSD	51	\$71.00
SIF-700-01	99	\$80.00	SIF-832-01	99	\$50.00	SYR-STSD	51	\$71.00
SIF-793-01	124	\$2.00	SIF-834-01	99	\$221.50	SYR-STSD-AC	51	\$153.00
SIF-793-02	124	\$2.00	SIF-844-01	80	\$3.00	SYR-STSD-AS	51	\$153.00
SIF-793-03	124	\$2.00	SIF-844-02	80	\$3.00	SYR-STSD-N	51	\$71.00
SIF-793-04	124	\$3.50	SIF-844-03	80	\$7.00	SYR-STST	51	\$50.00
SIF-794-01	124	\$9.00	SIF-844-04	117	\$25.00	SYR-VCSD	51	\$75.00
SIF-795-01	124	\$80.00	SIF-860-02	117	\$12.00	TB-100	124	\$8.20
SIF-795-02	124	\$80.00	SLS-101-01	92	\$240.00	TB-101	124	\$4.60
SIF-795-03	124	\$90.00	SLS-101-02	92	\$218.00	TB-102	124	\$11.80
SIF-801-01	115	\$20.00	SLS-101-03	92	\$60.00	TB-103	124	\$12.80
SIF-801-02	115	\$20.00	SYR-006-xx	54	\$69.00	TB-104	124	\$7.30
SIF-801-03	115	\$20.00	SYR-007-06	54	\$320.00	TB-200	124	\$8.20
SIF-802-01	116	\$25.00	SYR-007-07	54	\$320.00	TB-201	124	\$4.60
SIF-803-01	125	\$26.00	SYR-007-08	54	\$320.00	TB-202	124	\$11.80
SIF-803-02	125	\$90.00	SYR-007-09	54	\$320.00	TB-203	124	\$12.80
SIF-811-01	116	\$26.00	SYR-007-10	54	\$202.00	TB-204	124	\$7.30
SIF-813-01	116	\$47.00	SYR-007-10P	54	\$63.40	TC3-10	139	\$2.22
SIF-813-02	116	\$19.00	SYR-031-00	52	\$8.80	TC4-10	139	\$2.26
SIF-813-03	116	\$30.00	SYR-031-01	52	\$9.80	TC5-10	139	\$2.31
SIF-813-04	116	\$29.00	SYR-031-02	52	\$8.00	TUB-057-05	139	\$31.36
SIF-814-01	116	\$32.00	SYR-031-03	52	\$8.20	TUB-057-B	139	\$8.32

DENTAL COMPONENTS & TUBING

Index

Catalog No.	Page	Retail Price	Catalog No.	Page	Retail Price	Catalog No.	Page	Retail Price
TUB-058	139	\$9.26	VAF-482-08	80	\$19.30	VAF-521-12	86	\$112.30
TUB-058B	139	\$3.43	VAF-483-01	81	\$14.50	VAF-521-13	86	\$2.60
TUB-059	139	\$1.80	VAF-483-02	81	\$14.00	VAF-521-15	86	\$115.60
TUB-060	139	\$34.00	VAF-483-03	81	\$14.00	VAF-522-01	86	\$66.30
TUB-065	139	\$2.29	VAF-483-04	81	\$20.10	VAF-523-00A	86	\$105.00
VAF-144-017	85	\$400.00	VAF-483-05	81	\$74.70	VAF-523-00W	86	\$105.00
VAF-144-019	85	\$280.00	VAF-483-06	81	\$38.00	VAF-523-02	86	\$80.00
VAF-144-120	85	\$120.00	VAF-484-01	81	\$6.90	VAF-523-03	86	\$47.30
VAF-144-121A	85	\$130.00	VAF-484-02	81	\$9.00	VAF-523-04	86	\$8.00
VAF-144-121W	85	\$130.00	VAF-484-03	81	\$14.70	VAF-524-00	87	\$25.00
VAF-144-122	85	\$280.00	VAF-485-01	82	\$12.25	VAF-524-01	87	\$30.00
VAF-451-01	77	\$13.40	VAF-490-00	82	\$80.00	VAF-524-02	87	\$17.50
VAF-451-03	77	\$13.40	VAF-490-00D	82	\$2.60	VAF-524-03	87	\$3.00
VAF-451-05	77	\$1.25	VAF-490-01	82	\$90.00	VAF-524-03S	87	\$24.00
VAF-451-05	77	\$1.25	VAF-490-01D	82	\$3.20	VAF-524-04	87	\$2.50
VAF-452-01	77	\$13.40	VAF-490-02	82	\$125.00	VAF-524-05	87	\$12.40
VAF-452-03	77	\$13.40	VAF-490-02D	82	\$3.60	VAF-524-06	87	\$20.70
VAF-453-01	77	\$17.00	VAF-490-03	82	\$210.00	VAF-524-15	87	\$12.40
VAF-453-03	77	\$17.00	VAF-490-03D	82	\$2.20	VAF-525-01	87	\$60.00
VAF-453-04	77	\$17.50	VAF-490-04	82	\$5.70	VAF-525-02	87	\$25.50
VAF-454-01	77	\$17.35	VAF-490-05	82	\$9.50	VAF-525-03	87	\$28.00
VAF-454-02	77	\$17.35	VAF-491-01	82	\$150.00	VAF-525-05	87	\$2.95
VAF-454-03	77	\$17.35	VAF-491-02	82	\$130.00	VAF-527-01	87	\$120.00
VAF-454-04	77	\$17.35	VAF-491-03	82	\$180.00	VAF-527-02	87	\$28.00
VAF-455-01	77	\$18.00	VAF-492-02	82	\$13.70	VAF-527-03	87	\$42.00
VAF-455-02	77	\$18.00	VAF-492-05	83	\$24.30	VAF-527-04	87	\$18.00
VAF-461-01	78	\$15.50	VAF-492-06	83	\$40.00	VAF-528-01	87	\$100.00
VAF-462-01	78	\$13.00	VAF-492-07	83	\$5.20	VAF-528-02	87	\$25.50
VAF-463-01	78	\$3.50	VAF-492-08M	83	\$32.00	VAF-528-03	87	\$20.00
VAF-464-01	78	\$2.00	VAF-493-01	83	\$15.30	VAF-528-04	87	\$2.90
VAF-464-03	78	\$2.00	VAF-493-02	83	\$15.30	VAF-528-05	87	\$18.00
VAF-466-01	78	\$22.70	VAF-493-04	83	\$13.70	VAF-529-01	88	\$170.00
VAF-468-01	79	\$12.50	VAF-493-05	84	\$19.00	VAF-532-01	88	\$25.00
VAF-471-01	79	\$21.00	VAF-494-01	84	\$30.00	VAF-532-02	88	\$26.50
VAF-471-02	79	\$21.00	VAF-494-03	84	\$15.00	VAF-535-01	88	\$9.50
VAF-471-03	79	\$21.00	VAF-494-04	84	\$14.00	VAF-535-02	88	\$9.50
VAF-472-01	79	\$22.00	VAF-495-01	84	\$21.00	VAF-536-01	88	\$26.50
VAF-473-01	79	\$11.30	VAF-501-01	84	\$17.00	VAF-536-02	88	\$26.50
VAF-473-02	79	\$11.30	VAF-502-01	84	\$12.80	VAF-536-03	88	\$37.00
VAF-473-03	79	\$20.50	VAF-502-03	84	\$12.80	VAF-536-04	88	\$51.80
VAF-474-01	79	\$10.40	VAF-511-01N	85	\$384.00	VAF-537-01	88	\$302.80
VAF-474-02	79	\$10.40	VAF-511-03N	85	\$283.50	VAF-537-02	88	\$366.70
VAF-475-01	80	\$50.90	VAF-511-04N	85	\$82.70	VAF-537-03	89	\$84.00
VAF-475-02	80	\$42.50	VAF-511-06N	85	\$18.40	VAF-537-04	89	\$90.00
VAF-475-03	80	\$9.00	VAF-512-01N	85	\$304.00	VAF-537-05	89	\$40.00
VAF-476-01	80	\$34.00	VAF-512-02N	85	\$218.25	VAF-537-06	89	\$16.00
VAF-482-01	80	\$55.00	VAF-512-03N	85	\$82.70	VAF-537-07	89	\$96.00
VAF-482-02	80	\$19.90	VAF-512-04N	85	\$18.40	VAF-538-01	89	\$170.00
VAF-482-02E	80	\$20.60	VAF-514-01	85	\$63.30	VAF-538-02	89	\$170.00
VAF-482-03	80	\$35.30	VAF-521-01	86	\$42.00	VAF-538-03	89	\$20.00
VAF-482-04	80	\$18.20	VAF-521-01N	86	\$35.00	VAF-539-01	89	\$980.00
VAF-482-05	80	\$19.00	VAF-521-02	86	\$18.00	VAF-539-02	89	\$1,300.00
VAF-482-06	80	\$23.60	VAF-521-10	86	\$42.00	VAF-540-01	89	\$41.20
VAF-482-07	81	\$54.00	VAF-521-11	86	\$42.00	VAF-540-02	89	\$82.40

DENTAL COMPONENTS & TUBING

Index

Catalog No.	Page	Retail Price
VCF-200-01	59	\$290.00
VCF-200-01P	59	\$595.00
VCF-200-01S	59	\$620.00
VCF-200-02	59	\$355.00
VCF-200-03	59	\$190.00
VCF-201-01	59	\$290.00
VCF-201-02	59	\$355.00
VCF-201-03	59	\$190.00
VCF-202-01	60	\$280.00
VCF-202-02	60	\$345.00
VCF-202-03	60	\$180.00
VCF-300-01D	59	\$255.00
VCF-300-01E	59	\$245.00
VCF-301-01	61	\$34.00
VCF-301-02	61	\$34.00
VCF-302-00	60	\$68.00
VCF-302-01	60	\$45.00
VCF-302-02	60	\$47.00
VCF-303-00AL	60	\$40.00
VCF-303-00DL	60	\$30.00
VCF-303-01	60	\$42.00
VCF-305-00	60	\$22.30
VCF-305-01	60	\$26.50
VCF-305-02	60	\$30.20
VCF-305-03	60	\$30.20
VCF-306-01	61	\$49.00
VCF-306-02	61	\$13.00
VCF-306-03	61	\$20.00
VCF-307-01	61	\$46.00
VCF-307-02	61	\$13.00
VCF-307-03	61	\$20.00
VCF-308-01	61	\$43.00
VCF-308-02	61	\$11.00
VCF-308-03	61	\$20.00
VCF-309-00	61	\$40.00
VCF-309-01	61	\$36.00
VCF-311-01	61	\$25.00
VCF-311-02	61	\$15.00
VCF-311-03	62	\$24.00
VCF-311-04	62	\$24.00
VCF-312-01	62	\$18.00
VCF-312-02	62	\$18.00
VCF-313-01B	62	\$11.90
VCF-313-01G	62	\$19.00
VCF-314-01	62	\$7.00
VCF-314-02	62	\$8.00
VCF-314-03	62	\$4.00
VCF-314-04	62	\$5.00
VCF-315-01	62	\$70.00
VCF-315-02	62	\$64.00
VCF-315-03	62	\$64.00
VCF-316-01	62	\$20.60
VCF-316-02	62	\$20.60
VCF-316-03	62	\$24.70

Catalog No.	Page	Retail Price
VCF-316-04	62	\$20.60
VCF-316-05	62	\$24.70
VCF-316-06	62	\$24.70
VCF-320-01	63	\$8.80
VCF-320-02	63	\$5.20
VCF-321-01	63	\$2.90
VCF-321-02	63	\$2.90
VCF-322-01	63	\$14.00
VCF-322-02	63	\$4.50
VCF-323-01	63	\$0.60
VCF-323-02	63	\$9.00
VCF-325-01	63	\$3.50
VCF-325-02	63	\$4.00
VCF-326-01	64	\$12.00
VCF-331-04	65	\$12.50
VCF-331-04B	65	\$34.00
VCF-332-01	64	\$24.00
VCF-332-02	65	\$11.45
VCF-333-01	64	\$24.00
VCF-333-02	65	\$12.50
VCF-333-02B	65	\$34.00
VCF-333-05	65	\$49.45
VCF-333-06	65	\$50.00
VCF-333-07	65	\$50.00
VCF-334-01	65	\$0.26
VCF-334-02	65	\$0.34
VCF-335-01	64	\$22.00
VCF-336-01	64	\$22.00
VCF-337-01	65	\$11.45
VCF-338-01	65	\$3.00
VCF-338-02	65	\$3.00
VCF-339-00	65	\$25.00
VCF-339-01	65	\$16.50
VCF-339-02	65	\$14.00
VCF-342-01	66	\$64.00
VCF-342-02	66	\$64.00
VCF-345-01	66	\$19.00
VCF-345-02	66	\$24.00
VCF-345-03	66	\$24.00
VCF-351-01	66	\$9.30
VCF-351-02	66	\$170.00
VCF-352-01	66	\$14.00
VCF-352-02	66	\$14.00
VCF-353-01	66	\$8.70
VCF-353-02	66	\$8.70
VCF-354-01	67	\$190.00
VCF-354-02	67	\$189.50
VCF-354-03	67	\$9.30
VCF-354-04	67	\$6.70
VCF-355-01	67	\$215.00
VCF-355-02	67	\$235.00
VCF-361-01	67	\$160.70
VCF-362-01	68	\$49.50
VCF-362-02	68	\$32.50

Catalog No.	Page	Retail Price
VCF-362-03	68	\$91.70
VCF-372-01	68	\$370.00
VCF-372-01B	68	\$30.00
VCF-372-01B	68	\$30.00
VCF-372-02	68	\$370.00
VCF-374-01	68	\$10.00
VCF-374-02	68	\$11.00



tel: 818.502.9954

fax: 818.502.0501

sales@dentalpartsusa.com

www.DentalPartsUSA.com

418 North Glendale Ave, Unit J
Glendale, CA 91206, USA